EXODUS

EXODUS

AUTOMATIC EXOSOME ISOLATION SYSTEM





Automatic System for Exosomes Isolation



EXODUS is an automatic, label-free, and highly efficient exosome isolation system. With EXODUS, you can easily and quickly isolate high-quality, intact exosomes with excellent yield and purity from a variety of bio-fluids and sample volumes.


Experience the efficiency of EXODUS for yourself and take your research to the next level.




Isolation Principles

EXODUS has been developed using a dual-membrane nanofiltration system that integrates periodic negative pressure oscillation (NPO) and double-coupled ultrasonic harmonic oscillations (HO).



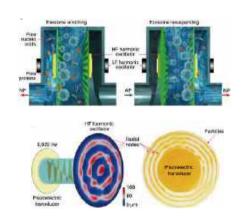

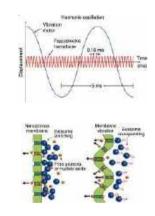


Nature Methods, 2021, 18(2):212-218.

EXODUS can rapidly remove free nucleic acid and protein impurities from the sample, resulting in the efficient purification and enrichment of exosomes. The exosomes are precisely intercepted by nanoporous membrane, allowing for a highly targeted isolation process.

EXODUS has great potential to revolutionize exosome isolation and drive new discoveries in biomedical research and translation.

Automatic

EXODUS is designed to automatically isolate high yield and purity exosomes from different biofluid sample volumes.

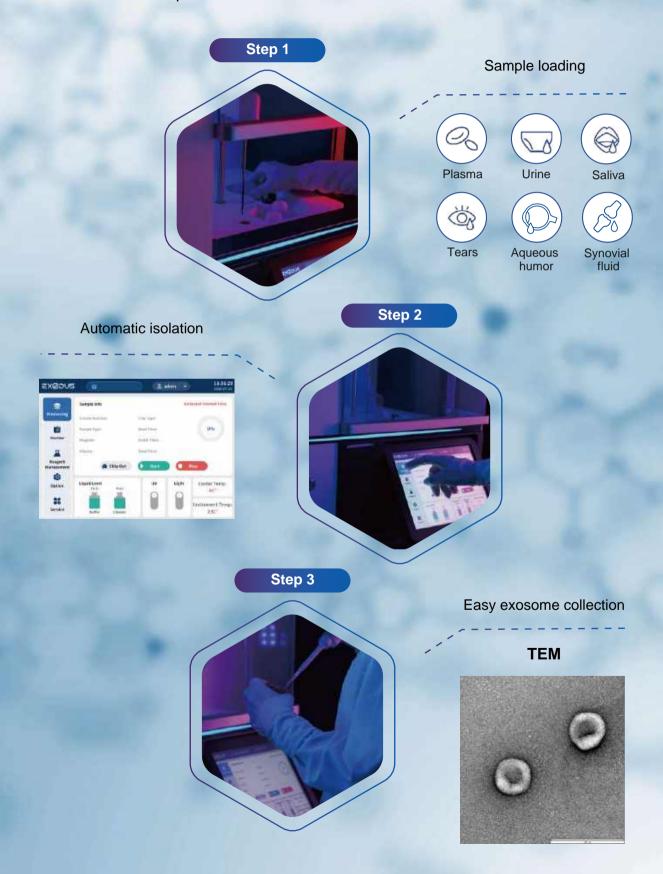







Rapid isolation

Maximum isolation speed: 200 mL/h

High purity and high yield

Purity ~ 99 %; Yield ~ 90 %



Wide application

Sample types	Sample volumes	Sample types	Sample volumes
Urine		Plasma	0.01 - 2 mL
Plant		Saliva	0.5 - 10 mL
Cell culture medium	1 - 250 mL	Tears	0.005 - 1 mL
Cell-derived vesicle		Aqueous humor	0.005 - 1 mL
Bacterial culture medium		Cerebrospinal fluid	0.5 - 25 mL

80


Label-free

Only need PBS buffer



10 mL urine




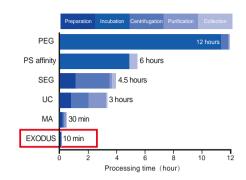

Nutaure Methods, 2021, 18(2):212-218.



High purity and high yield

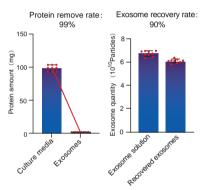

Nutaure Methods, 2021, 18(2):212-218.

Other small amount sample types



Wide application





Plasma

Cerebrospinal

fluid

Tears

Saliva

Aqueous

humor




Synovial fluid

TEM image of exosome



Label-free

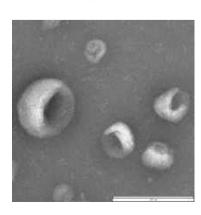



Various Sample Types



Plasma



Urine



Saliva



Cerebrospinal fluid



Tears



Aqueous humor



Synovial fluid



Tissue



Cell culture medium



Bacterial culture medium



Cell-derived vesicle



Plant

Applications


Early diagnosis

- Drug delivery
- Exosome therapeutics
- Regenerative medicine

Urine EV Isolation with EXODUS

Figure 1. Characterization of exosomes by transmission electron microscopy

TEM of exosomes harvested from urine by EXODUS, showing the characteristic cup and plate shape of exosomes. Scale bar = 600 nm

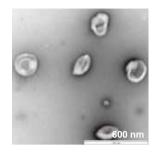


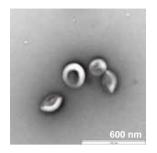


Figure 2. Method comparison for exosome isolation.

Western blot analysis of exosomal markers (Alix, CD63, TSG101, and CD81) and impurity Uromodulin (UMOD) of exosomes isolated using methods of EXODUS, ultracentrifugation (UC), PEG precipitation, Membrane affinity (MA), Size exclusion chromatography (SEC), and phosphatidylserine (PS) affinity, respectively. Compared to other methods, exosome isolated by EXODUS shows higher relative yield and purity in WB gel graph (left) and bar graph of normalized band intensity (right).

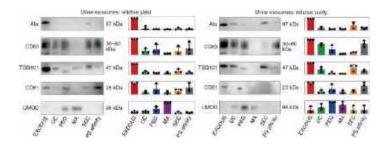


Nature Methods, 2021, 18(2):212-218.

2 Cell culture medium EV Isolation with EXODUS

Figure 1. Characterization of exosomes by TEM.

Electron micrographs of exosomes isolated from umbilical cord mesenchymal stem cell supernatants using EXODUS H (left) and UC (Right). TEM of exosomes obtained from EXODUS shows such typical cupped structure, more exosomal vesicular particles, and less particle aggregation compared to that from UC. Scale bar = 600 nm.



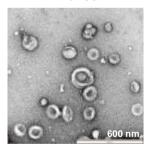

UC


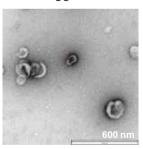


Figure 2. Scratched wound assay for the effects evaluation of exosomes on

Cells began to exhibit significantly migration 12h after treatment with exosomes isolated either from EXODUS or EXO-kit compared to that from UC and the control group. After 24h, enhanced cell migration was observed from the group with the addition of EXODUS-isolated exosome. Scale bar = 40 µm.

migration capacity of cells.

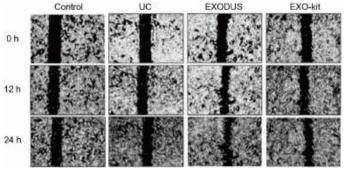



Electrophoresis, 2024, 0173-0835.

3 Plasma EV Isolation with EXODUS

Figure 1. Characterization of exosomes by TEM.

TEM analysis of exosomes isolated from the plasma of patients with esophageal squamous cell carcinoma using the EXODUS reveals their distinctive cup-shaped and plate-like morphology. Scale bar = 400 nm.

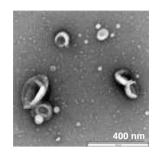


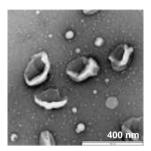

Technol Cancer Res Treat, 2024, 1533-0338.

Figure 2. The yield and purity analysis of sEVs isolated by methods of EXODUS, UC, and PEG, respectively.

sEVs isolated using EXODUS or PEG from 20 μ L of plasma (approximately 2 x 10 10 particles/mL) showed over a 50-fold increase compared to those isolated using UC. The purity of sEVs isolated using EXODUS was more than three times higher than that achieved with UC or PEG.

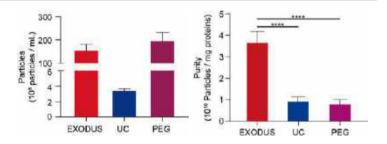


Int J Nanomedicine, 2024, 19:1351-1362.


4 Trace sample EV Isolation with EXODUS

Figure 1. Characterization of tears EVs and aqueous humor EVs by TEM.

TEM images of exosomes derived from two trace samples of 50 μ L tear fluid (left, scale bar = 100 nm) and 150 μ L aqueous humor (right, scale bar = 400 nm) using the EXODUS, reveals their characteristic cup-shaped and plate-like morphological features.

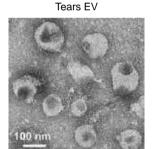

Aqueous Humor EV

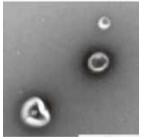


ACS Nano, 2022, 16(8): e11720.

Figure 2. The yield and purity analysis of tear exosome isolated by methods of EXODUS, UC, SEC, and PEG.

A. Total particles of exosomes isolated using EXODUS or PEG from 50 µL of tear (approximately 2 x 10¹⁰ particles/mL) showed over a 20-fold and 2-fold increase compared to those isolated from UC and SEC, respectively.

B. Equal-protein-mass (3 µg) western blot analysis of the exosomal markers (Alix, CD63, CD9, and CD81) and negative marker (Calnexin) of tear exosomes prepared by iTEARS EXODUS and other methods.

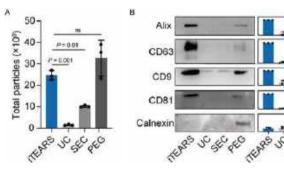

ACS Nano, 2022, 16(8): e11720.

5 Bacterial EV Isolation using with EXODUS


Figure 1. Characterization of bacterial EVs by TEM.


TEM images of EVs derived from Helicobacter pylori, Escherichia coli, Mycobacterium avium, and intestinal microbiota isolated by EXODUS H 600. The EXODUS H 600 is capable of isolating EVs from a variety of bacterial sources, including both Gram-positive and Gram-negative bacteria.

Helicobacter pylori EVs 1000 nm 600 nm 200 nm 1000 nm 600 nm 200 nm Mycobacterium avium EVs Gut microbiota EVs 1000 nm 400 nm 200 nm 200 nm

Figure 2. Comparison of yield and purity by using EXODUS and UC.

EXODUS achieves higher particle yield and EV purity compared to UC of *E.coli*

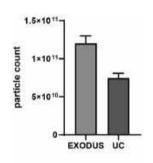


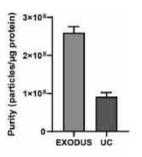


6 Virus and virus-like particles Isolation with EXODUS

Figure 1. TEM characterization of virus and virus-like particles.

Using EXODUS, high-resolution TEM images were obtained to highlight typical morphology of bacteriophages, phage protein self-assembly particles and Rift Valley fever virus. Gentle isolation of EXODUS is able to preserves the viral or particle integrity, which is crucial for accurate structural analysis. Scale bar = 100 nm.

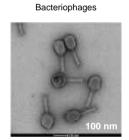

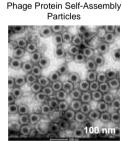



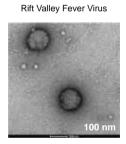


Figure 2.Lentiviral titer recovery rate results


After purification with the EXODUS, the titer of lentiviral infection increased from 9.86×10^7 to 1.56×10^9 , with a recovery rate of 63.29%, and the physical titer increased from 2.5×10^{10} to 3.95×10^{11} , with a recovery rate of 63.0%

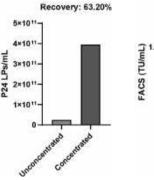

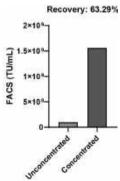

Publications

NO.	Title	Journal	Five-year average IF
01	Exosome detection via the ultrafast-isolation system: EXODUS	Nature Methods	45.6
02	UPCARE:Urinary Extracellular Vesicles-Derived Prostate Cancer Assessment for Risk Evaluation	Journal of Extracellular Vesicles	19.6
03	Discovering the Secret of Diseases by Integrated Tear Exosomes Analysis via Rapid-isolation System: iTEARS	ACS Nano	16.2
04	Robust Acute Pancreatitis Identification and Diagnosis: RAPIDx	ACS Nano	16.2
05	Interaction network of extracellular vesicles building universal analysis via eye tears: iNEBULA	Science Advances	13.7
06	Sensitive small extracellular vesicles associated circRNAs analysis combined with machine learning for precision identification of gastric cancer	Chemical Engineering Journal	13.2
07	Quantitative metabolic analysis of plasma extracellular vesicles for the diagnosis of severe acute pancreatitis	Journal of Nanobiotechnology	11.4
08	Lipidomic identification of urinary extracellular vesicles for non-alcoholic steatohepatitis diagnosis	Journal of Nanobiotechnology	11.4
09	Metabolomic investigation of urinary extracellular vesicles for early detection and screening of lung cancer	Journal of Nanobiotechnology	11.4
10	The genetic source tracking of human urinary exosomes	PNAS	10.8
11	Identification and detection of plasma extracellular vesicles-derived biomarkers for esophageal squamous cell carcinoma diagnosis	Biosensors & Bioelectronics	9.9
12	Audible Acoustic Wave Promotes EV Formation and Secretion from Adherent Cancer Cells via Mechanical Stimulation	ACS Applied Materials & Interfaces	8.7
13	Identification of circulating extracellular vesicle long RNAs as diagnostic biomarkers for patients with severe acute pancreatitis	Clinical and Translational Medicine	8.0
14	Prediction of Response to Chemoradiotherapy by Dynamic Changes of Circulating Exosome Levels in Patients with Esophageal Squamous Cell Carcinoma	International Journal of Nanomedicine	7.5
15	Investigating the proliferative inhibition of HepG2 cells by exosome-like nanovesicles derived from Centella asiatica extract through metabolomics	Biomedicine & Pharmacotherapy	6.8
16	Gut Subdoligranulum variabile ameliorates rheumatoid arthritis by promoting TSG-6 synthesis from joint cells	Frontiers in Immunology	6.8
17	Metabolomic analysis of exosomal-markers in esophageal squamous cell carcinoma	Nanoscale	6.1
18	Metabolic signatures of tear extracellular vesicles caused by herpes simplex keratitis	Ocular Surface	5.9
19	Sensitive electrochemical biosensor for rapid detection of sEV-miRNA based turbo-like localized catalytic hairpin assembly	Analytica Chimica Acta	5.5
20	Isolation of Exosome Nanoparticles from Human Cerebrospi- nal Fluid for Proteomic Analysis	ACS Applied Nano Materials	5.4
21	Human umbilical cord mesenchymal stem cells inhibit liver fibrosis via the microRNA-148a-5p/SLIT3 axis	International Immunopharmacology	5.0
22	Patient-derived induced pluripotent stem cells with a MERTK mutation exhibit cell junction abnormalities and aberrant cellular differentiation potential	World Journal of Stem Cells	4.2
23	Assessing Alzheimer's disease via plasma extracellular vesicle–derived mRNA	Alzheimer's & Dementia: Diagnosis, Assessment & Disease Monitoring	4.0
24	Efficient preparation of high-purity and intact mesenchymal stem cell–derived extracellular vesicles	Analytical and Bioanalytical Chemistry	3.8
25	Plasma-derived Exosomal miR-25-3p and miR-23b-3p as Predictors of Response to Chemoradiotherapy in Esophageal Squamous Cell Carcinoma	Technology in Cancer Research & Treatment	2.8
26	Comparative investigation of exosome extraction from rat bone marrow mesenchymal stem cells using three different methodologies	Electrorhoresis	2.8
27	Isolation of small extracellular vesicles from a drop of plasma via EXODUS and their fingerprint proteomics profiling by MALDI-TOF MS	Biosensors & Bioelectronics:X	/
28	Application of EXODUS system combined with allosteric DNA nanoswitches in the detection of miR-107 among plasma exosomes of Parkinson's disease patients	Chinese Journal of Preventive Medicine	/

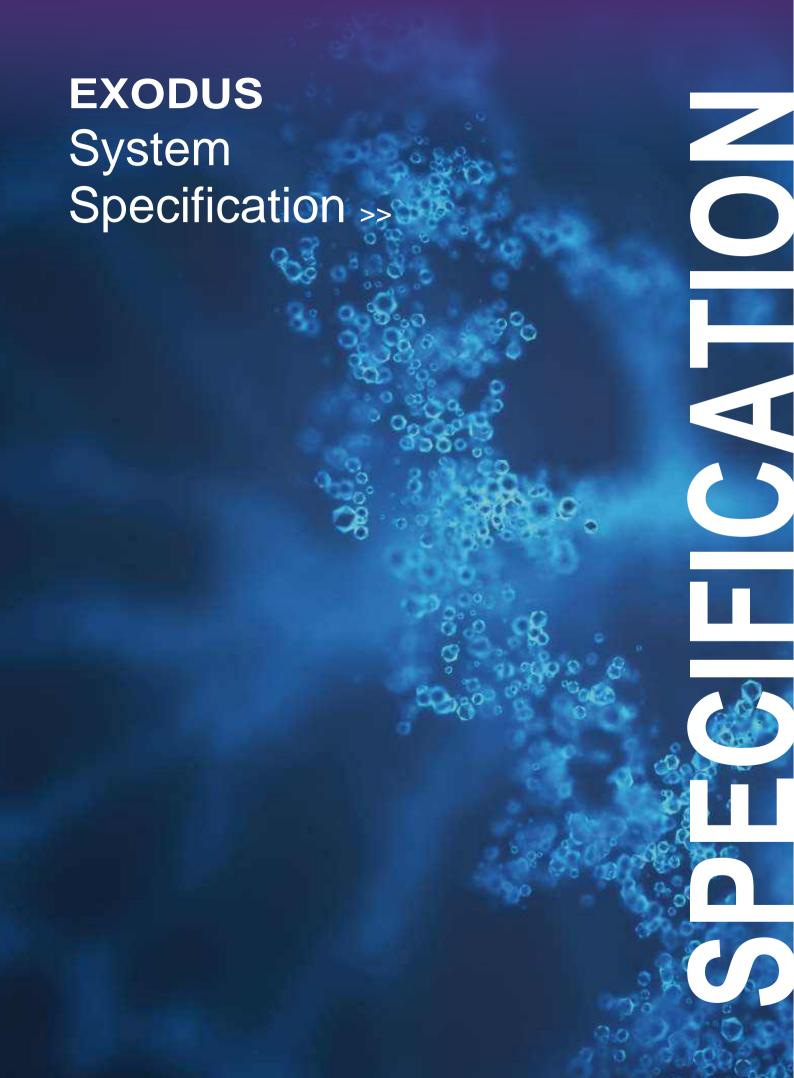

Model	EXODUS H-300	EXODUS H-600	
Isolation principles	Combination of the negative pressure oscillations (NPO) and double coupled harmonic oscillations (HO) on nanoporous membrane		
Sample types	Plasma, urine, saliva, cerebrospinal fluid, tears, aqueous humor, synovial fluid, tissue, cell culture medium, bacterial culture medium, cell-derived vesicle, plant, ect.		
Isolation device size	S/M	S/M/L	
Temperature of sample reservoir	2 - 8 °C		
Sample volumes	10 μL - 50 mL	10 μL - 250 mL	
Processing speed	Max speed 50 mL/h	Max speed 200 mL/h	
Isolation data saving	2000	20000	
Exosome recovery volumes	100 - 400 μL	100 - 1000 μL	
Ultraviolet sterilization	Internal UV lamp, turn off automatically after 30 min		
Display	10.4 inch touch screen, real time display with sample type, time, processing information ect. Supporting the operation without computer		
Dimension	535 x 510 x 475 mm (H x W x D)		
Net weight	40 kg (88 lbs)		
System interfaces	4 USB ports, 1 network port, 1 serial port		

EXODUS

Product specifications may change without notice, based on the latest technical data and test results.

★ service@exodus-bio.com

• 160 E Tasman Dr., San Jose, CA 95134, United States