


BG-gdsAUTO710 Mini Chemiluminescence Imaging Analysis System User Manual

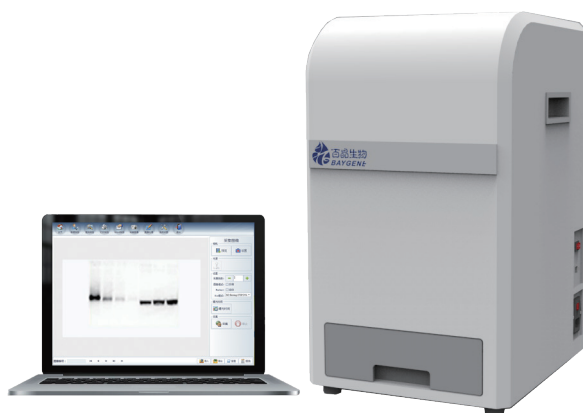

Please read this manual carefully before operating the instrument.
Please keep this manual properly so that it can be used as needed.

Revision: September 2019 V1.1

© 2018 Beijing Baygene Biotechnology Co., Ltd. All rights reserved.

Beijing Baygene Biotech Co., Ltd. ensures that the information contained in its supporting documents is accurate and clear, but assumes no responsibility for any errors or omissions. Beijing Baygene Biotech Co., Ltd.'s products and services are constantly developing. Ensure that any published information used for reference is current and relevant to the current condition of the product. If necessary, please contact Beijing Baygene Biotech Co., Ltd.

This manual may not be copied, photocopied, reproduced, translated or converted into other electronic or machine-readable versions without the written approval of Beijing Baygene Biotech Co., Ltd. in advance.

All information contained in this manual is proprietary, confidential and belongs to Beijing Baygene Biotech Co., Ltd. This manual is protected by copyright and reproduction is prohibited. This manual is for personal use only approved by Beijing Baygene Biotech Co., Ltd.

Contact information:

Address: Building 7, Science and Technology Pioneer Park, Zone B, Tianzhu Airport Industrial Zone, Beijing

Production address: No. 28, Yuhua Road, Shunyi District, Beijing (No. 7-8, Airport Area B)

Tel: 010-80483100/80483200

Fax: 010-80482859

Website: www.baygenebiotech.com.cn

E-mail: info@baygenebiotech.com

In this document, the BG-gdsAUTO710 MINI Chemiluminescence Imaging Analysis System is abbreviated as BG-gdsAUTO710.

CONTENT

Chapter 1 System Introduction	1
1.1 Instrument Introduction	1
1.2 Performance Introduction	1
1.3 Instrument Parameters	1
1.4 Software Installation Requirements	1
Chapter 2 Software Installation	2
2.1 Connect The Instrument	2
2.2 Install The Driver	2
2.3 Install The Software	2
2.4 Software Registration	3
Chapter 3 System Operation	4
3.1 Instrument Operation	4
3.2 Software Startup	4
Chapter 4 Report Management	5
4.1 New Report	5
4.2 Save Report	6
4.3 Query Report	6
4.4 Print Report	6
4.5 Word Report Output	7
4.6 Close Report	7
Chapter 5 Image Capture	8
5.1 Camera Settings	8
5.2 Light Source Operation	9
5.3 Capture Parameter Settings	9
5.4 Capture Image	10
5.5 Save And Recall Parameters	10

CONTENT

Chapter 6 Image	11
6.1 Selection Tool	11
6.2 Zoom, Rotate, Move Tools	12
6.3 Editing Tools	12
6.4 Mark Tools	12
6.5 Adjustment Tools	13
6.6 Image Merge Function	16
6.7 File Toolbar	16
Chapter 7 System Settings	17
7.1 User Management	17
7.2 Change Password	18
7.3 Report Title	18
7.4 Sample Type	18
7.5 Device Settings	19
7.6 System Information	19
7.7 Other	19
Chapter 8 Other Settings	20
8.1 Image Import	20
8.2 Image Export	20
8.3 System Backup And Recovery	20
Chapter 9 Precautions For Use	21
Daily Maintenance And Care	21
Chapter 10 Frequent Q & A For Use	21
Chapter 11 After Sales	22
11.1 Product Warranty Terms	22
11.2 After-sales Service Procedures	22

Chapter 1 System Introduction

1.1 Instrument Introduction

BG-gdsAUTO710MINI chemiluminescence imaging analysis system combines high-definition image capture, processing, database management and other functions. It is a chemiluminescence imaging analysis system specifically designed for Western Blot protein electrophoresis experiments to help researchers quickly capture the accurate protein images and analyze the results.

This system includes BG-gdsAUTO710MINI chemiluminescence imager and ChemiImaging imaging software (hereinafter referred to this software).

1.2 Performance Introduction

1. Efficient and automatic image capture and processing system;
2. Simple information entry method and fast data management system;
3. Friendly interface to generate report form automatically;
4. Large data storage capacity to realize the scientific and perfect image and information management for users.

1.3 Instrument Parameters

- > High-resolution professional cooled CCD;
- > Large target surface CCD chip;
- > Large aperture and fixed focus lens;
- > High sensitivity based on high light utilization efficiency;
- > Built-in LED light source;
- > ChemiImaging professional imaging software, switch between Chinese / English, easy operation; Free upgrade for life.
- > Wide voltage operation: 85-265V, suitable for different countries and regions;
- > Dimensions: 320x360x550mm (LxWxH);
- > Instrument power: 30W.

1.4 Software Installation Requirements

- > The requirements of the computer for normal system operation are as follows:
- > CPU: Core2 Duo or above is recommended
- > Memory: 2GB or more is recommended
- > Hard disk: more than 120G
- > Hi-Speed USB2.0 or USB3.0
- > Display resolution: 1024x768
- > Operating system: WinXP 32 / 64-bit, Win7 32 / 64-bit, Win10

Note: Specifications and designs may change at any time without notice.

Chapter 2 Software Installation

2.1 Connect The Instrument

When confirming that the power button of the instrument is off, insert the USB cable and power cable in the interface, and then turn on the power switch.

2.2 Install The Driver

2.2.1 Camera driver:

Open path

: \ Chemilming Software \ Drives \ Cammera
\ Camera Driver.exe, double-click the installer,
the installation interface pops up, and follow the
program guides to install. (Figure 2-1)

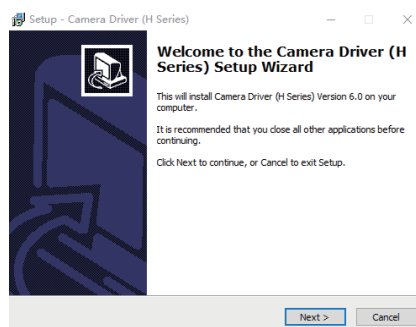


Fig 2-1

2.2.2 Serial driver:

Open path

: \ Chemilming Software \ Drives \ USBtoCom\
HL-340, double-click the installer, the installation
interface pops up, and click the <INSTALL> button to
install. (Figure 2-2)

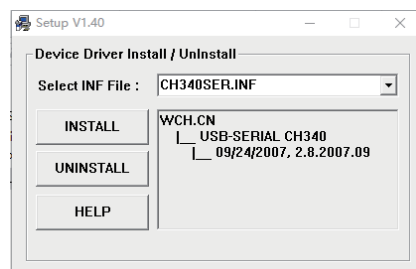


Fig 2-2

2.3 Install The Software

Open the path: \ Chemilming Software \ Setup \
Chemilming.exe, double-click the installation
software, pop up the installation interface, and follow
the program to guide the installation. (Figure 2-3)

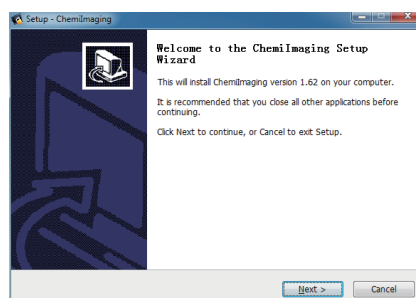


Fig 2-3

Note: To avoid accidental data loss, it is best not to install it in the system disk when selecting the installation path. Please click "Browse" to change the installation path. The change page is shown in Figure 2-4.

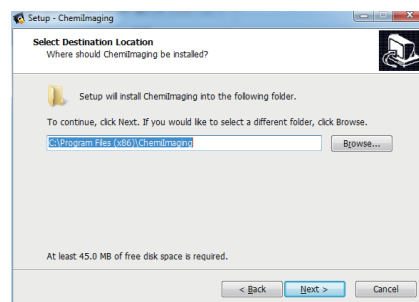


Fig 2-4

2.4 Software Registration

Start the software for the first time and ask for registration (see Figure 2-5).



Fig 2-5

Click the <Register> button, a window will pop up (see Figure 2-6). Copy the machine code.

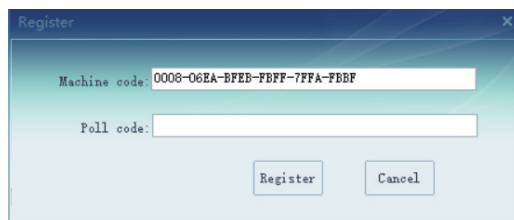


Fig 2-6

Open path
: \ ChemImaging software \ Keygen \
Keygen.exe, double-click to pop up the
interface, as shown in Figure 2-7.

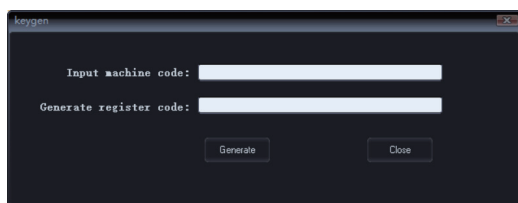


Fig 2-7

In the Enter Machine Code box, enter the
machine code in Figure 2-6, and click the
<Generate Registration Code> button (see
Figure 2-8).

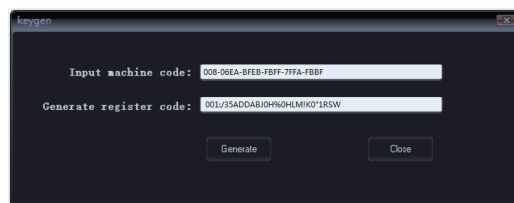


Fig 2-8

Copy the registration code in Figure 2-8 into the registration code box in Figure 2-6, and click the <Register> button to complete the registration.

Chapter 3 System Operation

3.1 Instrument Operation

This instrument is a dark box imaging system. When you use it, you only need to put in experimental samples, and you can perform exposure imaging on the built-in computer.

The steps are as follows:

1. Plug the instrument into a power source and turn on the main power switch on the back of the case.
2. Connect the computer USB data cable.
3. Open Chemilming.exe on the computer and log in, create a new report (database management of this software), enter the relevant information, capture the image (enter the capture image module), set the corresponding parameters, click capture to start imaging.
4. The software has the functions of image capture and image processing, which can process the images captured or imported from other external source, and print reports;
5. When finish, please remove the sample and clean the box;
6. Turn off the main power on the back of the case.

3.2 Software Startup

Double click the desktop or click the start menu icon  to enter the login interface, as shown in Figure 3-1.

Initial value: user name: admin; Password: admin

Once the user has correctly selected or entered a username and password, check the <Remember> item, the system will automatically remember the password.

Click <Login> to enter the system. The interface after login, (as shown in Figure 3-2):



Fig 3-1

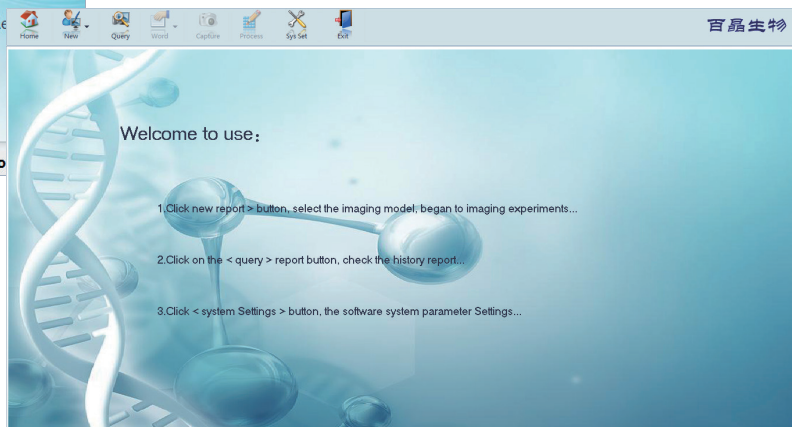


Fig 3-2

Chapter 4 Report Management

This chapter mainly describes how users manage the reporting process after logging in to the system.

4.1 New Report

After entering the system, click the <New Report> button. Select the shooting type (see Figure 4-1).

The system provides three modes:

Singel mode: 1 sample image + 1 marker image + 1 merge image
 Series mode: 5 sample images + 1 marker image

Manual mode:

- a. n sample images
- b. 1 sample image + 1 marker image + 1 merge image

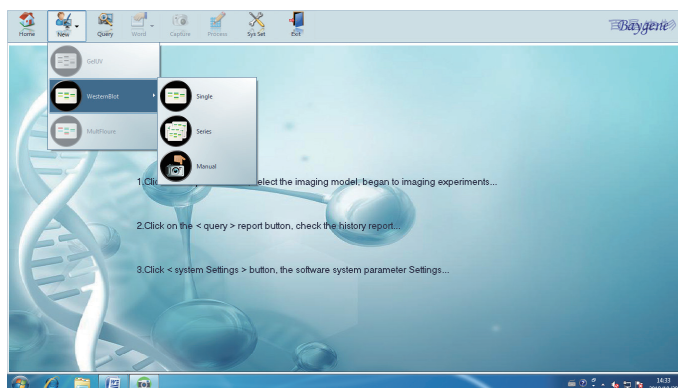


Fig 4-1

After selecting the mode, enter the interface to capture the image (see Figure 4-2).

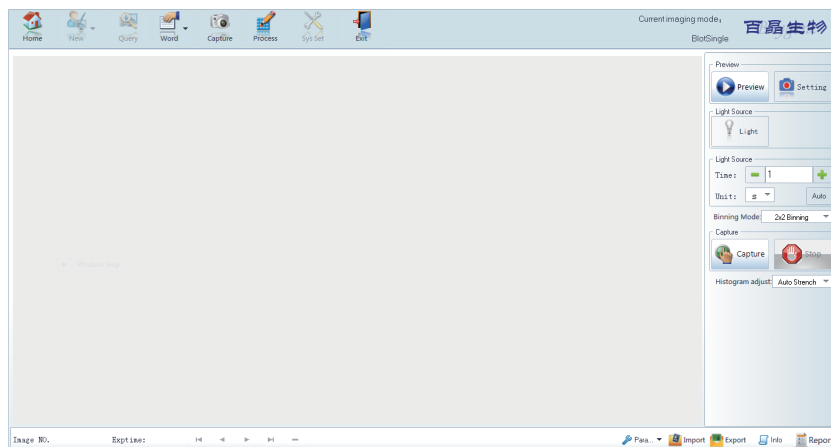


Fig 4-2

Click on the report in the lower right corner, the interface for adding report information will appear, Users can enter relevant information on this interface (see Figure 4-3). In the report information such as the date and user is automatically generated by the system. Users can modify as need.

Fig 4-3

4.2 Save Report

After entering the information, the user can click the <Save> button to save the report, which is convenient for later report search, printing, and output.

4.3 Query Report

Click the <Query Report> button to pop up the query report window (as shown in Figure 4-4), and query according to the relevant search conditions.

Set the query conditions

Sample NO.

Sample Type

Operator

Date

2019/10/19

To

2019/10/29

clear

CaptureModel	Samplpe NO.	Sample Type	Operator	Record Date
▶ BlotSingle	1	westblot	admin	2019/10/29
BlotSingle			admin	2019/9/24
BlotSingle			admin	2019/9/24
BlotSingle			admin	2019/9/24

1/4

Back

delete

open

Query

Fig 4-4

4.4 Print Report

Click the <Print Report> button under the <Output Report> button to pop up the interface, as shown in Figure 4-5.

The preview print effect is displayed. Click <Print> to print the image report.

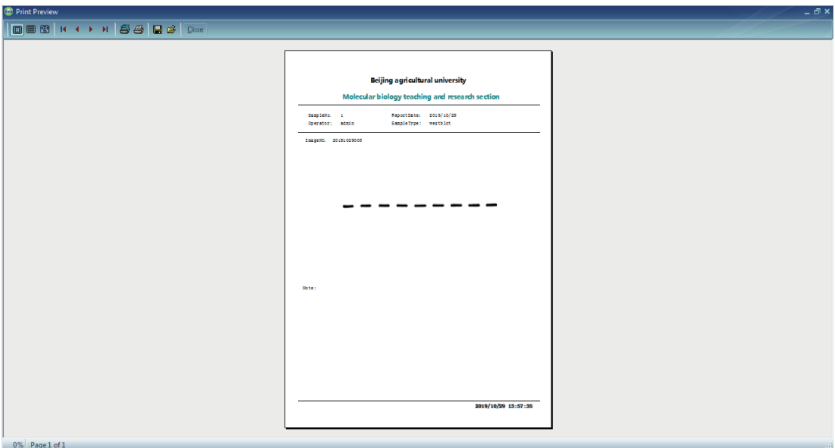


Fig 4-5

6

4.5 Word Report Output

Click the <Word Report> button under the Export Report button to output the report to a Word document.

(Figure 4-6)

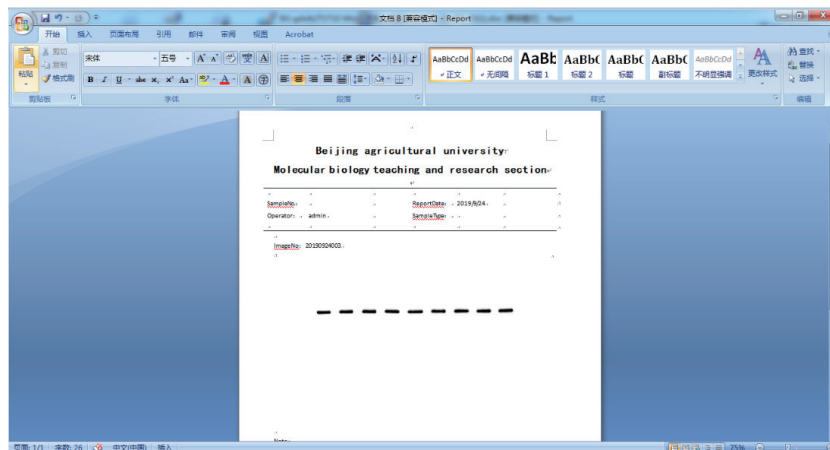


Fig 4-6

4.6 Close Report

Click the <Home> Close this report.

Chapter 5 Image Capture

This chapter mainly introduces the use and information of the image capture module. Place the prepared sample membrane in the center of the sample stage.

After creating a report, go directly to the selected capture interface.

Note: If a report has been created, open the report and click <Capture > to enter the interface.

The system provides three image capture modes (see Figure 5-1, Figure 5-2, and Figure 5-3).

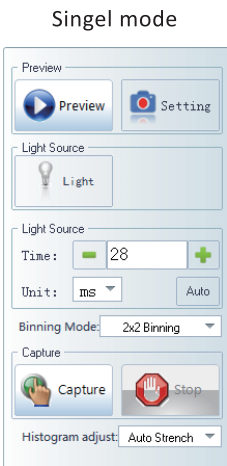


Fig 5-1

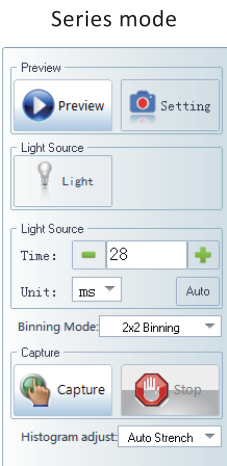


Fig 5-2

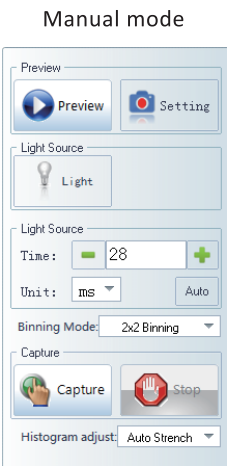


Fig 5-3

5.1 Camera Settings

Click the <Preview> button:

When the button changes to  , it indicates that the image is being previewed;

When the button changes to  , the preview stopped.

Click the <Setting> button: the camera parameter setting interface pops up (see Figure 5-4).

After the user sets the parameters as required, please click the <Save> button.

Click the <Call> button to call the factory default parameters.

Click the <Cover> button to override the factory default parameters.



Fig 5-4

5.2 Light Source Operation

The instrument has a built-in LED light source. Click the <Light> button. The button color turns green, which indicates that the light source is turned on, while the button color turns gray which indicate that the light source is turned off.

5.3 Capture Parameter Settings

The instrument has a built-in LED light source. Click the <Light> button. The button color turns green, which indicates that the light source is turned on, while the button color turns gray which indicate that the light source is turned off.

5.3.1 Singel mode: (Figure 5-5)

Time: <-> and <+> buttons to set the exposure time

Unit: Click to switch the unit of exposure time to m, s or ms

Auto: The software automatically calculates the exposure time and recommended to users.

A window will pop up after clicking (see Figure 5-6).

Binning mode: This is an image output pattern. The charges from adjacent pixels are added together and output in one-pixel pattern. Although lower image resolution, binning mode provides higher image sensitivity and output speed. The user can choose how many pixels binning by this mode to get the ideal sensitivity.

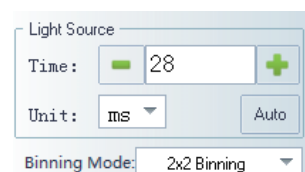


Fig 5-5

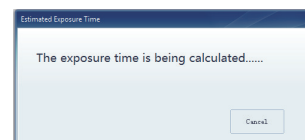


Fig 5-6

5.3.2 Series mode:(Figure5-7)

Select exposure time program ExpTime P4

Click the time value to change the exposure time:

Click <s> to switch the exposure time unit:

10 40 60 80 100 s

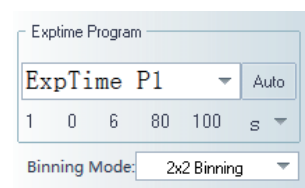


Fig 5-7

5.3.3 Manual mode:(Figure 5-8)

Number of captures: <-> and <+> buttons to set the number

Marker: ☒ Auto Whether to capture a Marker image.

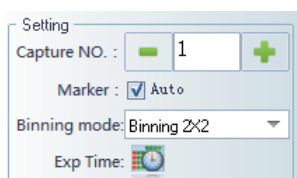


Fig 5-8

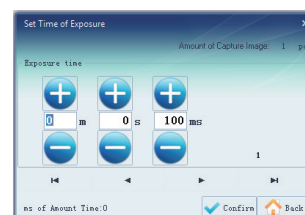


Fig 5-9

Exp Time: A window for setting the exposure time (see Figure 5-9):

5.4 Capture Image

Capture image: Click the <Capture> button, the instrument will capture the image, and the button will turn gray.

A pop-up window, as shown in Figure 5-10, displays the progress of image capture and the estimated time required.

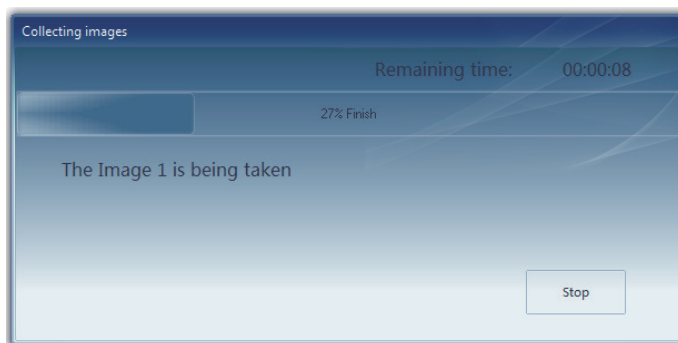


Fig 5-10


Stop the capture: If the user needs to stop the capture midway, please click the <Stop> button.

Note: When the image is being collected, please do not operate other parts because the camera consumes a lot of calculation memory. If necessary, stop the capture first.

5.5 Save And Call Parameters

5.5.1 Saving the parameters: the user can save the current experiment Settings if necessary.

the < save > button under < parameter >, and the popup window will appear (see figure 5-11).

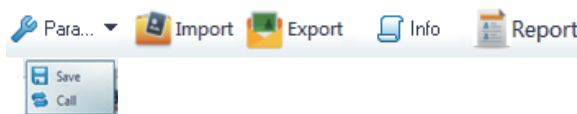

< add > : add new experimental parameters to the database

< delete > : delete unwanted experimental parameters

< save > : save the current experimental parameters to the database:

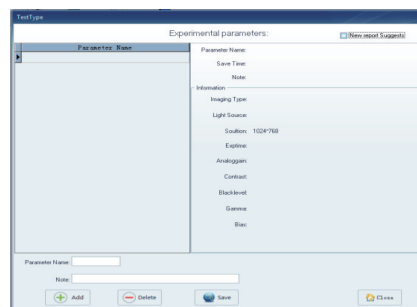


Fig 5-11



Whether to select when creating a new report

5.5.2 Calling the parameters: The user needs to call the experimental setting parameters before, click the < call > button under < parameter >, and the save window pops up (see figure 5-12).

< apply > : applies the selected experimental parameters to the current experiment.

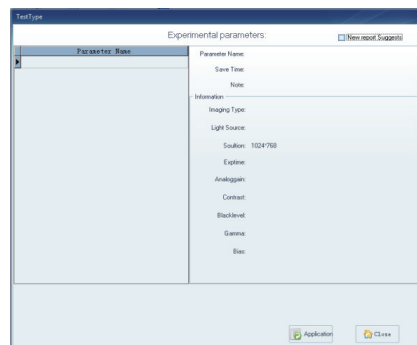


Fig 5-12

Chapter 6 Image Processing

This system provides more complete image processing functions, including image selection, copy and paste cutting, cropping, mirroring, adding text arrows, image scaling, rotation, color adjustment, brightness contrast adjustment, negative film effects, filtering, image merging and other functions.

In the shooting interface, click the <Image Processing> button on the toolbar to enter the image processing interface (see Figure 6-1).

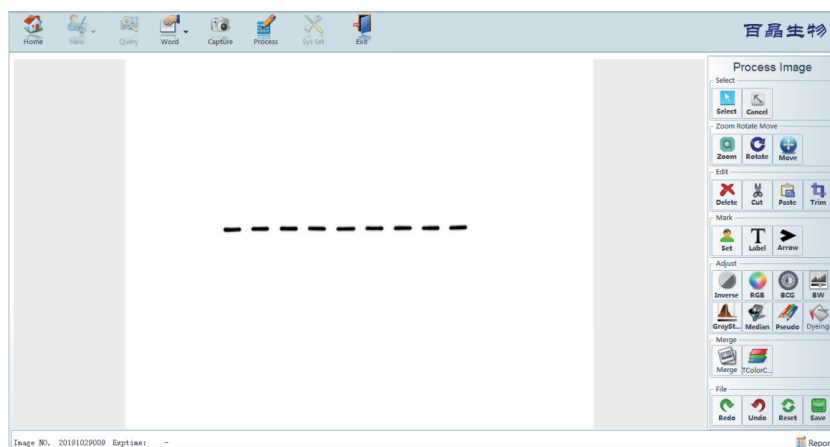


Fig 6-1

The user can realize each image processing function by clicking each tool button on the left side of the image processing window, and prompting the tool button function.

6.1 Selection Tool

Click the <Select> button, place the mouse in the selected state, and pop up the selection toolbar window (see Figure 6-2).

This software provides selection methods such as ellipse, rectangle, polygon, lasso, etc. After clicking the corresponding button, drag the image with the mouse to select the image.

Click the <Cancel> button and leave the mouse unselected.



Fig 6-2

6.2 Zoom, Rotate, Move Tools

6.2.1 Zoom:

Click the <Zoom> button to bring up the zoom toolbar window (see Figure 6-3).

Provides zoom mode: range zoom in, zoom in, zoom out, original image, fit window.

6.2.2 Rotate

Click the <Rotate> button to bring up the rotation toolbar window (see Figure 6-4).

Provide rotation mode: 90 degrees left rotation, 90 degrees right rotation, horizontal mirror, vertical mirror, rotate at any angle.

6.2.3 Move:

When the image is zoomed beyond the image window, you can click the <Move> button to drag the image.



Fig 6-3



Fig 6-4

6.3 Editing Tools

This software provides delete, cut, paste, and cut functions.

Methods as below.

6.3.1 Delete: Select the target and click the <Delete> button.

6.3.2 Cut: Select the area first, and then click the <Cut> button.

The selection one is cut and can be pasted elsewhere in the computer.

6.3.3 Paste: Select the area first, and then click the <Paste> button, which completes the copy and paste operations.

6.3.4 Crop: Select the area first, and then click the <Crop> button.

A window will pop up (see Figure 6-5). You can crop the image to the appropriate size as needed.

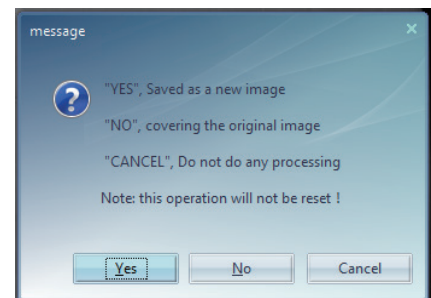


Fig 6-5

6.4 Mark Tools

This software provides text, arrow marks, and text and arrow settings. Methods as below.

6.4.1 Text: After clicking the <Text> button, drag the text window size in the image area with the mouse, and then double-click to enter text.

6.4.2 Arrow: After clicking the <Arrow> button, the mouse inserts an arrow in the image area.

6.4.3 Setting: After clicking the <Settings> button, the annotation setting window will pop up (see Figure 6-6).


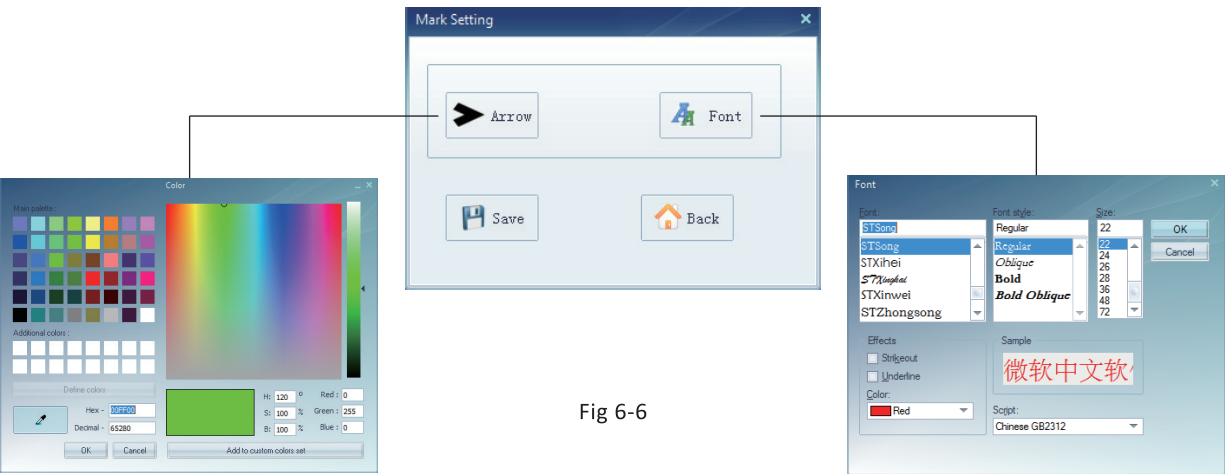


Fig 6-6

6.5 Adjustment tools

This software provides inverse color, RGB, BCG, black and white, averaging, filtering, false color, merge, and histogram adjustment.

6.5.1 Inverse Color: Click the <Invert Color> button to reverse the image.

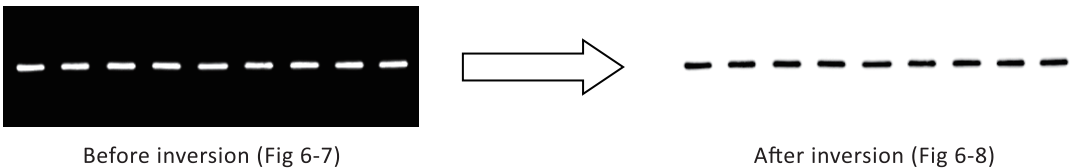

6.5.2 RGB adjustment:

Click the <RGB> button, a window will pop up (see Figure 6-9).

- a.RGB adjustment is to adjust the red, green, and blue component parameters of the image.
 - b.Method: keyboard input, up / down direction key adjustment, drag bar adjustment.
- Click the <Default> button to restore the default parameters.

6.5.3 BCG adjustment:

- Click the <BCG> button, a window will pop up, as shown in Figure 6-10.
- a.BCG adjustment is to adjust the brightness, ontrast, and gamma of the image
 - b.Method: keyboard input, up / down direction key adjustment, drag bar adjustment.
- Click the <Default> button to restore the default parameters.

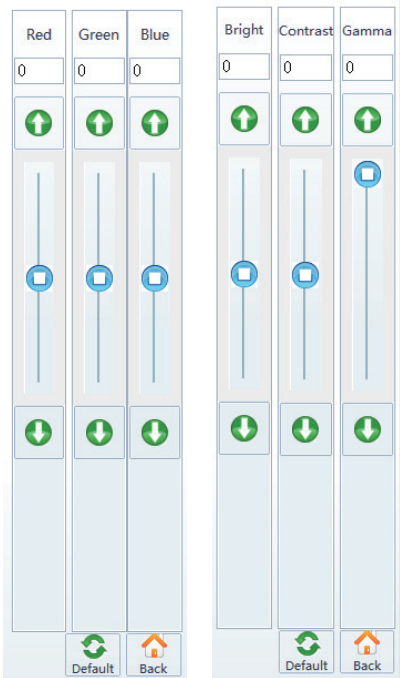



Fig 6-9

Fig 6-10

6.5.4 Black and white: The image can be converted from color mode to black and white mode.



Color image (Fig 6-11)

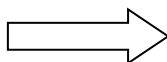

Black and white image (Fig 6-12)

6.5.5 Histogram:

Click the <Histogram> button to bring up the interface (see Figure 6-13).

Drag the scroll bar to adjust the histogram changes, and the image will follow the changes.



: Saves an image of the current histogram changes.



: Cancels the image with the current histogram change.


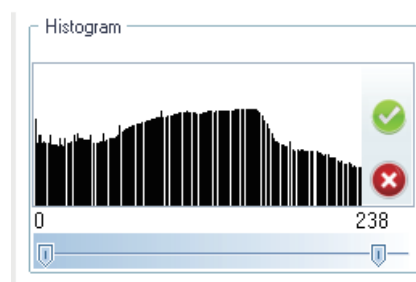


Fig 6-13

6.5.6 Filter: Click the <Filter> button to remove acnode noise and maintain the edge characteristics of the image.



Before filtering (Fig 6-14)

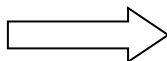


After filtering (Fig 6-15)

6.5.7 Pseudo-color: Click the <Pseudo-color> button to convert the grey level image to color image which help to find the difference among the sample bands.



Original image (Fig 6-16)

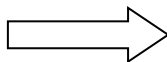

Pseudo-color image (Fig 6-17)

6.5.8 Dyeing: Click the <Dye > button to pop up the stain window (see Figure 6-18), you can stain the image.

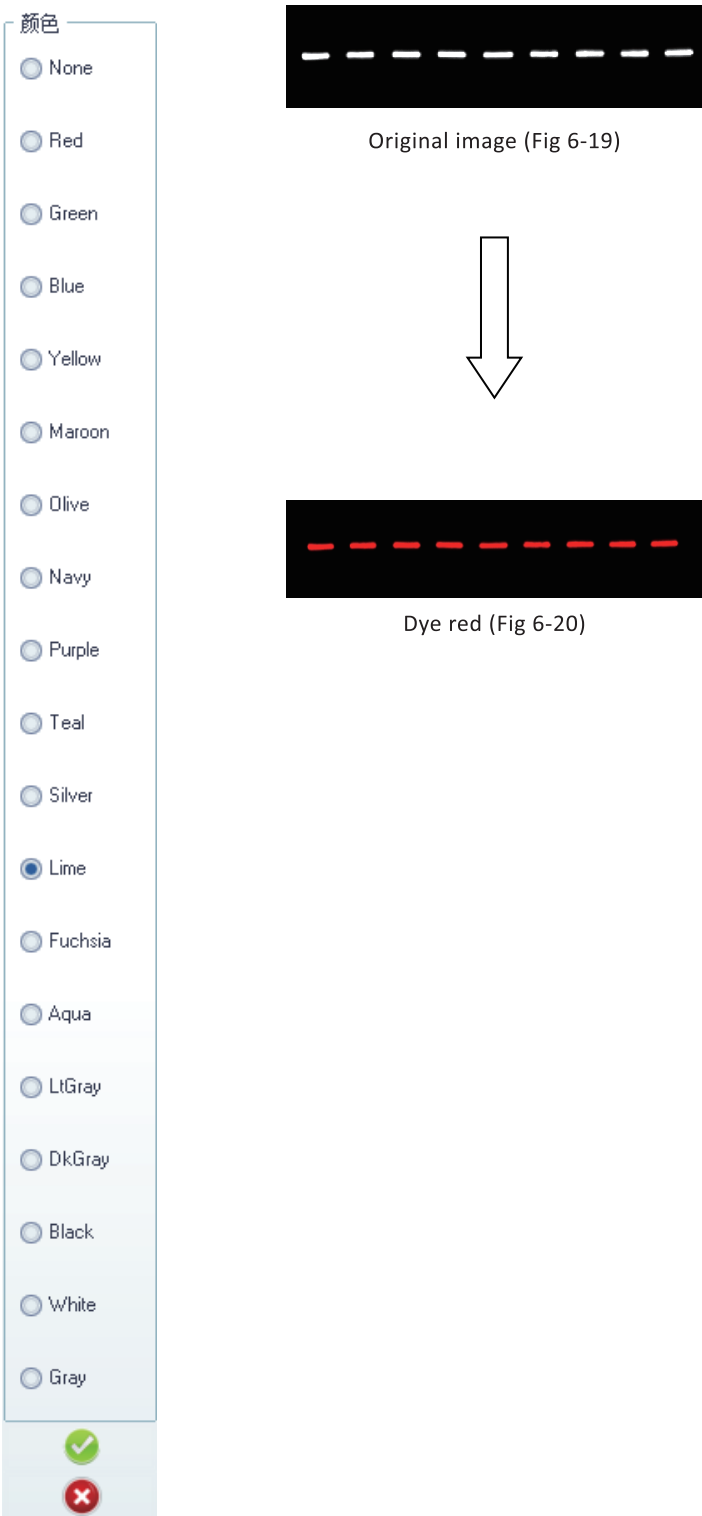


Fig 6-18

6.6 Image Merge Function

6.6.1Merge:
Click the <Merge> button to enter the image merge interface, as shown in Figure 6-21.
This function can merge the two images in different modes.The operation method is as follows:
Click the <Open> button to add the corresponding image.
Check the <Reverse> button to reverse the added image.
You can preview the effect by choosing different merging methods of Model.
Click the <Merge> button to merge the two images into one image.
And save it to the report.

6.6.2 Three-color merge:
Click the <Three Color Merge> button to enter the image merge interface. (Figure 6-22).
This function can merge the three images in different modes.
The operation method is as follows:
Check the color buttons (Red, Green, Blue) and click the <Open> button to add an image.
You can preview the effect by choosing different merging methods of Model.
Click the <Merge> button to merge the two images into one image and save it to the report.

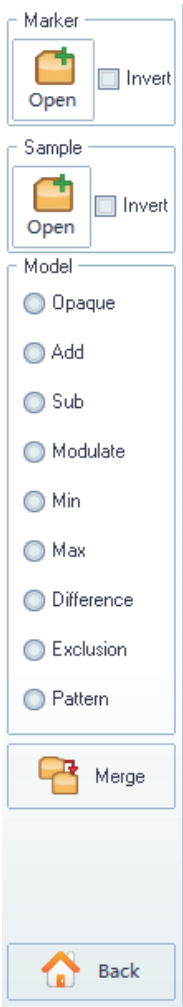


Fig 6-21



Fig 6-22

6.7 File Toolbar

This software provides functions for image operations: repeat, undo, reset, save, etc.

- 6.7.1 Repeat: Click the <Repeat> button to restore the undo operation.
- 6.7.2 Undo: Click the <Undo> button to undo the previous operation.
- 6.7.3 Reset: Click the <Reset> button to cancel all operations.
- 6.7.4 Save: Click the <Save> button, a message box will pop up (see Figure 6-23), and the user will follow the prompts as need.

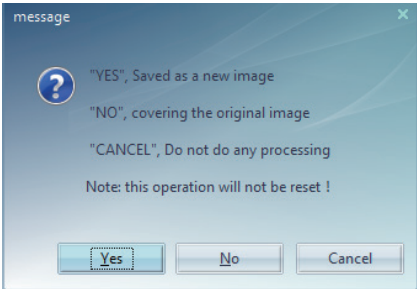


Fig 6-23

Chapter 7 System Settings

Click the <Sys Set> button to enter the system settings interface (see Figure 7-1).



Fig 7-1

7.1 User Management

The system provides flexible system user management functions. Different users can set different operation permissions.

This module can perform operations such as adding users, deleting users, and setting user permissions. After the user settings are completed, click the "Save Setting" button to save the current settings.

Click the <User Manage> button to enter the system user management interface (see Figure 7-2).

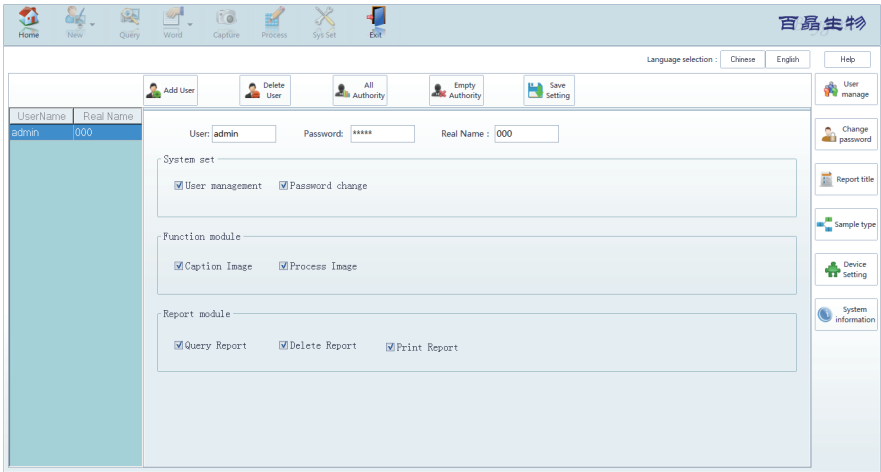


Fig 7-2

7.2 Change Password

Click the < Change Password > button to enter the password modification interface (see Figure 7-3), and the user can modify his own login password.

Original Password:

New Password:

Again Password:

Save

Fig 7-3

7.3 Report Title

Click the <Report Title> button to enter the report title modification interface (see Figure 7-4). The user can modify it or set the font format of the title.

Main Title:Beijing agricultural university

Sub Title:Molecular biology teaching and research section

Save

Fig 7-4

7.4 Sample Type

Click the <Sample Type> button to enter the sample type interface (see Figure 7-5). The user can add, edit, and modify the commonly used sample types for easy selection when entering report information.

DNA
RNA
westblot

Sample Type:

中文

Delete

Edit

Save

Fig 7-5

7.5 Device Settings

Click the <Device Settings> button to enter the device setting interface (see Figure 7-6). Set the communication port between the computer and the instrument. If the user does not know the specific serial port used by the computer, select <Automatically detect serial port at startup>.

The light source shuts down automatically

Regular time: 15 +

Communication serial port choice

Select COM 12 ☐ Automatic detection of serial port at startup

WesternBlot Marker

☒ Automatic shooting Marker-Image

Save

Fig 7-6

7.6 System Information

Click the <System Information> button to enter the system information interface and see software information (see Figure 7-7).

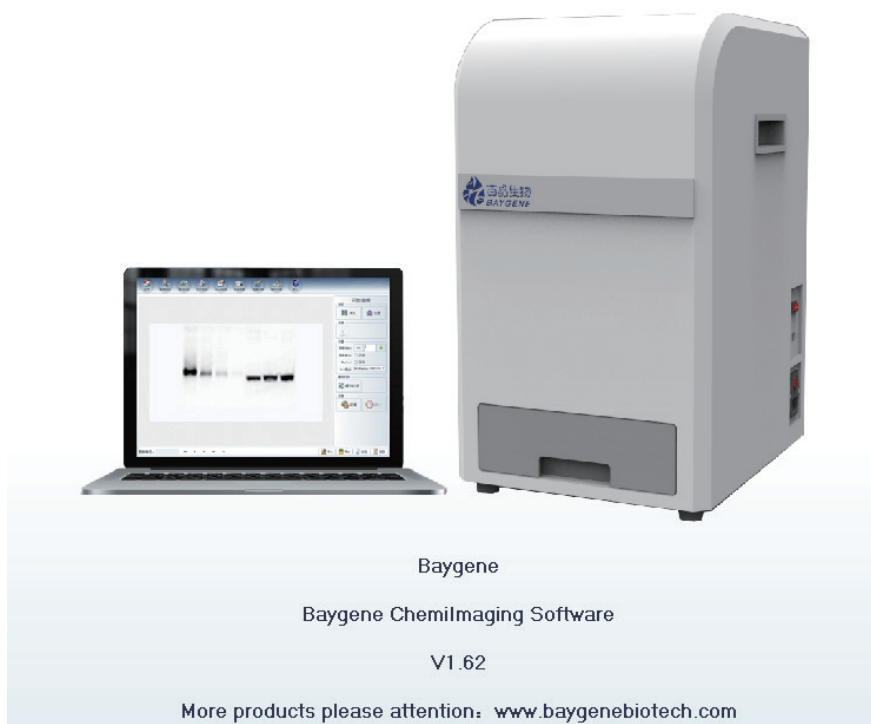


Fig 7-7

7.7 Other

7.7.1 Language selection: This software provides Chinese and English switching. Click the <Chinese> button or click the <English> button.

7.7.2 Help: Click the <Help> button, and the operation manual will pop up.

Chapter 8 Other Settings

In addition to the functions provided above, this software has image import and export functions.

This system supports common formats such as *.tif, *.jpg, *.bmp, *.png.

8.1 Image Import

Click the <Import> button at the bottom right of the window to pop up the window (see Figure 8-1).

Select the relevant image file and click the “Open” button to transfer the image into this system for image processing.

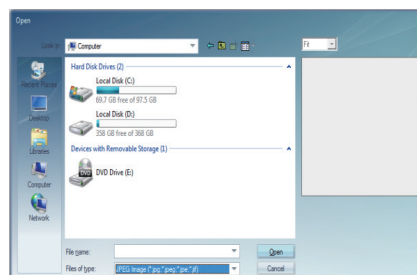


Fig 8-1

8.2 Image Export

Image output is to save the current image to the specified file path.

Click the <Export> button, a window will pop up, as shown in Figure 8-2.

Select the file format, enter the file name to be saved, and click the <Save> button to export and save the image.

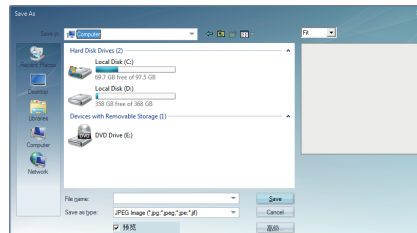


Fig 8-2

8.3 System Backup And Recovery

Copy the installation root folder completely to the backup address. (Figure 8-3).

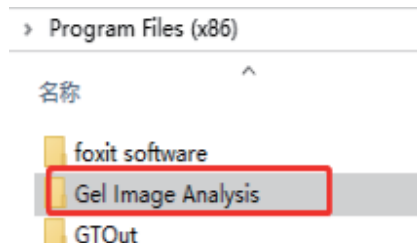


Fig 8-3

Chapter 9 Precautions for Use

Daily maintenance and care

1. The built-in camera of the instrument is fragile and valuable (not guaranteed), please protect it carefully.
2. Never disassemble the instrument at will, or you will not be responsible for the warranty.
3. After use, turn off the power and clean the sample stage.
4. The instrument should be placed in a well-ventilated, dry and dust-free place. When not in use, turn off the main power switch on the back of the cabinet.
5. The instrument should be placed in a packing box during packaging, storage and transportation. The upper and lower sides are lined with shock-proof materials. The outer packaging shall be fastened with an external cable tie and shall have the required marks. Each instrument must have a product certificate and instructions for use. Carton boxes should be handled carefully during transportation, and the storage warehouse should be well ventilated and free of corrosive gases.

Chapter 10 Frequent Q & A for Use

Common Error	Solution
Q1、Start the software and a pop-up window of "Camera not connected" will pop up.	A1、 (1)Please confirm whether the instrument is connected to the direct data line of the computer correctly; (2)Check whether the camera driver is installed in the device manager;
Q2、Start the software, and a prompt window "Communication Connection Error" pops up.	A2、 (1) Detect the communication port number (2) Check whether the device serial port driver is installed; (3) Please confirm the data connection between the instrument and the computer.
Q3、Click the LED light button, the light is off.	A3、 (1)Same as A2; (2)Check whether the power switch is on.
Q4、Image is not clear	A4、 Please contact the after-sales engineer to check if the lens focus is clear.
Q5、Other problems.	A5、 Please contact the after-sales engineer.

Chapter 11 After Sales

11.1 Product Warranty Terms

In order to enable consumers to use our products with confidence and satisfaction, our company will provide after-sales service in strict accordance with the relevant laws and regulations issued by the state.

1. Consumers should fill out the warranty card and stamp the dealer's seal when purchasing our products.
2. Service period: If our product fails during normal use within 7 days from the date of sale, consumers can choose services such as refund, replacement, and warranty. Products are warranted for one full year from the date of purchase which products are not being damaged by any man-made factors. For consumers who do not meet the criteria of free replacement or warranty service, our company still provides technical services, and only charges the material cost when need repair parts.
3. The purchase time is based on the date of the invoice or receipt issued by the dealer.
4. The "Three Guarantees" service cannot be provided under any of the following conditions:
 - 1) All damage caused by human factors and use in abnormal working s, failures and damage caused by using the product or using in the environment not following the user manual;
 - 2) The user disassemble, repairs, or modify the product without permission of the company;
 - 3) Damage caused by bad transportation after purchasing our products;
 - 4) Damage caused by other force majeure (such as floods, lightning, earthquakes, abnormal voltages);
 - 5) Instrument aging problems such as worn, abrasion, cracking and dipping;
 - 6) Products that do not belong to the company (such as fakes);
 - 7) Cannot present a valid shopping voucher, no warranty card, etc .;
 - 8) The serial No. (barcode) of the instrument is damaged.

11.2 After-sales service procedures

If there is a problem when use our product, bring the invoices or related receipts to the dealer for repair. If the problem is not properly resolved, please call or write to our customer service department, and we will help you solve the problem in the shortest time.

Product Warranty Card (Customer Page)

Product information	Customer information
Product name&model:	Institutions:
Serial number:	Contacts:
Installation date:	Address:
Installation engineer:	Tel:
Invoice number:	E-mail:

(Customer page is kept by the client)

Address: Building7, No.28 Yuhua Road, Airport Industrial Zone,
Area B, Beijing 101300, PR China.
Tel: 010-80483100 80483200 Fax: 010-80482859
web:www.baygenebiotech.com E-mail:info@baygenebiotech.com

Product Warranty Card (Corporation Page)

Product information	Customer information
Product name&model:	Institutions:
Serial number:	Contacts:
Installation date:	Address:
Installation engineer:	Tel:
Invoice number:	E-mail:

(Corporation page will be filled out by the client and handed over to the company)

Address: Building7, No.28 Yuhua Road, Airport Industrial Zone,
Area B, Beijing 101300, PR China.
Tel: 010-80483100 80483200 Fax: 010-80482859
web:www.baygenebiotech.com E-mail:info@baygenebiotech.com

BEIJING BAYGENE BIOTECHNOLOGIES CO., LTD.

Address: No. 28, Yuhua Road, Shunyi District, Beijing (No. 7-8, Airport Area B)

Postcode: 101300

Tel: 010-80483100/80483200 010-80483456/80483457

Fax: 010-80482859

Website: www.baygenebiotech.com

Email: info@baygenebiotech.com