

Fast and Gentle Multiplex Imaging for your LSM 800 Airyscan



Multiplex Imaging Modes available for LSM 800 Airyscan

Get better data faster.


Introduced with Airyscan 2 for LSM 900, Multiplex Imaging Modes opened a new dimension in acquisition performance. Multiplex Imaging Modes can now enhance the performance of your LSM 800 without purchasing a new system. The performance boost is easily achieved with an onsite upgrade.

With Multiplex Imaging Modes you can combine

- superresolution imaging
- high SNR
- larger field-of-views

leading to fast results.

Multiplex modes use knowledge of the shape of the excitation laser spot and the location of the individual detector array elements to extract more spatial information, even during the parallel pixel readout. This allows for larger steps when sweeping the excitation laser over the field of view, improving acquisition speed. Capturing more spatial information in the pinhole plane across the detector array enables final image reconstruction with better resolution than the acquisition sampling.

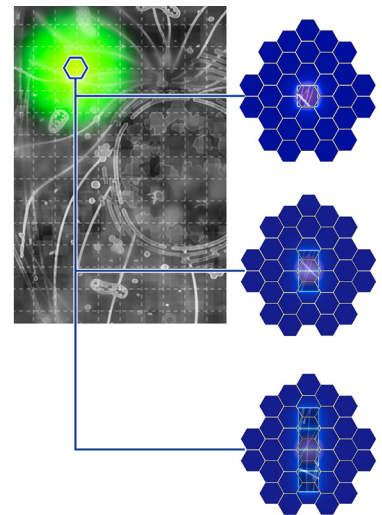

Superresolution

- Higher acquisition speed
 - 4 fps @ Airyscan SR
 - 8.4 fps @ MPLX SR-2Y
 - 18.9 fps @ MPLX SR-4Y
- High SNR
- Rapid tiling of large areas
- Efficient live cell imaging
- Fast volumetric imaging



Continuity

- Apply your existing LSM 800 expertise and proven workflows to new applications
- Better software experience
- ZEISS Service and Support



Sustainability

- Continue using your LSM 800 Airyscan and achieve a whole new performance level



Seeing beyond

Multiplex Imaging Modes on LSM 800 Airyscan

Fast and gentle Multiplex imaging

ZEISS LSM 800 with Airyscan after MPLX upgrade				
	Airyscan SR	Multiplex SR-2Y	Multiplex C0-2Y	Multiplex SR-4Y
ZEN Module	Airyscan	Airyscan 2 MPLX basic	Airyscan 2 MPLX basic	Airyscan 2 MPLX standard
Paralellization	1	2	2	4
Resolution x/y [nm]	120/120	140/140	180/180	140/140
Resolution z [nm]	350	450	550	450
FPS at 512x512 pixels	4	8.4	8.3	18.9
FPS at max FOV (Zoom 1.3)	0.4	0.8	3.5	3.5
Antibody labeling, fine structures	+++++	++++	++	++++
Antibody labeling, tiling	++	+++	+++	+++++
Live cell imaging	++	+++	+++++	++++
with C-Apo 63x/1.4 Oil @ 488				

What's in the box?

- Upgrade to ZEN 3.7
- ZEISS workstation*
- ZEN Module Airyscan 2 MPLX basic
- ZEN Module Airyscan 2 MPLX standard*
- Airyscan 2 detector*

*optional

What is required?

- LSM 800
- Airyscan detector
- ZEISS workstation
 - hp Z2G4 or
 - hp Z2G5 or
 - hp Z6G4
 - min. 64 GB RAM

What options are available?

- Airyscan 2 detector
- GaAsP detectors
- Sample Navigator
- ZEN Module LSM Plus
- ZEN Module Airyscan jDCV
- Autoimmersion Module**
- AI Sample Finder**
- Predictive Service
- ... and more

**only available for Axio Observer 7

Upgrade your ZEISS microscope

Upgrade your ZEISS microscope and benefit from a performance boost and compatibility with other ZEISS imaging systems.

Please contact us to learn more about the suitability of your LSM 800 for this upgrade opportunity and how your research processes will benefit from an upgrade:

microscopy@zeiss.com



microscopy@zeiss.com
www.zeiss.com/microscopy