

Nano ITC, Affinity ITC, Affinity ITC Auto



Site Preparation Guide

Table of Contents

Table of Contents	2
Ideal Setup	3
System Components.....	4
Instrument Measurements	5-6
Nano ITC, Affinity ITC	5
Affinity ITC Auto, Autosampler	6
Utility Requirements.....	7
Computer Requirements	8-9
Hardware.....	8
Software	9
Site Preparation Checklist	10
TA Instrument Offices.....	11

Ideal Setup



IDEAL PLACEMENT AND BENCH MEASUREMENTS

Select a location with adequate floor space and a rigid laboratory bench that is level and is in a vibration-free environment.



Bench length: 2.4 m (8 ft)

Bench depth: 76 cm (30 in) min.

Distance from the wall: 20 cm (8 in) min.



An **Autosampler** requires an additional 30 cm (1 ft) bench space to accommodate up to 6 bottles (0.5 L-4 L) to the left of the Autosampler and 1-2 waste outlet collection bottles (1 L-4 L) at the back. Additional floor space beneath the Autosampler is required for 4 L spent rinse fluid collection bottle.

System Components



MAIN SYSTEM COMPONENTS



- A. Computer (Controller)
- B. Autosampler (Optional)
- C. Instrument

Instrument Measurements



NANO ITC



Height: 28 cm (11 in)

Width: 35 cm (14 in)

Depth: 53 cm (21 in)

Weight: 17 kg (37 lbs)



AFFINITY ITC

Height with arm LOWERED: 53 cm (21 in)

Height with arm RAISED: 71 cm (28 in)

Width: 35 cm (14 in)

Depth: 53 cm (21 in)

Weight: 17 kg (38 lbs)



Instrument Measurements



AFFINITY ITC AUTO – INSTRUMENT



Height with arm LOWERED: 53 cm (21 in)

Height with arm RAISED: 71 cm (28 in)

Width: 35 cm (14 in)

Depth: 53 cm (21 in)

Weight: 20 kg (44 lbs)




AFFINITY ITC AUTO – AUTOSAMPLER

Height: 35 cm (14 in)

Width: 30 cm (12 in)

Depth: 56 cm (22 in)

Weight: 21 kg (46 lbs)

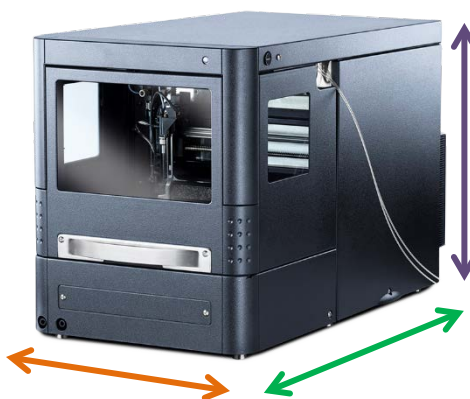

Utility Requirements



LABORATORY

Item	Requirement
Temperature	15–30°C
Relative Humidity	5–80% (non-condensing)
Temperature Stability	±1°C with changes in temperature being gradual shifts instead of fast changes
Instrument Location Environment	• Dust-free • Vibration-free • Away from exposure to direct sunlight and direct air drafts • Pollution Degree 2 Environment • Maximum altitude: 2000 m (6560 ft)



POWER

Item	Requirement
Instrument power	100–240 VAC, 3A, 50/60 Hz
Power line	Grounded, single-phase line for instrument and computer, not shared with motors, heaters, or compressors • 15 A for voltages near 120 VAC • 10 A for voltages near 230 VAC
Electrical power cord	• The plug of the cord must be rated to carry at least 125% of the product current rating. • The cord length must be less than 4.5 meters and must be UL or CSA approved.
Customer-supplied	Surge suppressor power strip



Computer Requirements



HARDWARE REQUIREMENTS

Item	Requirement
Processor	• Intel® Core™ i3 or better • 2.93 GHz with 3 MB L2 cache
Memory	Required: ≥ 8 GB RAM (if the computer is only acting as the controller) Recommended: ≥ 16 GB RAM or greater (for best analysis performance)
Hard drive	≥ 300 GB free space
DVD	≥ 16 x CD-ROM optional. Obtain the latest software updates at http://www.tainstruments.com/support/software-downloads-support/downloads
Screen resolution	Minimum: 1280 x 1024 with 24-bit colors Recommended: 1920 x 1080 with 24-bit colors
Graphic card	128 MB DirectX10 or higher recommended
Screen (LCD) size	Recommended: 24" wide screen
USB Ports (1 for instrument and 1 for Autosampler)	2.0

Computer Requirements



SOFTWARE REQUIREMENTS

Item	Requirement
Operating System	- Windows 7, 8, 10 Ultimate, Enterprise & Professional - Home version not supported
	Required: 32-bit or 64-bit version Recommended: 64-bit version
Browser	Internet Explorer
Other	Microsoft Operating System Service Pack
Other	Windows Operating System and associated Microsoft updates must be up-to-date
Network	<i>TA Instruments is not responsible for resolving issues associated with connections to your corporate network.</i>
Conflicts	<i>TA Instruments is not responsible for resolving hardware/software conflicts created by the addition of third party hardware or software to the computer.</i>

Site Preparation Checklist



Nano ITC, Affinity ITC, Affinity ITC Auto



Sufficient bench space for instrument, computer, sample prep space and tools, and Autosampler (if applicable):

- ☐ Length: 2.4 m (8 ft) for instrument, computer, and sample prep and tools. Add an additional 30 cm (1 ft) for Autosampler.
- ☐ Depth: 76 cm (30 in) minimum.
- ☐ Distance from the wall: 20 cm (8 in) minimum

Laboratory conditions meet the following requirements:

- ☐ Dust-free, vibration-free, away from exposure to direct sunlight/air drafts, in a pollution degree 2 environment
- ☐ Temperature is 15–30°C with a stability of $\pm 1^\circ\text{C}$.
- ☐ Relative humidity is 5–80% non-condensing
- ☐ Maximum altitude is 2000 m (6560 ft)



- ☐ Instrument power: 100–240 VAC, 3A, 50/60 Hz

I hereby acknowledge that all utility requirements have been met per the checklist above and that they will be ready at the agreed time of installation.

If all utility requirements are not met at the agreed time of installation, additional charges may be incurred for a return Service trip.

Customer DD / MM / YYYY

Company City State Country

Please send a signed copy of the completed checklist to your local Service representative.

TA Instruments Offices

For information on our latest products, contact information, and more, see our website at:
<http://www.tainstruments.com>.

To find your local TA Instruments office and contact information, visit
<http://www.tainstruments.com/contact/ta-directory/>

TA Instruments – Waters LLC
Corporate Headquarters
159 Lukens Drive
New Castle, DE 19720
USA

Telephone: 302-427-4000
Fax: 302-427-4001
Email: info@tainstruments.com