

TOUCH IMAGER

Simple Operation Instruction

I . Powering On the Instrument

1. Turn on the power switch located on the right side of the instrument.
2. Wait for the indicator light to change from red to deep green. Then, double-click the **"Touch Viewer"** software icon on the desktop. You may log in with your personal account, or simply use the default username "Default" without a password.
3. After logging in, please wait for the instrument to complete the automatic calibration before proceeding with any operations.

Note: During the automatic connection and calibration process, **please keep the top lid closed** and **ensure that no objects are present inside the device**.

II . Image Acquisition

Chemiluminescence Imaging: Automatic Continuous Exposure and Manual Exposure.

1. Automatic Mode:

Click **"Auto Capture"** to begin. The software will automatically capture 8 images over 30 seconds, at the following exposure times in sequence: 1s, 2s, 3s, 5s, 8s, 10s, 20s, 30s. If an overexposure warning (highlighted in pink) appears during the process, you may click **"Stop"** at any time to halt the capture.


2. Manual Mode:

a. Sample Capture: Click **"Capture"** to take a single image with the default exposure time of 1 second.


You may also manually customize the settings by selecting or entering the initial exposure time, interval time, and number of images to be captured. Then click **"Capture"** to start a customized continuous exposure sequence.

Note:

1. To set an exposure time shorter than 1 second, use the dropdown menu and select the desired exposure time.

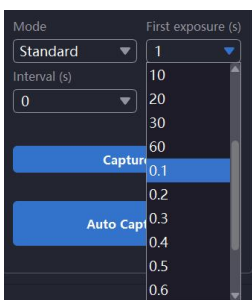

2. If the signal is too strong and overexposure still occurs even at 0.1s, please switch from Standard Mode to High Full Well Mode.

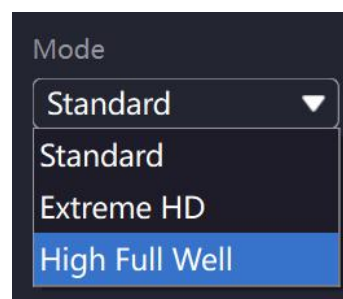

b. Marker Capture: Click the “**Marker**” button to capture the marker image.

c. Merge Images: Select one **Sample** image and one **Marker** image, then click “**Merge**” at the bottom right of the software interface to generate the merged image.

Note: The brightness and contrast of all images can be manually adjusted at the bottom of the acquisition interface.



d. Image Export: Click the checkbox at the top right corner of each image thumbnail you want, then click “**Export**” to save the selected images.

III. Precautions

1. Use the designated anti-static, silicone-tipped tweezers only. **Do not** use metal tweezers or other sharp plastic tweezers.
2. During operation, clean the imaging chip surface with the provided soft cloth only. **Do not** use rough paper or rigid cloths.
3. Avoid direct exposure of the imaging chip surface to sunlight or other strong light sources.
4. This instrument is designed for immediate use without pre-cooling. Turn off the device immediately after use.
5. If the device will not be used for an extended period, please disconnect the power supply.

Technical Support Email: **support@e-blot.com**

Official Website: **www.e-blot.com**