


# QIAgility™ — Pure Precision




Sample & Assay Technologies

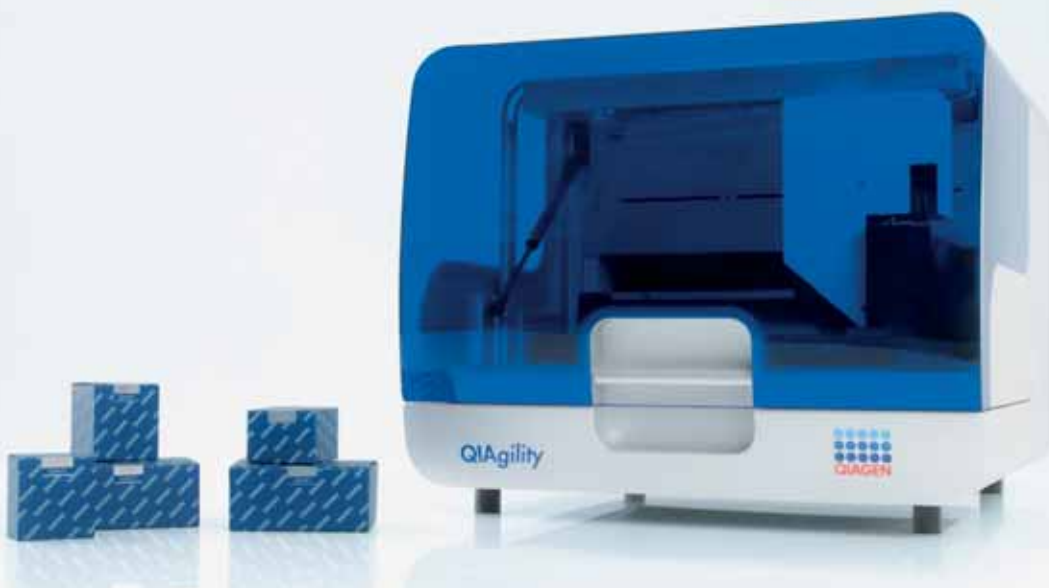

## Fully automated sample and assay technologies

The QIAgility is a compact benchtop instrument that enables automated PCR setup in a wide range of formats, including PCR assays for the Rotor-Gene™ Q. The high precision of the QIAgility delivers the reproducible results you need in your end-point and real-time PCR assays, from sample to sample and laboratory to laboratory. Automated PCR setup is rapid and reliable and eliminates manual pipetting steps that can be prone to human error. The unmatched versatility of the QIAgility means that almost all tube and plate formats are supported as well as Rotor-Discs for the Rotor-Gene Q.

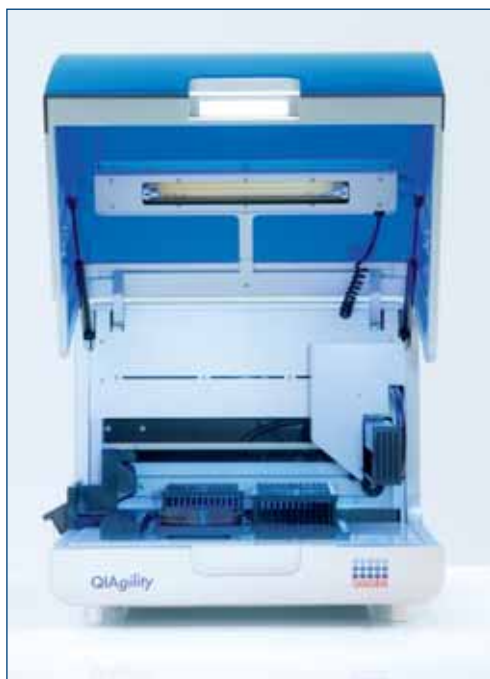

The QIAgility provides:

- Automated PCR setup in all formats
- Convenient, easy-to-use software
- Standardized results and increased productivity
- Elimination of manual pipetting steps
- Seamless integration with QIAGEN® sample technologies


The QIAgility is operated via a computer with easy-to-use software. Simply select a program, set up the worktable, and walk away during PCR setup. Samples are ready for immediate use in downstream applications.

## Advanced technology for pure precision

Dynamic liquid-level sensing enables the QIAgility to deliver high-precision, standardized pipetting (Table 1 and Figure 1). The instrument senses the meniscus during each pipetting event and retracts the tip while dispensing, enhancing accuracy and minimizing liquid carry-over. If the liquid is viscous, the QIAgility can adapt its pipetting speed accordingly.

## Exceptional safety features

The advanced design of the QIAgility delivers improved process safety, giving you peace of mind. An in-built sensor ensures that the instrument pauses when the hood is opened during a run, protecting users from moving parts. An optional UV light ensures effective worktable decontamination and an optional HEPA filter provides positive clean air pressure throughout PCR setup. These safety features are complemented by intelligent safety warnings from the instrument software, should an error occur.


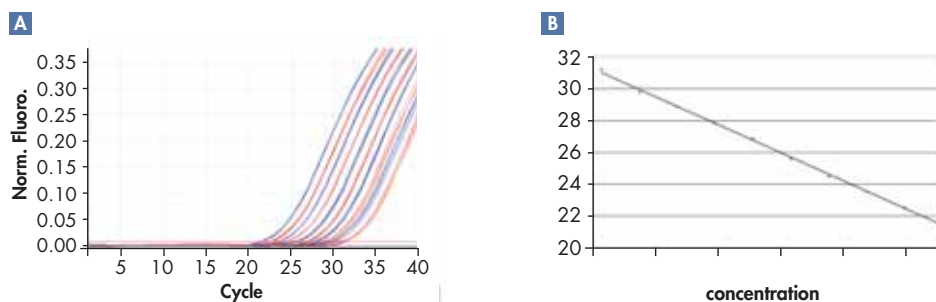


**Figure 1. Exceptional reproducibility and accurate series dilutions.** A 2-fold series dilution of human genomic DNA (30 ng to 0.6 ng) was performed using the QIAgility. IL1R2 was amplified from these dilutions, using the Rotor-Gene Q and the QuantiFast® Probe PCR Kit. Reactions were set up in replicates of 5. DNA yields were quantified by quantitative, real-time PCR, where the Rotor-Gene Q was also used for analysis. Plots show **A** highly precise and reproducible replicates with a 1 cycle difference between each 2-fold dilution and **B** the standard curve.

**Table 1. Unrivalled precision**

| Volume (µl) | Mean (mg) | S.D. (mg) | CV (%) |
|-------------|-----------|-----------|--------|
| 1           | 0.99      | 0.06      | 5.64   |
| 2           | 2.03      | 0.05      | 2.20   |
| 5           | 5.02      | 0.03      | 0.56   |
| 10          | 9.99      | 0.02      | 0.24   |
| 20          | 19.97     | 0.04      | 0.22   |
| 50          | 49.99     | 0.05      | 0.10   |
| 100         | 100.84    | 0.14      | 0.14   |
| 150         | 150.79    | 0.25      | 0.17   |
| 200         | 198.82    | 0.26      | 0.13   |


The QIAgility was used to pipet 1, 2, 5, 10, 20, 50, 100, 150, and 200 µl water, in replicates of 16. The actual volume pipetted was determined for all aliquots based on weight (i.e., 1 g=1000 µl). The mean, standard deviation (S.D.), and coefficient of variation (CV) was calculated for each volume.

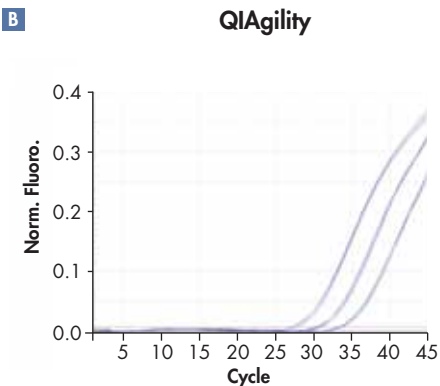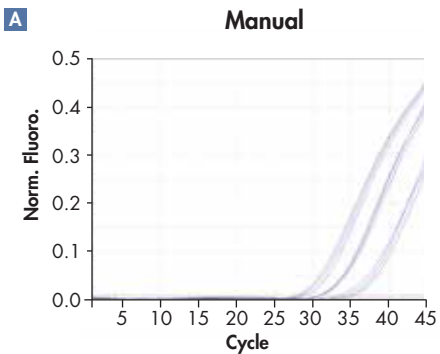

**Figure 3. Automated reaction setup is more accurate than manual reaction setup.** IL12RB1 was amplified from leukocyte RNA (1, 10, and 100 ng), using the Rotor-Gene Q and the QuantiFast Probe RT-PCR Kit. Replicates of 3 were set up **A** manually and **B** using the QIAgility. DNA yields were quantified by quantitative, real-time RT-PCR, where the Rotor-Gene Q was also used for analysis.

### Absolute reliability and reproducibility

The QIAgility is provided with optimized protocols for unmatched reliability and reproducibility. In the future, more protocols will become available. Manual pipetting steps, which are typically prone to human error, are eliminated (Figure 2 and Figure 3). Using the software in “Virtual Mode” enables user-developed protocols to be tested in the absence of reagents, avoiding costly repetition of experiments. To reduce your work load, there is no need for time-consuming manual documentation of runs since comprehensive run reports, containing all run information, are generated before and after each protocol.

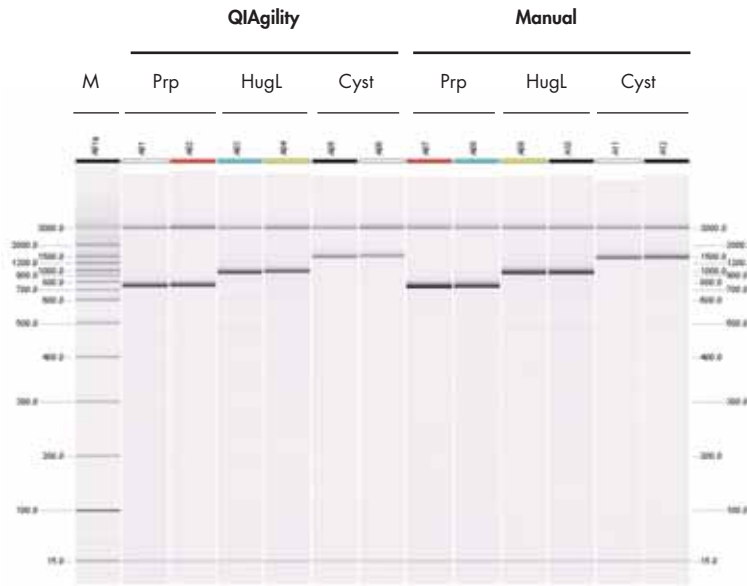

**Figure 2. Eliminate manual pipetting steps.** Prp, HugL and Cyst were amplified from human genomic DNA (10 ng). PCR reactions were set up using the TopTaq™ Master Mix Kit, in replicates of 2, using the QIAgility and manually. After cycling, reactions were analyzed by capillary electrophoresis, using the QIAxcel™.

# QIAgility

## Flexible workflow with built-in versatility

The QIAgility is compatible with all your workflow needs. The worktable is easy to set up and, for optimal versatility, is configured to support almost all PCR tube types and plate formats, including 96-well and 384-well plates as well as Rotor-Discs for the Rotor-Gene Q. The QIAgility can set up multiple master mixes and can process multiple PCR setups in parallel.



In addition to PCR setup, the flexibility of the QIAgility allows the user to program an unrivaled range of high-precision pipetting applications, which includes:

- Normalization of DNA and RNA concentration
- Transfer of liquid samples from one tube format to another
- Serial dilutions with variable dilution ratios
- Restriction digest setup
- Sample pooling
- Selective pipetting from archived sample banks





## Complete your fully automated workflow

The QIAgility fully automates proven QIAGEN PCR kits (Table 2). Visit [www.qiagen.com](http://www.qiagen.com) for up-to-date product information. For increased convenience, the QIAgility provides a broad range of output formats, enabling compatibility with all real-time cyclers, including the Rotor-Gene Q. When used in combination with the EZ1® Advanced automated solutions, the QIAgility can complete and enhance your automated, real-time PCR workflow (see next page).

**Table 2. Use the QIAgility to fully automate proven QIAGEN PCR kits.**

| Application                               | Method                   | Kit                                         | Description                                                                            |
|-------------------------------------------|--------------------------|---------------------------------------------|----------------------------------------------------------------------------------------|
| Gene expression analysis                  | Real-time RT-PCR         | Rotor-Gene Kits                             | Rapid and reliable RT-PCR on the Rotor-Gene Q.                                         |
|                                           | Real-time RT-PCR         | QuantiFast Kits                             | Fast, sensitive RT-PCR on any real-time cycler.                                        |
|                                           | Real-time RT-PCR         | QuantiTect® Kits                            | Sensitive RT-PCR on any real-time cycler.                                              |
|                                           | End-point RT-PCR         | QIAGEN OneStep RT-PCR Kit                   | Sensitive RT-PCR, with simple setup, using any RNA template.                           |
| Virus detection                           | Real-time PCR and RT-PCR | QuantiTect Kits                             | Highly sensitive detection of multiple viral DNA and RNA targets in the same tube.     |
|                                           | End-point RT-PCR         | QIAGEN OneStep RT-PCR Kit                   | Highly sensitive detection of viral RNA.                                               |
| Genotyping                                | End-point PCR            | Type-it™ PCR Kits                           | Accurate and reliable SNP genotyping, mutation detection, and microsatellite analysis. |
| Cloning                                   | End-point PCR and RT-PCR | QIAGEN LongRange Kits                       | Sensitive and accurate amplification of extremely long PCR products.                   |
|                                           | End-point PCR            | HotStar HiFidelity Polymerase Kit           | High-fidelity, hot-start PCR.                                                          |
| Pathogen detection                        | Real-time PCR            | <i>cador</i> ™ veterinary PCR Kits          | Ready-to-use, validated real-time PCR assays to detect veterinary pathogens.           |
| Standard and specialized PCR applications | End-point PCR            | HotStarTaq® and HotStarTaq <i>Plus</i> Kits | Highly specific hot-start PCR without optimization.                                    |
|                                           | End-point PCR            | QIAGEN Taq Kits                             | Highly reliable PCR at first attempt.                                                  |
|                                           | End-point PCR            | QIAGEN Fast Cycling PCR Kit                 | Rapid and specific PCR on any thermal cycler.                                          |
|                                           | End-point PCR            | TopTaq Kits                                 | Unmatched reliability with unrivaled ease-of-use.                                      |



# Automated solutions from sample to result — free up your time!

## Purification

**EZ1 Advanced  
EZ1 Advanced XL**  
Nucleic acid purification  
solution using magnetic-  
particle technology


## Assay setup

**QIAgility**  
PCR setup system

## Detection and analysis

**Rotor-Gene Q**  
Real-time PCR cyclers

**Real-time PCR**

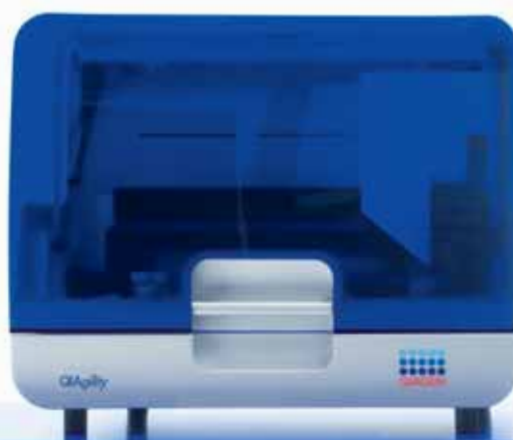



## Efficient, walkaway automation

The QIAgility enables rapid and efficient processing of complex reaction setups, freeing up your valuable laboratory time (Table 3). The single-channel pipet can process up to 96 samples in just 17 minutes, providing enormous time savings and eliminating tedious manual pipetting which can cause repetitive strain injury. While the QIAgility is performing reaction setup, the user can concentrate on more important laboratory tasks.

**Table 3. Speed and versatility.**

| Application              | Description                                  | Number of reactions | Tip reuse?    | Tip consumption                                | Run time (minutes) |
|--------------------------|----------------------------------------------|---------------------|---------------|------------------------------------------------|--------------------|
| Normalization            | One-step dilution                            | 96                  | No            | 108 x T <sup>50</sup> *                        | 28                 |
| Reformatting             | 1:1 transfer                                 | 96                  | No            | 96 x T <sup>50</sup>                           | 17                 |
| Plate replication        | 1:4 transfer                                 | 96–384              | No            | 154 x T <sup>50</sup>                          | 71                 |
| Plate replication        | 1:4 transfer                                 | 96–384              | Yes           | 96 x T <sup>50</sup>                           | 46                 |
| Serial dilution          | 16 concentrations                            | 16 dilutions        | No            | 32 x T <sup>50</sup>                           | 7                  |
| Repetitive aliquoting    | Primer solution (2 ml) into aliquots (20 µl) | 96 aliquots         | Yes (8 times) | 12 x T <sup>50</sup>                           | 10                 |
| Repetitive aliquoting    | Primer solution (2 ml) into aliquots (20 µl) | 96 aliquots         | No            | 96 x T <sup>50</sup>                           | 17                 |
| Preparation of standards | 6 concentrations                             | 6 dilutions         | No            | 11 x T <sup>50</sup><br>5 x T <sup>200</sup> † | 3                  |
| Master mix preparation   | 8 components for 96 samples                  | 96                  | Yes (8 times) | 10 x T <sup>200</sup>                          | 6                  |
| PCR setup                | Premixed master mix                          | 96                  | Yes (8 times) | 108 x T <sup>50</sup>                          | 36                 |
| PCR setup                | Master mix preparation                       | 96                  | Yes (8 times) | 110 x T <sup>50</sup><br>10 x T <sup>200</sup> | 39                 |
| PCR setup                | Premixed master mix and prediluted           | 96                  | Yes (8 times) | 108 x T <sup>50</sup>                          | 36                 |
| PCR setup                | Master mix preparation and standard dilution | 96                  | Yes (8 times) | 115 x T <sup>50</sup><br>15 x T <sup>200</sup> | 41                 |

\* T<sup>50</sup>=50 µl Tips. † T<sup>200</sup>=200 µl Tips.

The QIAgility can be used in a wide range of liquid-handling operations. The user can control tip usage, enabling cost-efficient use of consumables. Reactions can be set up using master mix that has already been prepared, or the QIAgility can prepare the master mix.





## Intuitive software for ease of use


With the QIAgility, specialist training is not required. QIAgility software simulates the instrument worktable, making it exceptionally easy to use by the automation novice or more experienced laboratory personnel. Fast setup and development of user protocols are facilitated by wizards, which guide the user through software procedures in a stepwise fashion. The QIAgility is preinstalled with optimized protocols for convenient and effective PCR setup. The software automatically calculates reagent volumes for the selected number of samples, giving you complete cost control. In addition, the user can develop custom protocols tailored to specific research requirements.

## Comprehensive service and support

QIAGEN Instrument Service provides flexible service support agreements tailored to your needs. We also provide Warranty PLUS extended warranties, giving you complete cost control and enhanced coverage through priority response times.

### Specifications of the QIAgility

|                         |                                                                                                           |
|-------------------------|-----------------------------------------------------------------------------------------------------------|
| Dimensions              | Width: 55 cm (21.6 in.)<br>Depth: 63 cm (24.8 in.) without cables<br>Height: 45 cm (17.7 in.) hood closed |
| Weight                  | 41 kg (90.2 lb.) standard configuration                                                                   |
| Electrical requirements | 100–240 V, 250 VA, 50–60 Hz, 5 A fuse                                                                     |
| Sample throughput       | 1–96 samples per run                                                                                      |
| Warranty                | 1 year warranty on parts and labor                                                                        |

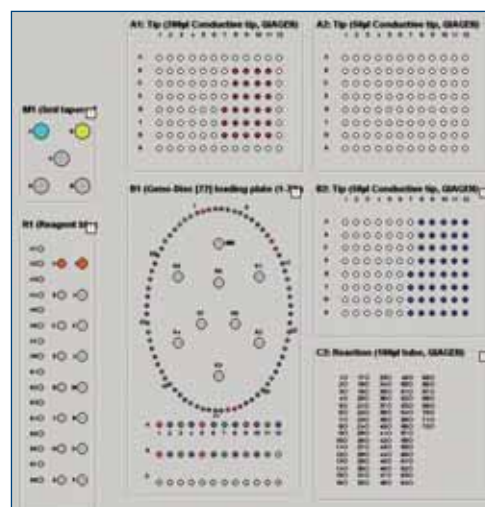

Visit [www.qiagen.com/goto/QIAgility](http://www.qiagen.com/goto/QIAgility) to find out more!

## Ordering Information

| Product                                  | Contents                                                                                                                                                                                        | Cat. no. |
|------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------|
| QIAgility System<br>(incl. PC)           | Robotic workstation for automated PCR setup (without UV light and HEPA filter),<br>1 year warranty on parts and labor, installation and training included                                       | 9001531  |
| QIAgility System<br>HEPA/UV (incl. PC)   | Robotic workstation for automated PCR setup (with UV light and HEPA filter),<br>1 year warranty on parts and labor, installation and training included                                          | 9001532  |
| Rotor-Disc 100<br>Starter Kit            | Kit includes: 2 Rotor-Disc 100 packs, Rotor-Disc Heat Sealer, Rotor-Disc Heat<br>Sealing Film, Rotor-Disc 100 Rotor and Locking Ring, Rotor-Disc 100 Loading<br>Block, Rotor-Disc Pipetting Aid | Inquire  |
| Rotor-Disc 72<br>Starter Kit             | Kit includes: 3 Rotor-Disc 72 packs, Rotor-Disc Heat Sealer, Rotor-Disc Heat<br>Sealing Film, Rotor-Disc 72 Rotor and Locking Ring, Rotor-Disc 72 Loading Block,<br>Rotor-Disc Pipetting Aid    | Inquire  |
| 50 µl Conductive<br>Filtered Tips (960)  | Conductive Filtered Tips for use with the QIAgility (960 x 50 µl)                                                                                                                               | 990512   |
| 200 µl Conductive<br>Filtered Tips (960) | Conductive Filtered Tips for use with the QIAgility (960 x 200 µl)                                                                                                                              | 990522   |
| Tip Receptacle<br>Box (10)               | Tip Receptacle Boxes for disposal of used tips from the QIAgility (10)                                                                                                                          | 990550   |
| 5 ml Tube; Graduated,<br>Flat-Base (50)  | Graduated, Flat-Base Tubes for use with the QIAgility (50 x 5 ml)                                                                                                                               | 990552   |
| Reagent Trough<br>Lid 70 ml (10)         | Reagent Trough Lid for use with the QIAgility (10 x 70 ml)                                                                                                                                      | 990554   |
| Reagent Trough<br>Lid 170 ml (20)        | Reagent Trough Lid for use with the QIAgility (20 x 170 ml)                                                                                                                                     | 990556   |
| Reagent Trough<br>Lid 270 ml (15)        | Reagent Trough Lid for use with the QIAgility (15 x 270 ml)                                                                                                                                     | 990558   |

The QIAgility is intended for molecular biology applications. This product is neither intended for the diagnosis, prevention, or treatment of a disease, nor has it been validated for such use either alone or in combination with other products. The Rotor-Gene Q, if used in combination with QIAGEN kits indicated for use with the Rotor-Gene Q instrument is intended for the applications described in the respective QIAGEN kit handbooks. If the Rotor-Gene Q instrument is used with kits other than QIAGEN kits, it is the user's responsibility to validate the performance of such product combination for any particular application. Rotor-Gene Kits, QuantiFast Kits, QuantiTect Kits, QIAGEN OneStep RT-PCR Kit, Type-it PCR Kits, QIAGEN LongRange Kits, HotStar HiFidelity Polymerase Kit, *cador* veterinary PCR Kits, HotStarTaq Kits, HotStarTaq Plus Kits, QIAGEN Taq Kits, QIAGEN Fast Cycling PCR Kits and TopTaq Master Mix Kit are intended for research use. No claim or representation is intended to provide information for the diagnosis, prevention, or treatment of a disease. The *cador* BVDV RT-PCR Kit is approved by the Friedrich-Loeffler-Institute in Germany for veterinary diagnostic use. The *cador* BVDV Type 1/2 RT-PCR Kit and the *cador* T. equigenitalis PCR Kit is under review by the Friedrich-Loeffler-Institute in Germany for veterinary diagnostic use. TopTaq DNA Polymerase is intended for general laboratory use. No claim or representation is intended to provide information for the diagnosis, prevention, or treatment of a disease.

Purchase of this product (Rotor-Gene, QuantiFast, and QuantiTect SYBR Green Kits) is accompanied by a limited, non-transferable immunity from suit to use it with detection by a dsDNA-binding dye as described in U.S. Patents Nos. 5,994,056 and 6,171,785 and corresponding patent claims outside the United States for the purchaser's own internal research. No real-time apparatus or system patent rights or any other patent rights, and no right to use this product for any other purpose are conveyed expressly, by implication or by estoppel.

This product (Rotor-Gene, QuantiFast and QuantiTect Probe and Multiplex Kits, and Type-it Probe Kits) is an Authorized 5' Nuclease Core Kit without Licensed Probe. Its purchase price includes a limited, non-transferable immunity from suit under certain patents owned by Roche Molecular Systems, Inc. or F. Hoffmann-La Roche Ltd, for using only this amount of the product in the practice of the 5' nuclease process solely for the purchaser's own internal research when used in conjunction with Licensed Probe. No right under any other patent claims (such as apparatus or system claims) and no right to use this product for any other purpose is hereby granted expressly, by implication or by estoppel. Further information on purchasing licenses may be obtained by contacting the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.

Purchase of this product (Rotor-Gene, QuantiFast, and QuantiTect RT-PCR Kits) is accompanied by a limited, non-transferable license under RT and Reverse Transcription-PCR patents owned by Roche Molecular Systems, Inc. and F. Hoffmann-La Roche Ltd to use it for the purchaser's own internal research. No real-time patent rights of any kind, no right under any other patent claims (such as apparatus or system claims), and no right to use this product for any other purpose is hereby granted expressly, by implication or by estoppel.

Purchase of this product (QuantiTect Kits) is accompanied by a limited license under U.S. Patent Numbers 5,035,996; 5,945,313; 6,518,026 and 6,287,823 and corresponding foreign patents.

The purchase of this product (QuantiTect SYBR Green Kits) includes a limited, non-transferable license under U.S. Patent No. 5,871,908 and all continuations and divisionals, and corresponding claims in patents and patent applications outside the United States, owned by Roche Diagnostics GmbH, for internal research use or for non-in vitro diagnostics applications with authorized reagents with regard to Melting Curve Analysis. No right is conveyed, expressly, by implication or estoppel, under any other patent or patent claims owned by Roche Diagnostics GmbH, or by any other Party.

Purchase of this product (TopTaq DNA Polymerase) is not accompanied by a license under patents owned by Roche Molecular Systems, Inc. and F. Hoffmann-La Roche Ltd. No rights to real-time PCR or to any primers, probes or pathogens are conveyed expressly, by implication or by estoppel.

The QIAGEN silver logo is exclusively licensed to Corbett Research.

#### NOTICE TO PURCHASER:

The purchase of this product (Rotor-Gene Q, Rotor-Disc) includes a limited, non-transferable license to certain patents (see details below) surrounding rapid polymerase chain reaction (PCR) methods and instrumentation, the use of SYBR® Green I in PCR reactions, melting curve analysis, analysis methods of DNA melting data, specifically high resolution melting (HRM) and others.

The purchase of this product (Rotor-Gene Q, Rotor-Disc) includes a limited, non-transferable license to one or more of US Patents Nos 6,787,338; 7,238,321; 7,081,226; 6,174,670; 6,245,514; 6,569,627; 6,303,305; 6,503,720; 5,871,908; 6,691,041; 7,387,887; and U.S. Patent Applications Nos. 2003-0224434 and 2006-0019253 and all continuations and divisionals, and corresponding claims in patents and patent applications outside the United States, owned by the University of Utah Research Foundation, Idaho Technology, Inc., and/or Roche Diagnostics GmbH, for internal research use or for non-in vitro diagnostics applications. No right is conveyed, expressly, by implication or estoppel, for any reagent or kit, or under any other patent or patent claims owned by the University of Utah Research Foundation, Idaho Technology, Inc., and/or Roche Diagnostics GmbH, or by any other Party. For information on purchasing licenses for in-vitro diagnostics applications or reagents, contact Roche Molecular Systems, 4300 Hacienda Drive, Pleasanton, CA 94588, USA.

#### ADDITIONAL NOTICE TO PURCHASER IN USA, EUROPE, AND AUSTRALIA ONLY:

Practice of the polymerase chain reaction (PCR) process may require a license depending on the reagents and methods applied. Purchase of this instrument (Rotor-Gene Q) does not itself convey to the purchaser a complete license or right to perform the PCR process. The use of this instrument (Rotor-Gene Q) with Authorized Reagents such as Authorized QIAGEN kits and assays provides a limited PCR license in accordance with the label rights accompanying such reagents.

Further information on purchasing licenses to practice the PCR process may be obtained by contacting the Director of Licensing at Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA and/or Roche Molecular Systems, 4300 Hacienda Drive, Pleasanton, CA 94588, USA. No right is conveyed expressly, by implication or by estoppel under any other patent claim, such as claims to apparatus, reagents, kits, or methods such as 5' nuclease methods.

[www.qiagen.com](http://www.qiagen.com)

**Australia** ■ Orders 03-9840-9800 ■ Fax 03-9840-9888 ■ Technical 1-800-243-066

**Austria** ■ Orders 0800/28-10-10 ■ Fax 0800/28-10-19 ■ Technical 0800/28-10-11

**Belgium** ■ Orders 0800-79612 ■ Fax 0800-79611 ■ Technical 0800-79556

**Brazil** ■ Orders 0800-557779 ■ Fax 55-11-5079-4001 ■ Technical 0800-557779

**Canada** ■ Orders 800-572-9613 ■ Fax 800-713-5951 ■ Technical 800-DNA-PREP (800-362-7737)

**China** ■ Orders 0086-21-3865-3865 ■ Fax 0086-21-3865-3965 ■ Technical 800-988-0325

**Denmark** ■ Orders 80-885945 ■ Fax 80-885944 ■ Technical 80-885942

**Finland** ■ Orders 0800-914416 ■ Fax 0800-914415 ■ Technical 0800-914413

**France** ■ Orders 01-60-920-920 ■ Fax 01-60-920-925 ■ Technical 01-60-920-930

**Germany** ■ Orders 02103-29-12000 ■ Fax 02103-29-22000 ■ Technical 02103-29-12400

**Hong Kong** ■ Orders 800 933 965 ■ Fax 800 930 439 ■ Technical 800 930 425

**Ireland** ■ Orders 1800 555 049 ■ Fax 1800 555 048 ■ Technical 1800 555 061

**Italy** ■ Orders 02-33430-420 ■ Fax 02-33430-426 ■ Technical 800-787980

**Japan** ■ Telephone: 03-6890-7300 ■ Fax 03-5547-0818 ■ Technical 03-6890-7300

**Korea (South)** ■ Orders 1544 7145 ■ Fax 1544 7146 ■ Technical 1544 7145

**Luxembourg** ■ Orders 8002-2076 ■ Fax 8002-2073 ■ Technical 8002-2067

**Mexico** ■ Orders 01-800-7742-639 ■ Fax 01-800-1122-330 ■ Technical 01-800-7742-639

**The Netherlands** ■ Orders 0800-0229592 ■ Fax 0800-0229593 ■ Technical 0800-0229602

**Norway** ■ Orders 800-18859 ■ Fax 800-18817 ■ Technical 800-18712

**Singapore** ■ Orders 65-67775366 ■ Fax 65-67785177 ■ Technical 65-67775366

**Spain** ■ Orders 91-630-7050 ■ Fax 91-630-5145 ■ Technical 91-630-7050

**Sweden** ■ Orders 020-790282 ■ Fax 020-790582 ■ Technical 020-798328

**Switzerland** ■ Orders 055-254-22-11 ■ Fax 055-254-22-13 ■ Technical 055-254-22-12

**UK** ■ Orders 01293-422-911 ■ Fax 01293-422-922 ■ Technical 01293-422-999

**USA** ■ Orders 800-426-8157 ■ Fax 800-718-2056 ■ Technical 800-DNA-PREP (800-362-7737)

