



Superior sectioning results – **perfectly simple**

Epredia Rotary and Sliding Microtomes
For research and clinical laboratories





Section safely with **excellent quality**

Microtomes from EpreDia are designed to be precise and stable, helping to yield superior sectioning results for even the most difficult-to-cut paraffin sections.

Designed with the end user in mind

- With both automatic and manual options available, we have the right fit for you, your lab, and your paraffin sections

Specificity in action

- User-friendly control panels facilitate flexibility
- Fully-adjustable to ensure the correct thickness

Non-tiring comfort

- Ergonomic designs with integrated waste trays to save you time and energy
- Emergency stop features to keep you safe (EpreDia HM355S only)

The total solution

- Wide variety of steel low- and high-profile blades for cutting various tissue types
- EpreDia Cool-Cut and EpreDia™ Section Transfer System™ (STS) assist in placing cooled, wrinkle-free sections on glass microscope slides
- Get peace of mind with EpreDia service contracts that include preventative maintenance and thorough, expedited repairs, ensuring minimal downtime in your laboratory



Lean sectioning

Improve workflow at the microtome by transferring sections directly onto slides as they are generated. Combining an EpreDia microtome with a water bath and EpreDia™ SlideMate™ AS or SlideMate Pro slide printer offers straight-forward, cost-effective means of simultaneously streamlining workflow and reducing the risk of sample identification errors.

Rotary Microtomes

Designed for precision and stability to help labs achieve quality ribbons even with difficult-to-cut paraffin sections.



Manual

Epredia HM 325

A premier manual rotary microtome with additional capability to support rapid set-up and cleaning for improved workflow.

- 64 mm vertical stroke can accommodate Epredia™ Super Mega™ Cassettes in horizontal orientation
- Wrap-around, detachable waste tray for safe and speedy debris collection and disposal
- XY fine orientation and zero positioning for rapid re-orientation of pre-cut blocks
- Ultra-light touch, ergonomic flywheel
- Optional Epredia Section Transfer System supports workflow transition from ribbon to slide
- Compatible with Epredia Cool Cut to extend cutting periods for efficient step and serial sectioning

HM 325 manual microtome	902100
• With E blade holder	902100A
• With ER blade holder	902100ER


Semi-automated

Epredia HM 340E

A semi-automated rotary microtome combining mechanized feed with manual cutting.

- 72 mm vertical stroke for efficient sectioning of macro and Super Mega Cassettes
- Memory function for specimen positioning with XY fine orientation and wraparound waste tray to speed up preparation and clean-up
- Ergonomic removable control panel to suit user preferences
- Compatible with Epredia Section Transfer System and Epredia Cool Cut

HM 340E semi-automated microtome	905190
• With E blade holder	905190A
• With ER blade holder	905190ER

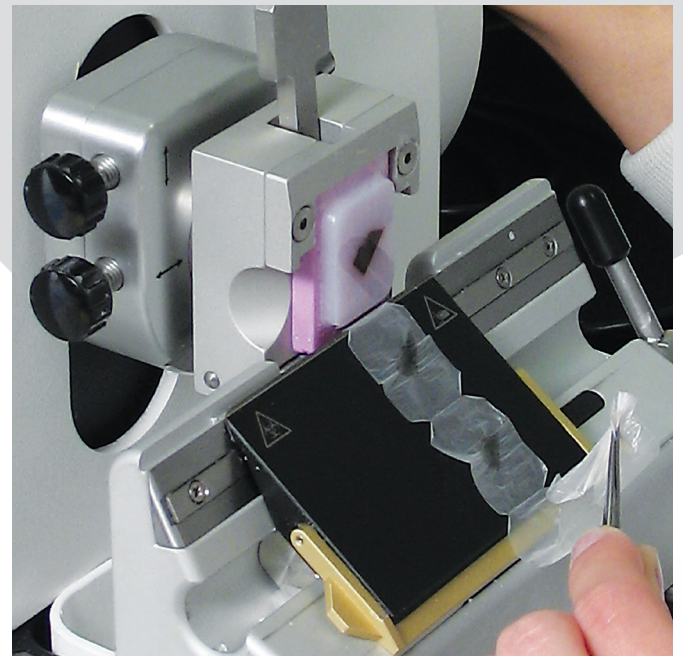

Fully automated

Epredia HM 355S

A fully automated microtome with four mechanized cutting modes to deliver exceptional sections across an extensive range of specimen types.

- 72 mm vertical stroke for efficient sectioning of Super Mega Cassettes and other large samples
- Choice of manual or mechanized sectioning
- Single, multiple, interval and continuous cutting modes for optimum control according to sectioning requirements
- Unique “double-tap” startup is designed to prevent accidental operation
- Compatible with Cool-Cut and Section Transfer System for improved workflow
- Programmable cutting window to define the sectioning area for more efficient sectioning
- Memory function for specimen positioning with XY fine orientation and wraparound waste tray to speed up preparation and clean-up
- Ergonomic removable control panel to suit user preferences
- Emergency stop, manual, and electronic brake

HM 355S fully automated microtome 905200

- With E blade holder 905200A
 - With ER blade holder 905200ER
-

Histological examination continues to be the gold standard for detection, diagnosis and characterization of many clinical conditions



Sliding Microtomes

Epredia sliding microtomes – smart choices for reliability, comfort and versatility.




Manual

Epredia HM 430

A manually operated microtome suitable for a range of botanical samples in addition to anatomical tissue.

- Variable stroke length to meet specimen requirements
- Accommodates specimens up to 80 mm x 60 mm
- Coaxial specimen orientation and memory function for rapid re-orientation of pre-cut blocks
- Choice of manual or automated advance

HM 430 manual microtome	910010
• With low-profile blade carrier	910010L
• With high-profile blade carrier	910010H

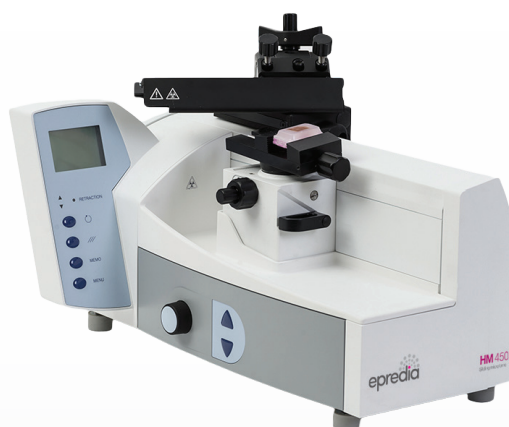

Fully automated

Epredia HM 450

An automated sliding microtome with additional capability for cutting large and hard specimens.

- Choice of manual or mechanized operation
- Coaxial specimen orientation and memory function for rapid re-orientation of pre-cut blocks
- User-friendly control panel for easy operation
- Optional retraction to protect specimen

HM 450 fully automated microtome	910020
• With low-profile blade carrier	910020L
• With high-profile blade carrier	910020H

Essential accessories

Epredia microtomes can be custom-configured for optimal application-specific performance.



Epredia Cool-Cut and Section Transfer System

Discover the efficiencies of a complete sectioning solution. Simply pair one of our Epredia rotary microtomes with the Epredia Cool-Cut paraffin block cooler, and Epredia Section Transfer System. These accessories are designed to help laboratories improve quality and efficiency – two key influencers of positive patient outcomes.

The Cool-Cut provides object cooling for high-quality sectioning and process efficiency. The Section Transfer System automatically transfers section ribbons from the blade to the water bath reducing manual handling. Their compact, integrated design saves valuable bench space.

Both are designed to be used exclusively with the Epredia HM 325, HM 340E, and HM 355S rotary microtomes.



Clamps and blade holders

Our variety of clamps and blade holders allow users to fine-tune the performance of their microtomes to precisely match specific specimen and media types, section thicknesses and cutting styles.



Microtome blades

Epredia high- and low-profile disposable blades deliver longevity, reliability and ribbon generation across a range of tissue types.

4481 Campus Drive
Kalamazoo, MI 49008
UNITED STATES
+1 (800) 522-7270

Tudor Road, Manor Park
Runcorn, WA7 1TA
UNITED KINGDOM
+44 (0) 800 018 9396
+44 (0) 1928 534 000

epredia Enhancing precision
cancer diagnostics

Find out more at [epredia.com](https://www.epredia.com)

POWERED BY SHANDON MICRON MENZEL-GLÄSER