

LABGARD®




NU-425 CLASS II, TYPE A2
BIOLOGICAL SAFETY CABINETS




Best Products. Best Performance. Best Protection.

The LabGard® Class II Type A2 Biological Safety Cabinet can be either exhausted back into the room or connected to a facility HVAC system. The cabinet's airflow is 30% exhausted / 70% recirculated to minimize cross-contamination of low to moderate risk biologicals in the absence of volatile toxic chemicals.

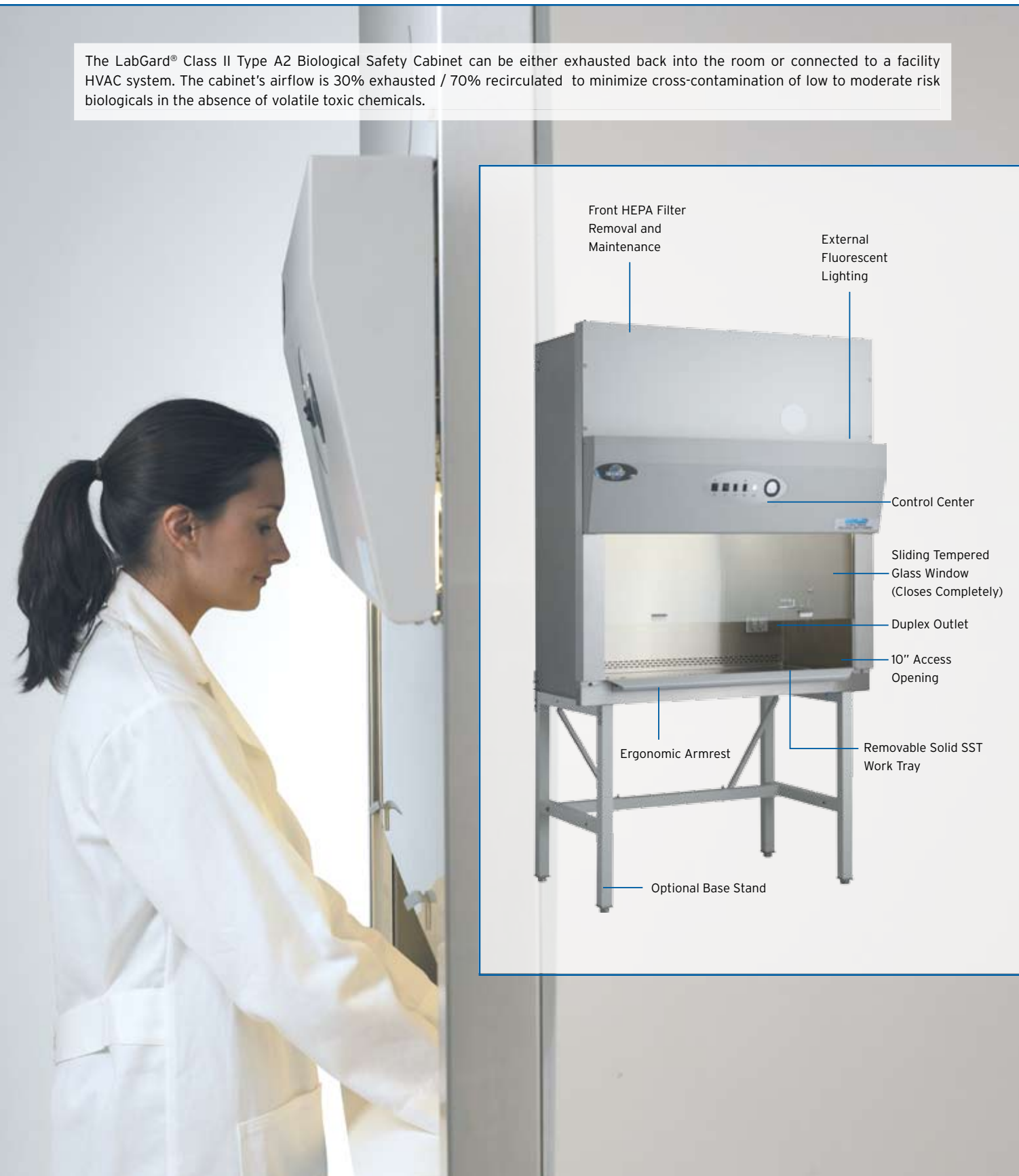

Improved Principles of Biological Safety Cabinet Design

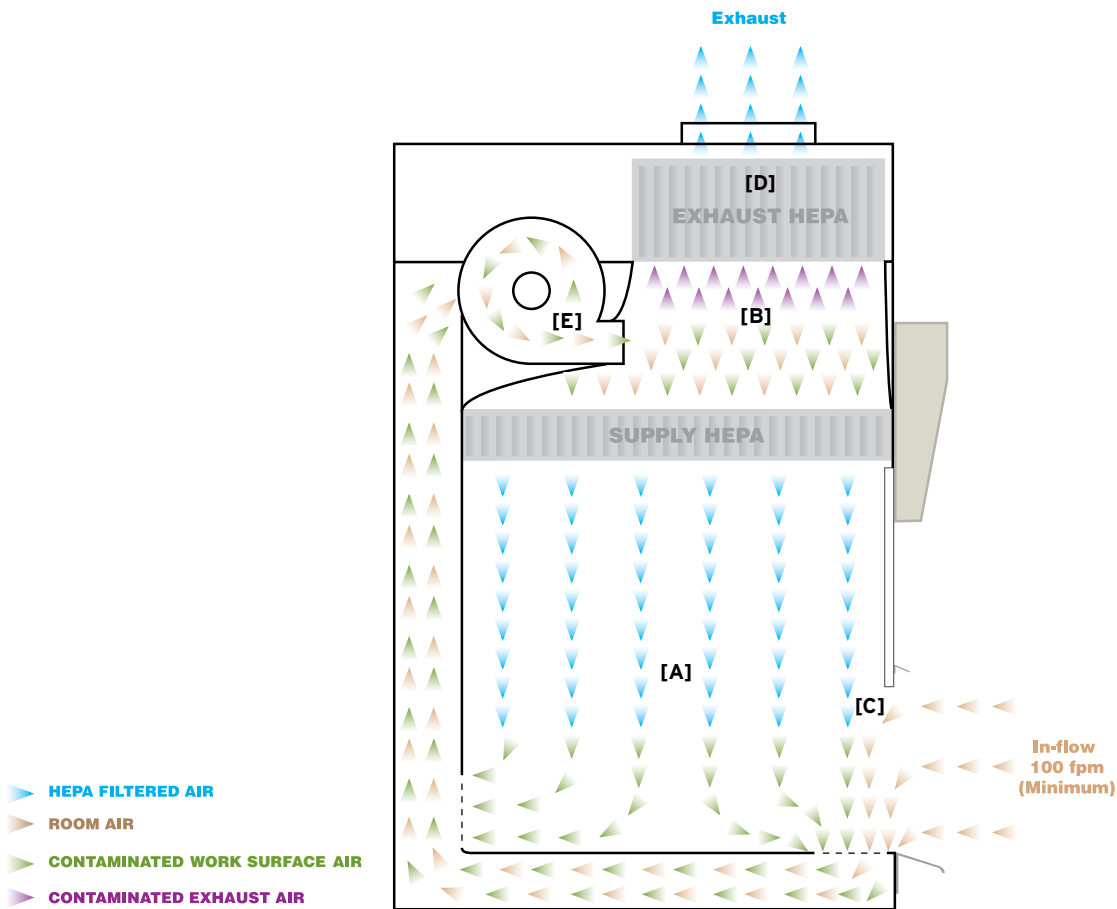

Class II Type A2 Biological Safety Cabinet Airflow Schematic

True Laminar Airflow [A]

Unidirectional airflow moving along parallel flow lines at a constant velocity minimizes air turbulence within the work zone.

HEPEX™ Zero Leak Airflow System [B]

The HEPEX™ Zero Leak Airflow System eliminates the possibility of cabinet and gasket leaks.

Greater Volume of Air Creates a Stronger Air Barrier [C]

Air exchanges within the Class II Type A2 Biological Safety Cabinet occur 25 times per minute.

Larger, Long Lasting HEPA Filters [D]

NuAire® utilizes HEPA filters that are typically 25% larger than those in competitive products.

Motor/Blower System Extends HEPA Filter Life [E]

As the HEPA filter loads the motor/blower system automatically compensates to assure performance and extended filter life.

Single Piece Stainless Steel Construction

Centrally Located Control Center

State-of-the-Art Electronic Safety System

Improved Quality, Reduced Noise

Biological Tolerance Testing Assures Optimal Performance



Registered to: ISO 9001:2000
ISO 13485:2003



For more information contact NuAire, Inc.® at 1.800.328.3352 or www.nuaire.com

NU-425 Control Center



NU-425 (Class II Type A2)

Cabinet controls include a minihelic pressure gauge, motor voltage regulator, audible window alarm, main and outlet power circuit breakers, and power switches for exterior mounted fluorescent lights and / or ultraviolet lights, interior outlets, and blower motor.

LABGARD® NU-425 Class II, Type A2 Biological Safety Cabinet

Standard Features

- Personnel, Product, and Environmental Protection
- True Laminar Airflow
- 30% Exhaust - 70% Recirculation
- HEPEX™ Zero Leak Airflow System
- Large HEPA Filters: 99.99% Efficiency Rating
- Metal Diffuser over Supply Filter
- Removable Control Center
- Spill Trough with Drain Valve
- 100% Stainless Steel Sealed Construction
- Minimum Air Turbulence
- 100 LFPM Air Barrier
- Sliding Window with 10" Work Access Opening Window is in Vertical Position, Minimizing Glare and Allows for Ergonomically Correct Posture (NU-425-300, 8" Opening)
- Frameless Polished Edge Window, 19" Maximum Opening
- Removable Work Tray
- External Non-Glare Fluorescent Lighting
- Electrical Duplex Outlets
- One Gas Service Petcock (Right Side Wall)
- Minihelic Gauge
- Ergonomic PVC Armrest
- NSF Standard 49 and cUL Approved

Optional Features and Accessories


- Germicidal Ultraviolet Lamp
- Prop-Up Removable Work Tray
- Ground Fault Interrupter
- IV Bar with Six (6) Hooks
- Base Support Stand
- Footrest
- Elbow Rest
- Microscope Viewing Window
- Additional Service Valves for Gas, Air, or Vacuum
- Exhaust Transitions
- Storage Pull-Out Trays
- Alarm Systems [exhaust]
- Adjustable Control for Support Stand or Storage Cabinet
- Decorative Side Panels
- Sinks with Hot / Cold or DI Water Faucets
- Heated Work Surface
- Seismic Studs and Brackets

Model	Size (Nominal)	Electrical (All Sizes)	Interior Dimensions (W x H x D) in. [mm]	Exterior Dimensions (W x H x D) in. [mm]
NU-425-300	3 ft.	115 VAC / 60 Hz	34 ³ / ₈ x 28 ¹ / ₂ x 23 ¹ / ₂ [873 x 724 x 597]	41 ⁵ / ₈ x 63 x 32 ⁷ / ₈ [1057 x 1600 x 835]
NU-425-400	4 ft.	D: 100 VAC / 50-60 Hz	46 ³ / ₈ x 28 ¹ / ₂ x 23 ¹ / ₂ [1178 x 724 x 597]	53 ⁵ / ₈ x 63 x 32 ⁷ / ₈ [1362 x 1600 x 835]
NU-425-500	5 ft.	E: 230 VAC / 50 Hz	58 ³ / ₈ x 28 ¹ / ₂ x 23 ¹ / ₂ [1483 x 724 x 597]	65 ⁵ / ₈ x 63 x 32 ⁷ / ₈ [1667 x 1600 x 835]
NU-425-600	6 ft.	G: 220 VAC / 60 Hz	70 ³ / ₈ x 28 ¹ / ₂ x 23 ¹ / ₂ [1788 x 724 x 597]	77 ⁵ / ₈ x 63 x 32 ⁷ / ₈ [1972 x 1600 x 835]




For more information contact NuAire, Inc.® at 1.800.328.3352 or www.nuaire.com

Options and Accessories

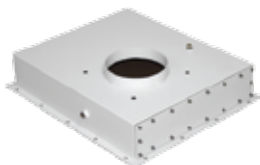

Gas Tight / Low Profile
Exhaust Transition

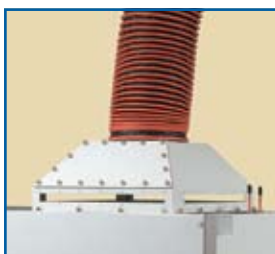

Canopy Exhaust Transition
and Flex Duct



Butterfly Valves



203 mm to 254 mm
Round Transition



Neoprene Connection Kit



Exhaust Airflow Monitors




Elbow Rest



IV Bar With Hooks



Remote Controlled
Service Valves

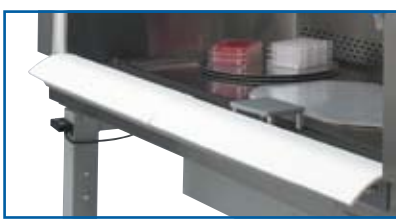

Foam Armrest Pad



Stainless Steel Turntable



Adjustable Footrest



Motorized Adjustable Base Stand



Base Stand with Storage Cabinet

Also Available:

Ultraviolet Germicidal
Lamps
Faucets and Sinks, Cup Sinks
Pre-Piped Plumbing
Ground Fault Interrupter
Side Panels

Seismic Brackets
Metal Frame HEPA Filter
Microscope Window
Heated Work Surface
Custom Configurations



For more information contact NuAire, Inc.® at 1.800.328.3352 or www.nuaire.com



NuAire, Inc.®

2100 Fernbrook Lane

Plymouth, MN 55447

U.S.A.

1.800.328.3352

www.nuaire.com