

Shakers

GFL®



- For gentle motion, vigorous mixing or intensive shaking - you can choose between 15 quality products and five different types of shaking motion.

Space-saving and versatile

What sets this robust, space-saving compact Shaker apart from the 3005 Shaker is its shaking motion. Its impressive features are both its gentle back-and-forth motion and its intensive shaking of liquids, especially in separating funnels. This unit is maintenance-free and carries the CE mark.

3006



3006 Analogue Reciprocating Shaker with Tray 3950 and Clamps for separating funnels 3955 - 3959

Technical Data

Overall dimensions (W x D x H):
380 x 510 x 140 mm

Moving platform: 330 x 330 mm

Load: max. 8 kg

Shaking motion: reciprocating
(back and forth)

Timer: up to 60 minutes, or continuous
operation

Shaking amplitude: 20 mm

Shaking frequency: 20 - 300 rpm

Electrical connection:

230 V* 50...60 Hz, 65 W

* Other voltages on request!

Net/gross weight 10/12 kg

Packing volume (cardboard box): 0.1 m³

Order No. 3006

Features

- ▶ compact, low-wear counterbalanced drive mechanism, ensuring high stability and dependability for continuous operation
- ▶ housing made of high impact strength polystyrene, off-white varnish. The base plate is made of electrolytically galvanised, powder-coated sheet steel, the shaking platform of anodised aluminium, equipped with four plastic pins for secure attachment of accessory equipment
- ▶ an a.c. motor with overload protection drives the unit
- ▶ clearly laid out control panel for easy operation

- ▶ electronic speed control, stepless, with gentle start-up
- ▶ constant speed during continuous operation, independent of load

Applications

This Shaker can be used in medical diagnostics or analysis, as well as for growing micro-organisms, cell and tissue cultures.

This unit is used particularly in laboratories, incubation rooms and moderating rooms at ambient temperatures between +10 °C and +50 °C.

Robust and variable

Versatile, powerful and exceptionally quiet - those are the descriptive qualities of this robust Shaker with reciprocating (back-and-forth) motion. Use of this device is recommended for vigorous and intensive shaking of clamped lying flasks, mixing flasks or separating funnels in continuous operation. Both of these units are maintenance-free and carry the CE mark.

3016 / 3018



Technical Data

Overall dimensions (W x D x H):
510 x 625 x 145 mm

Moving platform: 450 x 450 mm

Load: max. 15 kg

Shaking motion: reciprocating

Timer:

3016 / up to 60 minutes, or continuous operation

3018 / 1 min. - 99:59 hours, or continuous operation

Shaking amplitude: 30 mm

Shaking frequency: 20 - 300 rpm

Electrical connection: 230 V*

* Other voltages on request!

3016 / 50...60 Hz, 65 W

3018 / 50 or 60 Hz**, 65 W

** Please indicate when ordering!

Net/gross weight: 18/22 kg

Packing volume (cardboard box): 0.19 m³

■ **Order No. 3016**

■ **Order No. 3018**

with RS 232 port



Features

- ▶ compact, low-wear counterbalanced drive mechanism, ensuring high stability and dependability for reciprocal shaking motion
- ▶ housing made of high impact strength polystyrene, off-white varnish. The base plate is made of electrolytically galvanised, powder-coated sheet steel, the shaking platform of anodised aluminium, equipped with four plastic pins for secure attachment of accessory equipment
- ▶ an a.c. motor with overload protection drives the unit
- ▶ 3016 - clearly laid-out control panel for easy operation
3018 - speed and remaining time indicated by two LC displays

- ▶ 3016 - electronic speed control, stepless, gentle start-up
3018 - microprocessor-controlled, adjustable in steps of 1.0 rpm, gentle start-up
- ▶ constant speed during continuous operation, independent of load

Applications

The producers of suspensions and emulsions know this shaker well. The shaking platform's jerky back-and-forth motion creates high turbulence, thus thoroughly mixing the media. This is a vital prerequisite, for example, in medical diagnostics, or in foodstuffs or environmental analyses.

Both of these units are well-suited for use in laboratories, incubation rooms and moderating rooms in ambient temperatures between +10 °C and +50 °C.



3016 Analogue Back-and-forth Shaker
with Tray 3966 and Test Tube Racks 3953

Accessories



3005



3006



Shaking tray

made of stainless steel, 330 x 330 mm, with holes to accommodate clamps for Erlenmeyer flasks and other accessory equipment

Order No. 3950




Clamps

made of stainless steel for Erlenmeyer flasks, for screwing on to shaking tray 3950, complete with attachment materials

- Order No. 3983** for 25 ml flasks (45*)
- Order No. 3984** for 50 ml flasks (25*)
- Order No. 3985** for 100 ml flasks (16*)
- Order No. 3986** for 200 ml flasks (12*)
- Order No. 3987** for 250-300 ml flasks (9*)
- Order No. 3988** for 500 ml flasks (9*)
- Order No. 3989** for 1000 ml flasks (4*)
- Order No. 3990** for 2000 ml flasks (2*)

Other clamps on request

* max. number of clamps per shaking tray

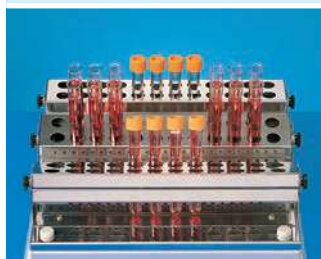


Test tube rack

made of stainless steel, for test tubes; perforated section can be tilted by 90°, with clamping springs for secure grip and silent shaking of glasses; can be screwed on to shaking tray 3950

Order No. 3953

for a max. of 24 tubes: 12 - 17 mm diameter
Capacity of tray: three test tube racks

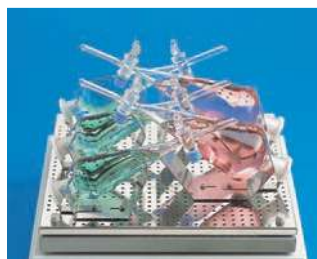

Clamps

made of stainless steel, for separating funnels, for screwing on to shaking tray 3950

Order No. 3957 for 50 ml separating funnels (6*)
(ISO/Squibb, neck 19/26)

Order No. 3958 for 100 ml separating funnels (6*)
(ISO/Squibb, neck 19/26)


Order No. 3959 for 250 ml separating funnels (4*)
(ISO/Squibb, neck 19/26)

Order No. 3955 for 250 ml separating funnels (4*)
(conical form, neck 29/32)

Order No. 3956 for 500 ml separating funnels (3*)
(ISO/Squibb, neck 29/32)

Other clamps on request

* Max. number of clamps per shaking tray

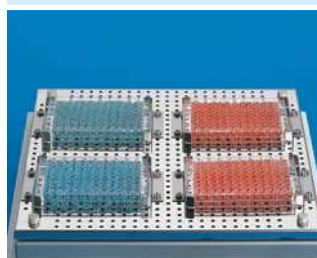


Holding frame

made of stainless steel, for attachment of a test plate, for screwing on to shaking tray 3950, complete with attachment materials

Order No. 3910

Capacity of tray: four test plates

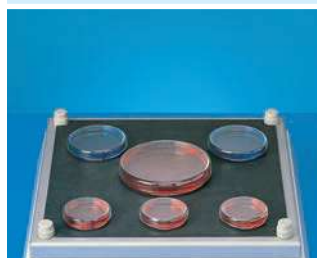

Non-slip mat

for the shaking platform, 300 x 300 mm size, for slow moving, e.g. nutrient solutions in Petri dishes and flasks

Order No. 3951



Universal mount

for secure attachment of different shaking objects between the four rubber-coated bars

Order No. 3952

Accessories



3011



3012



3013



3014



Order No. 3965

Non-slip mat


for the shaking platform, 420 x 420 mm size, for slow moving, e.g. nutrient solutions in Petri dishes



Order No. 3967

Universal mount

for secure attachment of different shaking objects between the six rubber-coated bars

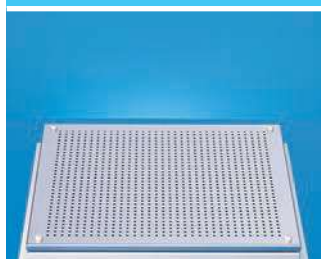

Order No. 3966

Shaking tray

made of stainless steel, 450 x 450 mm, with holes to accommodate clamps for Erlenmeyer flasks and other accessory equipment



Order No. 3968

Platform frame

with four levels, with top rims, three levels removable. Frame and trays made of stainless steel. Distance between trays - 60 mm. Grooved rubber mats on trays to ensure secure grip of objects.



Clamps

made of stainless steel for Erlenmeyer flasks, for screwing on to shaking tray 3966, complete with attachment materials

- Order No. 3983** for 25 ml flasks (79*)
- Order No. 3984** for 50 ml flasks (49*)
- Order No. 3985** for 100 ml flasks (36*)
- Order No. 3986** for 200 ml flasks (22*)
- Order No. 3987** for 250-300 ml flasks (16*)
- Order No. 3988** for 500 ml flasks (12*)
- Order No. 3989** for 1000 ml flasks (9*)
- Order No. 3990** for 2000 ml flasks (4*)

Other clamps on request

* Max. number of clamps per shaking tray

Accessories



3015




3017



3016



3018

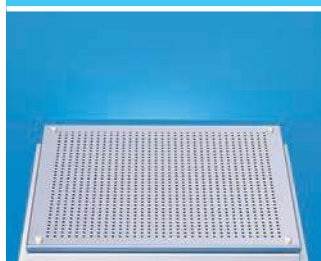

Shaking tray

made of stainless steel, 450 x 450 mm, with holes to accommodate clamps for Erlenmeyer flasks and other accessory equipment

Order No. 3966




Clamps

made of stainless steel for Erlenmeyer flasks, for screwing on to shaking tray 3966, complete with attachment materials

- Order No. 3983** for 25 ml flasks (79*)
- Order No. 3984** for 50 ml flasks (49*)
- Order No. 3985** for 100 ml flasks (36*)
- Order No. 3986** for 200 ml flasks (22*)
- Order No. 3987** for 250-300 ml flasks (16*)
- Order No. 3988** for 500 ml flasks (12*)
- Order No. 3989** for 1000 ml flasks (9*)
- Order No. 3990** for 2000 ml flasks (4*)

Other clamps on request

* Max. number of clamps per shaking tray

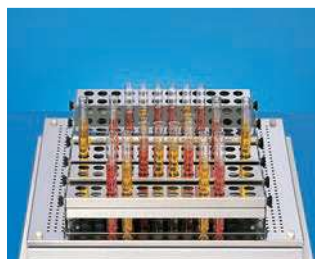

Test tube rack

made of stainless steel, for test tubes; perforated section can be tilted by 90°, with clamping springs for secure grip and silent shaking of glasses; can be screwed on to shaking tray 3966

Order No. 3953

for a max. of 24 tubes: 12 - 17 mm diameter
capacity of tray: six test tube racks



Clamps

made of stainless steel, for separating funnels, for screwing on to shaking tray 3966

Order No. 3957 for 50 ml separating funnels (11*)
(ISO/Squibb, neck 19/26)

Order No. 3958 for 100 ml separating funnels (11*)
(ISO/Squibb, neck 19/26)

Order No. 3959 for 250 ml separating funnels (8*)
(ISO/Squibb, neck 19/26)

Order No. 3955 for 250 ml separating funnels (8*)
(conical form, neck 29/32)

Order No. 3956 for 500 ml separating funnels (6*)
(ISO/Squibb, neck 29/32)

Other clamps on request

* Max. number of clamps per shaking tray

Holding Frame

made of stainless steel, for attachment of a test plate, for screwing on to shaking tray 3966, complete with attachment materials

pict. see Shakers 3019 - 3020

Order No. 3910

Capacity of tray: six test plates



Non-slip mat

for the shaking platform, 420 x 420 mm size, for slow moving, e.g. nutrient solutions in Petri dishes

Order No. 3965



Universal mount


for secure attachment of different shaking objects between the six rubber-coated bars

Order No. 3967

Accessories



3019

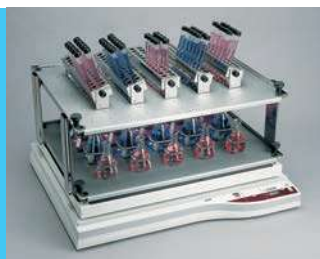

3020



Order No. 3980

Shaking tray

made of 3-mm-thick anodised aluminium, 670 x 537 mm, with holes to accommodate clamps for Erlenmeyer flasks and other accessory equipment




Clamps

made of stainless steel for Erlenmeyer flasks, for screwing on to shaking tray 3980, complete with attachment materials

- Order No. 3983** for 25 ml flasks (99*)
- Order No. 3984** for 50 ml flasks (99*)
- Order No. 3985** for 100 ml flasks (50*)
- Order No. 3986** for 200 ml flasks (26*)
- Order No. 3987** for 250-300 ml flasks (26*)
- Order No. 3988** for 500 ml flasks (26*)
- Order No. 3989** for 1000 ml flasks (12*)
- Order No. 3990** for 2000 ml flasks (9*)

Other clamps on request

* Max. number of clamps per shaking tray

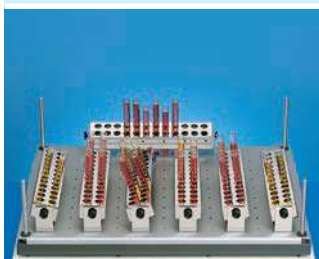

Order No. 3953

for a max. of 24 tubes: 12 - 17 mm diameter
Capacity of tray: nine test tube racks

Test tube rack

made of stainless steel for test tubes; perforated section can be tilted by 90°, with clamping springs for secure grip and silent shaking of glasses; can be screwed on to shaking tray 3980



Order No. 3910

Capacity of tray: fifteen test plates

Holding frame

made of stainless steel, for attachment of a test plate, for screwing on to shaking tray 3980, complete with attachment materials



Order No. 3982

Universal mount

for secure attachment of shaking objects. The mount consists of two longitudinal bars and six transverse bars. Screw-type clamps are used for adjusting the longitudinal bars vertically and transverse bars horizontally.



Order No. 3981

Rack frame

for doubling the capacity by accommodating two shaking trays 3980; distance between trays - 205 mm. The bottom tray can be loaded with flasks up to a max. of 500 ml.