



Europe Only

RUO model

NanoZoomer[®] S360

Digital slide scanner C13220-04

The latest high-throughput model, ideal for uses in pharmaceutical companies and academic research.

High-throughput
scanning

More than 82 slides/h
(40× mode)

Low-operational
workload

Automated assistant of
image quality check

High-capacity
scanning

360 slides
in one batch

HAMAMATSU
PHOTON IS OUR BUSINESS

High-throughput scanning of tissue slides with low operational workload

High-throughput and high-capacity scanning

By improving scan speed as well as other processes such as slide loading and data transfer.

High-throughput of 82 slides/h!

The drastic improvement of the scan speed, one half of conventional models, realizes throughput of 82 slides/hour for both 20× and 40× mode.

High-throughput

82 slides/h
(40× mode)

Automatic scan up to 360 slides!

Up to 30 tissue slides are mounted in a cassette, and up to 12 cassettes are mounted in a system. Total 360 slides are automatically scanned once you started.

Automatic scanning


360 slides

Daily yield of scans

1080 slides

Three runs of batches a day!

4.5 hours scanning time of 360 slides enable to start three batches within working hour.

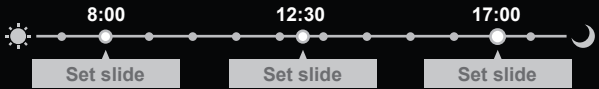


Automated assistant of image quality check

Greatly simplified image quality check process with automated focus evaluation.

Reduces the workload of checking image quality

Easy identification of slides that need to be visually checked and rescanned!

Focus quality of scanned images are scored, and it is presented on a display monitor. Users can identify slides need to be visually checked, then define slides need to be rescanned. This process improve scanning efficiency and greatly reduces operational workload.

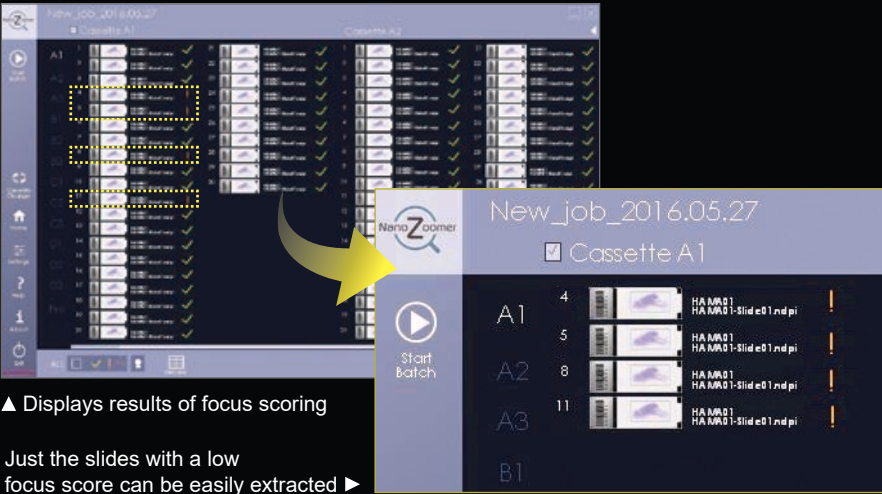


▲ Displays results of focus scoring

Just the slides with a low focus score can be easily extracted ►

Easy identification of areas in a slide that need to be visually checked!

The focus pass/fail results are superimposed and displayed over entire tissue on a slide, and users can easily identify areas that need to be checked.

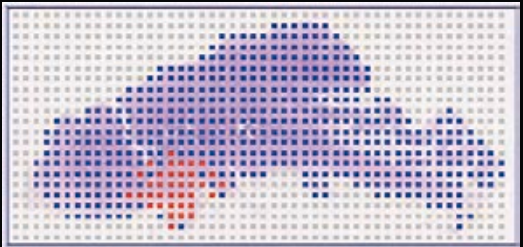

◀ Blue colored: in-focus area
Red colored: out-of-focus area

More productive and convenient

Scan process monitoring

Users can check progress of slide scans. Display panel shows status of each cassette as "Waiting for scan", "Scanning" and "Scan completed".



Cassette based management of slide scan mode

Slide scan mode is independently manageable for each cassette labeled with a unique barcode. It is useful when different kind of tissues or stains included in one batch of scan.




Barcode management for each cassette



*The NanoZoomer S360 is designed for use with a dedicated cart. Images above shown without the dedicated cart.

System configuration

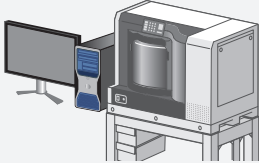

NanoZoomer S360
Digital slide scanner
C13220-04

+



NZAcquire
Image acquisition software
U14297-01

+



NDP®.view2
Image viewing software
U12388-01


+

Option software

- NDP.view2 Plus
Image viewing software
- NZConnect Serve Basic
Web slide server software
- Barcode option
Data matrix for NZAcquire

Dimensional outlines (Unit: mm)

● Weight
Main unit: Approx. 117 kg
Dedicated rack: Approx. 78 kg

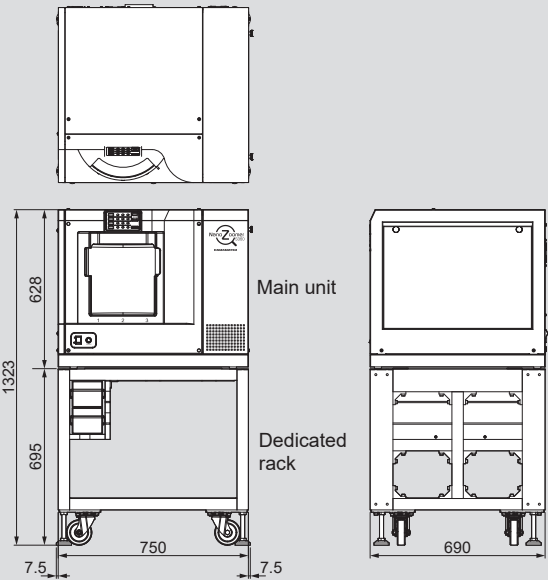

Main unit

Dedicated rack

* Excluding leveling feet.

Specifications

Product name		NanoZoomer S360
Product number		C13220-04
Scanning speed	20× mode (15 mm×15 mm)	Approx. 30 s
	40× mode (15 mm×15 mm)	Approx. 30 s
Throughput	20× mode (15 mm×15 mm)	More than 82 slides/h ^{*1}
	40× mode (15 mm×15 mm)	More than 82 slides/h ^{*1}
Objective lens		20× NA 0.75 User can select 20× or 40× mode at start of scanning
Compatible glass slide		26 mm×76 mm (Thickness 0.9 mm to 1.2 mm)
Slide loader	Standard size slide	360 slides (30 slides×12 cassettes)
Scanning resolution	20× mode	Approx. 0.46 μm/pixel
	40× mode	Approx. 0.23 μm/pixel
Focusing method		Pre-Focus map
Z-stack function		Included
Image compression		JPEG compression
Power supply		AC 100 V to AC 240 V
Power consumption (Scanner only)		Approx. 200 VA

*1 For the case of 5 focus points

• NanoZoomer and NDP are registered trademarks of Hamamatsu Photonics K.K. (EU, Japan, UK, USA).
• Subject to local technical requirements and regulations, availability of products included in this brochure may vary. Please consult your local sales representative.
• The product described in this brochure is designed to meet the written specifications, when used strictly in accordance with all instructions.
• Specifications and external appearance are subject to change without notice.
© 2026 Hamamatsu Photonics K.K.