


## UV Crosslinkers and Translinker

### UV Crosslinkers

UVP Crosslinkers assure consistent UV output for many applications including DNA bonding and UV curing. UV crosslinking for attaching nucleic acids to a membrane for example, takes seconds as compared to oven baking. The crosslinkers utilize a microprocessor to control and measure the dose of UV radiation.

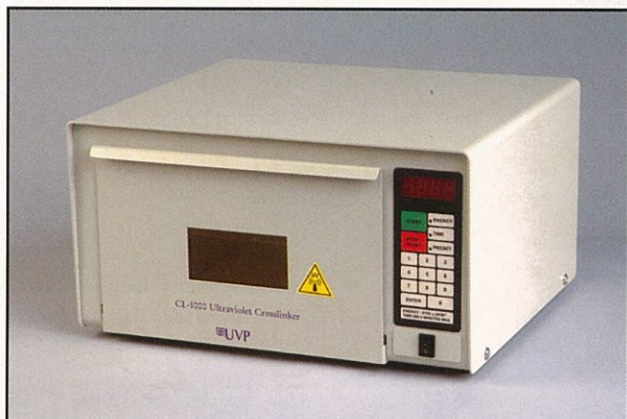


**CL-1000 Crosslinker** is an exposure instrument which uses shortwave UV energy to bond the DNA to a medium.

- Preset and manual controls for ultraviolet or time exposures
- Factory presets deliver 120,000 microjoules or 5 minutes of exposure; maximum UV energy setting of 999,900 microjoules per cm<sup>2</sup>
- Microprocessor optimizes the use of energy while ensuring consistent, maximum energy efficiency and output
- Large LED display
- Window on the door allows viewing of the process yet blocks the UV radiation
- Five 8-watt ultraviolet tubes for uniform overhead ultraviolet illumination



**Translinker** An ultraviolet blocking cover is included with the transilluminator.

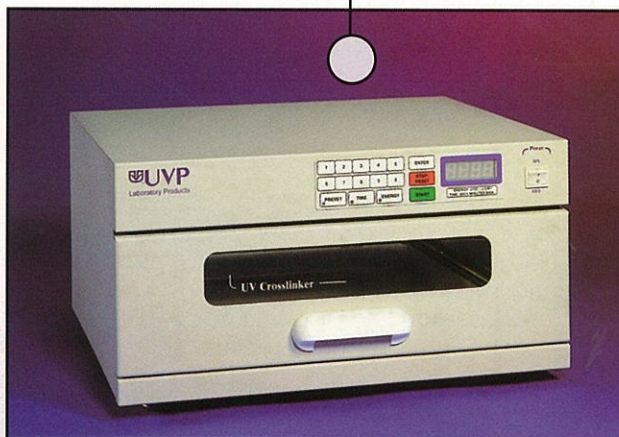

**CX-2000** Features a pull-out drawer and shortwave UV. CX-2000 operates with the same functions as the CL-1000 with microprocessor controls.

CX-2000 Dimensions:

Exterior: 15D x 17.5W x 9H in. (445 x 229 x 381mm)  
Interior: 11.25D x 11.5W in. (292 x 286mm)

CL-1000 Dimensions:

Exterior: 13.7D x 15.6W x 8.75H in. (248 x 396 x 222mm)  
Interior: 12D x 10W x 5H in. (305 x 254 x 127mm)

#### Ordering Information

|                                                                                                                                                         |                                                                        |
|---------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------|
| <b>CL-1000 Crosslinker</b><br>254nm UV<br>Replacement Tube                                                                                              | 95-0174-01 (115V)<br>95-0174-02 (230V)<br>34-0007-01                   |
| <b>CL-1000L Crosslinker</b><br>365nm UV<br>Replacement Tube                                                                                             | 95-0228-01 (115V)<br>95-0228-02 (230V)<br>34-0006-01                   |
| <b>CL-1000M Crosslinker</b><br>302nm UV<br>Replacement Tube                                                                                             | 95-0230-01 (115V)<br>95-0230-02 (230V)<br>34-0042-01                   |
| <b>CX-2000 Crosslinker</b><br>254nm UV<br>Replacement Tube                                                                                              | 95-0339-01 (115V)<br>95-0339-02 (230V)<br>34-0007-01                   |
| <b>TL-2000 Translinker</b><br>254nm Crosslinker/<br>302nm Transilluminator<br>Filter Size 20x20cm<br>Replacement Tube, 254nm<br>Replacement Tube, 302nm | 95-0300-01 (115V)<br>95-0300-02 (230V)<br><br>34-0007-01<br>34-0042-01 |

### TL-2000 Translinker

The TL-2000 Translinker components include the shortwave ultraviolet CL-1000 crosslinker and an M-20 transilluminator built onto the crosslinker cabinet.

- On/Off controls for the transilluminator are located on the Crosslinker panel
- The Transilluminator uses midrange 302nm UV

TL-2000 Dimensions:

13.7D x 19W x 8.75H in. (348 x 483 x 222mm)  
Weight: 36 lb (16.4 kg)