

EZ-Retriever® **IR System**

All-in-One | All-at-Once



SUPERIOR QUALITY



REAL FAST



REAL EZ

Features and Benefits:


- Infra-Red (IR) system maintains precise consistent processing temperature
- High boiling temperature EZ-AR elegance solution prevents boiling and spillovers
- Large touch screen with user friendly GUI for ease of protocol set-up
- Dewax, rehydration, antigen, nucleic acid and decal retrieval all in one
- High throughput standardized sample prep, 96 slides in 30 minutes
- Uses non-toxic, non-hazardous, and non-Xylene green solution
- One solution provides clean, crisp intense stain for all antibodies and Nucleic acids

Antigen retrieval technology pioneered by BioGenex disrupts formalin induced cross-linking bonds introduced during tissue processing. All FFPE (Formalin Fixed Paraffin Embedded), formalin fixed paraffin embedded tissues sections require retrieval process to make biomarker epitopes available for antibody binding or nucleic acid molecules for hybridization followed by visualization in IHC and ISH techniques.


Product	Xmatrx® Elite 200 Slides (RTU*)	Xmatrx® Infinity 200 Slides (RTU*)	Manual 1000 mL (RTU)
EZ-AR™ 1 Elegance	HX031-YCD	HX031-YCX	HK546-XAK
EZ-AR™ 2 Elegance	HX032-YCD	HX032-YCX	HK547-XAK

*80 µL/test for Xmatrx Elite, 70 µL/test for Xmatrx Infinity


Results with EZ-AR™ Elegance series

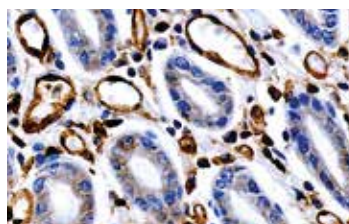


Gastrin

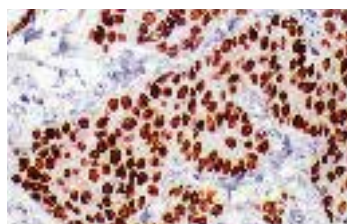


Estrogen Receptor

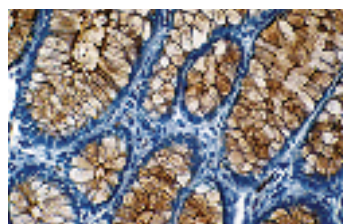


MUC4

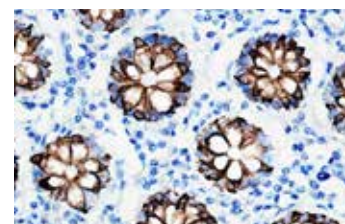


Cytokeratin 20

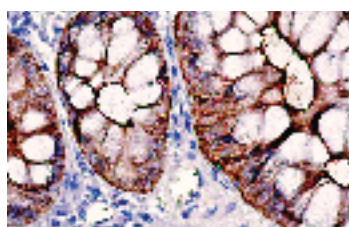


Villin

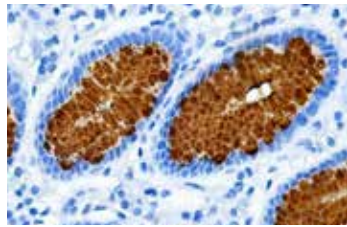

MUC5AC

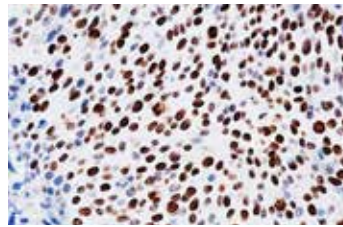

Ki67

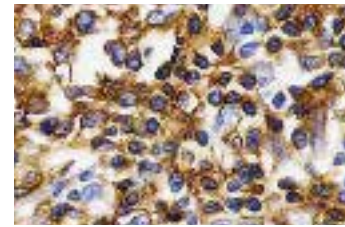

Alpha-1-Antitrypsin

PARAMETERS	EZ-RT (MW015-IR)	EZ-RT (MW016-IR)
Slide capacity	96 Slides	96 Slides
Oven capacity	30L	34L
Operating voltage	230 V AC, 50Hz	110 V AC, 60Hz
Microwave power	1500 W	1550 W
Dimensions (W/D/H)	24" X 16" X 12"	23" X 18" X 13"
Weight	46.3 lbs / 21kg	46.3 lbs / 21kg

USA : 48810 Kato Road, Suite 100E & 200E, Fremont, CA 94538 USA
 India : Aerospace & Precision Engineering, TSIC SEZ, Adibatla, Telangaa-501510
 China : Shanghai Exhibition Center, No.1000 Middle Yan An Road, Shang Hai Shi, China, 200041

Customer service - US : customerservice@biogenex.com
 Customer service - India : indiaccs@biogenex.com
 Customer service - Global : internationalcs@biogenex.com
 Instrument support : instrument.care@biogenex.com
 Reagent support : reagentsupport@biogenex.com

For US
 T: +1-510-824-1400
 T: +1-800-421-4149
 F: +1-510-824-1490

For India
 T: +91-6309448856

For China
 T: +86-21-2220-6680
 F: +86-21-2220-6679

For Other Countries
 T: +91-7702301043



www.biogenex.com