



Thermo Scientific
Heraeus Megafuge 8 Centrifuge Series

Fits in.
Stands out.

Thermo
SCIENTIFIC

Introducing the NEW!

Thermo Scientific™

Heraeus™ Megafuge™ 8 and 8R
small benchtop centrifuges



Fits in.



- ▶ to **your space**, with its compact footprint.
- ▶ to **your applications**, accommodating multi-laboratory settings with the flexibility to adapt with evolving clinical and research needs, in both ventilated and refrigerated options.

Heraeus Megafuge 8

- 580 mL capacity (4 x 145 mL)
- Simplified routine operating interface
- Performance drive
- 4-program memory
- 13 optional rotor systems

Heraeus Megafuge 8R


- 580 mL capacity (4 x 145 mL)
- Simplified routine operating interface
- High performance flexible shaft drive
- 99-program memory
- 15 optional rotor systems
- Powerful refrigeration system with pre-temperature function

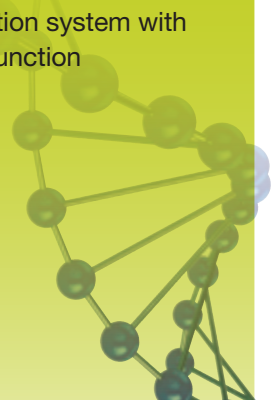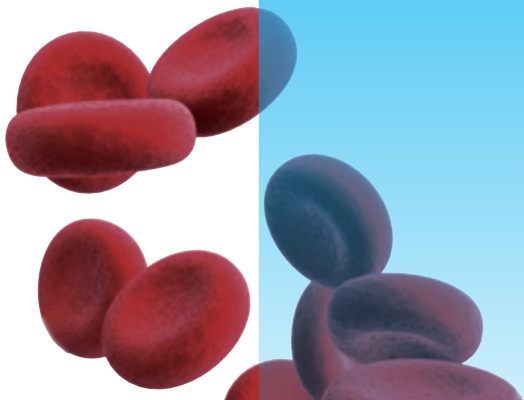



Stands out.

- ▶ designed for **solid reliability and consistent results.**
- ▶ featuring **innovative Thermo Scientific™ technologies:**
- ▶ anticipating and addressing your most pressing lab needs – even the ones you don't know you have yet.

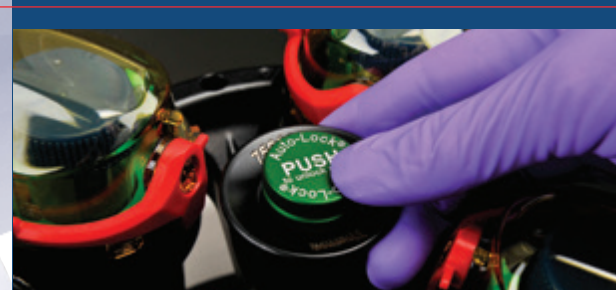

Auto-Lock™ Rotor Exchange

Secure push-button application versatility and cleaning convenience



ClickSeal™ Biocontainment Lids

One-handed, certified¹ sample protection



Maximized Swing-Out Capacity

Process up to 8 x 50 mL conical tubes, 24 x 5/7 mL blood tubes, microplates and microtubes, all in one versatile centrifuge



Fixed Angle Flexibility

Spin 50 mL conical tubes at over 12,000 x g or a wide range of microtubes at speeds up to 30,279 x g



Ventilated (left) and refrigerated (right)



¹ Biocontainment certification by the Public Health England, Porton Down, UK

Fits in. **Stands out.**



Simple solutions. **Flexible processing.**

Secure, push-button **Auto-Lock** rotor exchange in as little as 3 seconds delivers:

- Trouble-free rotor installation and removal
 - 1] No tools are required
 - 2] No complicated rotor changing procedures
- Clear chamber access for cleaning convenience
- Quick switch between applications – from conical tubes to blood tubes, and from microtubes to microplates – and the flexibility to evolve with the changing needs of your laboratory

Optimize bench space with **maximized capacity:**

- 8 x 50 mL swinging bucket rotor capacity or 4 microplates for research needs
- Micro-volume rotor portfolio supports 0.2, 0.25, 0.5, 1.5, 2 and 5 mL microtubes
- Clinical processing up to 24 x 5/7 mL blood tubes per run

Micro-volume Solution Rotors

Accommodates up to 30 spin columns, such as DNA and RNA mini-preps, including open lid designs



Reliable quality. **Safe operation.**

- Biocontainment sealing options, including certified¹ **ClickSeal** lids for glove-friendly, one-handed operation
- Manufactured with quality materials providing broad chemical resistance with exceptional strength for durability and reliability
- Low noise levels, including just 52 dBA at performance of 17,850 rpm/ 30,279 x g (with *Heraeus Megafuge 8R centrifuge*), to ensure a peaceful laboratory
- Conforms to the latest clinical and safety standards, such as UL, CE and IVD²

¹  Biocontainment certification by the Public Health England, Porton Down, UK

² See full list of standards & certifications on page 14



Sample-ready. User-friendly.

- One-touch operation with pre-saved protocols
- Highly visible backlit display for easy reading of parameters across the lab
- Optional indicators at end of run, including automatic lid opening, full flashing display and adjustable audible signal
- Glove- and detergent-friendly
- Multilingual instructions – English, Dutch, French, German, Italian, Russian and Spanish — on programming, run conditions, alerts and service messages



Fits in. **CLINICAL APPLICATIONS.**

Stands out.

**Outstanding capacity in a compact design
with a smart, simple interface.**




- Runs up to 24 x 5/7 mL blood tubes at a time in swinging bucket configuration
- Auto-Lock rotor exchange simplifies cleaning
- Conforms to the latest clinical and safety standards
- Fast acceleration and deceleration times for quicker spins

Expanded Programming Options

- **Alphanumeric program naming**, of up to 12 characters, promotes correct program selection every time.
- **Program-only mode** limits the control of the centrifuge to the set program and the control quadrant (Start, Stop, Pulse and Open), ideal for controlled environments



Rotors for Clinical Applications

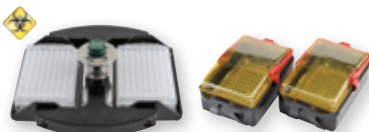Thermo Scientific Rotor	Description	Stand-out Features	Performance		
VENTILATED					REFRIGERATED
TX-150 Swinging Bucket Rotor					
	Designed for a broad range of high-speed, high-capacity swing-out processing needs, with a capacity of 24 x 5/7 mL blood tubes per run and flexibility with a wide range of adapters.	Certified ¹ sample protection with ClickSeal biocontainment lids.	4,500 rpm 3,260 x g	4,500 rpm 3,260 x g	
M10 Swinging Bucket Rotor					
	Centrifuge 4 microplates or 2 midi-deepwell plates, ideally suited for quick spins. Spin up to 4 standard microplates with optional certified ¹ ClickSeal biocontainment lids.	Transportable microplate carriers, complemented by Auto-Lock rotor exchange.	4,400 rpm 2,576 x g	4,400 rpm 2,576 x g	
8 x 50 mL Individually Sealed Fixed Angle Rotor					
	Designed for high-risk sample processing for a range of clinical and research tubes.	Individual sealed canisters for high-risk samples.	n/a	5,600 rpm 5,014 x g	
CLINIConic Fixed Angle Rotor					
	Spin 30 blood tubes or 15 mL conicals in a fixed angle configuration with fast acceleration to optimize blood separation applications.	Auto-Lock rotor exchange for ultimate ease-of-use.	4,400 rpm 3,030 x g	4,400 rpm 3,030 x g	
Hematocrit Rotor					
	Achieve red blood cell capillary separations.	Result reading chart provided.	13,300 rpm 16,810 x g	13,300 rpm 16,810 x g	

¹  Biocontainment certification by the Public Health England, Porton Down, UK

Dedicated Rotors for Routine Clinical Processing



TX-100S Clinical Swinging Bucket Rotor with Sealed Carriers

Spin 8 tubes with biocontainment sealing

Ventilated
4,500 rpm
3,260 x g

Refrigerated
4,500 rpm
3,260 x g



TX-100 Clinical Swinging Bucket Rotor with Carriers

Capacity of 16 tubes per run

Ventilated
4,500 rpm
3,260 x g

Refrigerated
4,500 rpm
3,260 x g

Fits in. **RESEARCH APPLICATIONS.**

Stands out.

Innovative design supports multiple users and applications to stay one step ahead of your evolving research needs.



- 8 x 50 mL swinging bucket rotor capacity, 4 microplates, or 30 filtration columns for research needs
- Auto-Lock rotor exchange for application flexibility – go from 50 mL tube to microplates to microtubes effortlessly



30,279 x g Performance


High-speed microliter performance at up to 30,279 x g enables

- quicker spins
- improved separations
- access to a wider range of applications

while maintaining full swing-out rotor flexibility, combined with quick, safe Auto-Lock rotor exchange.



Rotors for Research Applications

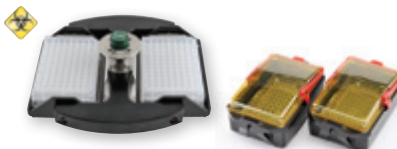Thermo Scientific Rotor		Description	Stand-out Features	Performance	
				VENTILATED	REFRIGERATED
TX-150 Swinging Bucket Rotor					
		Designed for a broad range of swing-out processing needs at high speed and high capacity with the flexibility of a wide range of adapters.	Dedicated bucket available for 8 conical tubes capacity or round buckets with certified ¹ ClickSeal biocontainment lids.	4,500 rpm 3,260 x g	4,500 rpm 3,260 x g
M10 Swinging Bucket Rotor					
		Centrifuge 4 standard microplates or 2 midi-deepwell plates, ideally suited for quick spins. Spin up to 4 standard microplates with optional certified ¹ ClickSeal biocontainment lids. Non-shielded design for quick access to microplates.	Transportable microplate carriers, complemented by Auto-Lock rotor exchange.	4,400 rpm 2,576 x g	4,400 rpm 2,576 x g
MT-12 Swinging Bucket Rotor					
		Design capable of spinning up to 12 x 1.5/2 mL microtubes and filtration columns, ideal for small-volume research applications.	Individual swinging bucket rotors for centralized pellets within tubes.	13,000 rpm 16,438 x g	13,000 rpm 16,438 x g
HIGHConic III Fixed Angle Rotor					
		Process 6 x 50 mL tubes in fixed angle configuration and at RCF above 10,000 x g in certified ¹ sealed conditions.	New transparent lid enables viewing for added safety and convenience, and Auto-Lock rotor exchange enhances ease-of-use.	8,700 rpm 10,155 x g	9,500 rpm 12,108 x g
CLINIConic Fixed Angle Rotor					
		Spin 30 x 15 mL round or conical tubes and ensure fast acceleration and throughput for cell culture applications.	Support conical tubes in a fixed angle rotor with the flexibility of Auto-Lock rotor exchange.	4,400 rpm 3,030 x g	4,400 rpm 3,030 x g
MicroClick 18 x 5 Rotor					
		18-place capacity for 5.0 mL microtubes for high throughput and over 22,000 x g to reduce run times.	Lightweight and ergonomic ClickSeal biocontainment lid allows sample viewing.	n/a	14,000 rpm 22,351 x g
Micro-Volume Solution Rotors					
		Capable of supporting micro-volume protocols, such as nucleic acid preparation, PCR reaction set-up and filtration columns, choose from the complete micro-volume rotor portfolio.	With Auto-Lock rotor exchange, switch from swing-out to a microtube rotor in as little as 3 seconds, coupled with security of certified ¹ ClickSeal biocontainment lids.	16,000 rpm 24,328 x g	17,850 rpm 30,279 x g
			Additionally, the MicroClick 30 x 2, is designed with an extended neck to accommodate 30 nucleic acid filter columns.	14,000 rpm 21,694 x g	14,000 rpm 21,694 x g
				15,000 rpm 17,860 x g	15,000 rpm 17,860 x g

¹  Biocontainment certification by the Public Health England, Porton Down, UK


Stands out.

Extensive Rotor
Portfolio for
Application Flexibility

Versatile Swinging Bucket Rotors



TX-150 Swinging
Bucket Rotor

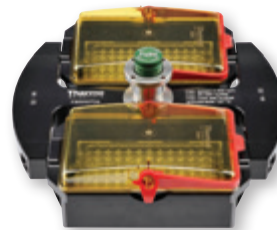

M10 Swinging
Bucket Rotor



MT-12 Swinging
Bucket Rotor



Cat. No.	Description	Rotor Capacity (places x volume, mL)	Max Tube Dimensions (O x L, mm)	Max Speed (rpm)	Max RCF (x g)	Max Speed (rpm)	Max RCF (x g)
Versatile Swinging Bucket Rotors				VENTILATED		REFRIGERATED	
75005701	TX-150 Swinging Bucket Rotor, 90°, R_{max} 144 mm (order buckets separately)	4 x 145	50 x 100	4,500	3,260	4,500	3,260
75005703	50 mL Conical Buckets (unsealed, no adapter needed) (set of 4)	8 x 50	29.5 x 120	4,500	3,260	4,500	3,260
75005702	Round Buckets (set of 4)	4 x 145	50 x 100	4,500	3,260	4,500	3,260
75005707	ClickSeal Biocontainment Lids for Round Buckets (set of 4)						
75005724	Replacement O-rings for Lids (set of 4)						
Adapters for 50 mL Conical Buckets (sets of 2)							
75005808	15 mL Conical Tube	8 x 15	17 x 123				
Adapters for Round Buckets (sets of 4)							
no adapter needed	145 mL Bottle (75005734)	4 x 145	50 x 100				
1 75005735	100 mL Round Bottom Open-Top Tube	4 x 100	45 x 117				
2 75005736	50 mL Conical or Skirted Tube	4 x 50	29.5 x 120				
3 75005744	25 mL Universal Flat/Skirted Tube	4 x 25	25 x 120				
4 75005737	15 mL Conical Tube	8 x 15	17 x 122				
4 75005737	11 mL IVF Tube	8 x 11	17 x 122				
5 75005738	14 mL Urine Tube	16 x 14	18 x 124				
5 75005738	12 mL Blood Collection Tube (Greiner™)	16 x 12	16 x 100				
5 75005738	10 mL Blood Collection or 15 mL Corex™/Kimble™ Tube	16 x 15	17 x 117				
6 75005739	5/7mL mL Blood Collection Tube	24 x 5/7	13 x 116				
7 75005740	3/5 mL Blood Collection Tube	28 x 3/5	12 x 124				
8 75005743	1.5/2 mL Microtube	40 x 2	11 x 65				
Rotor Packages							
75005760	Cell Culture Package, INCLUDES: TX-150 Rotor (75005701), Round Buckets (75005702), Adapters for 50 mL Conical Tubes (75005736)	4 x 50	29.5 x 120				
75005761	High Capacity Cell Culture Package, INCLUDES: TX-150 Rotor (75005701), Conical Buckets (75005703), Adapters for 15 mL Conical Tubes (75005808)	8 x 50	18 x 124				
75005762	Clinical Rotor Package, INCLUDES: TX-150 Rotor (75005701), Round Buckets (75005702), ClickSeal Biocontainment Lids (75005707), Adapters for Blood Collection Tubes: 5/7 mL (75005739) and 10 mL (75005738)	24 x 5/7	18 x 124				
75005706	M10 Swinging Bucket Rotor, 90°, R_{max} 119 mm (order buckets separately)		n/a	4,400	2,576	4,400	2,576
75005723	Unsealed Buckets (set of 2)	4 Standard or 2 Midi-Deepwell	Height < 33 mm				
75005721	Sealed Buckets (set of 2) R _{max} 109mm	4 Standard or 2 Midi-Deepwel	Height < 33 mm				
75005600	MT-12 Swinging Bucket Rotor, 90°, R_{max} 87 mm	12 x 2	11 x 50	13,000	16,438	13,000	16,438

Clinical Swinging Bucket Rotors



TX-100S Clinical Swinging
Bucket Rotor




TX-100 Clinical Swinging
Bucket Rotor

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Max Tube Dimensions (O x L, mm)	Max Speed (rpm)	Max RCF (x g)	Max Speed (rpm)	Max RCF (x g)
Clinical Swinging Bucket Rotors				VENTILATED		REFRIGERATED	
75005704	TX-100S Clinical Swinging Bucket Rotor with Sealed Carriers, 90°, Rmax 144 mm	8 x 15	17 x 106	4,500	3,260	4,500	3,260
75005705	TX-100 Clinical Swinging Bucket Rotor with Carriers, 90°, Rmax 144 mm	16 x 15	17 x 121	4,500	3,260	4,500	3,260

Adapters for TX-100S Clinical Rotor and TX-100 Clinical Rotor (each)

11203666	13.5 mL Urine Tube	16/8 x 13.5	16 x 82				
11172596	5/7 mL BD Hemogard™/ BD Vacutainer™ Tube	16/8 x 5/7	13 x 106				
11172595	5 mL BD Hemogard Tube	16/8 x 5	13 x 75				
11172287	3 mL Blood Collection Tube	16/8 x 3	11 x 70				
11172288	1.5/2 mL Microtube	16/8 x 1.5/2	10 x 41				

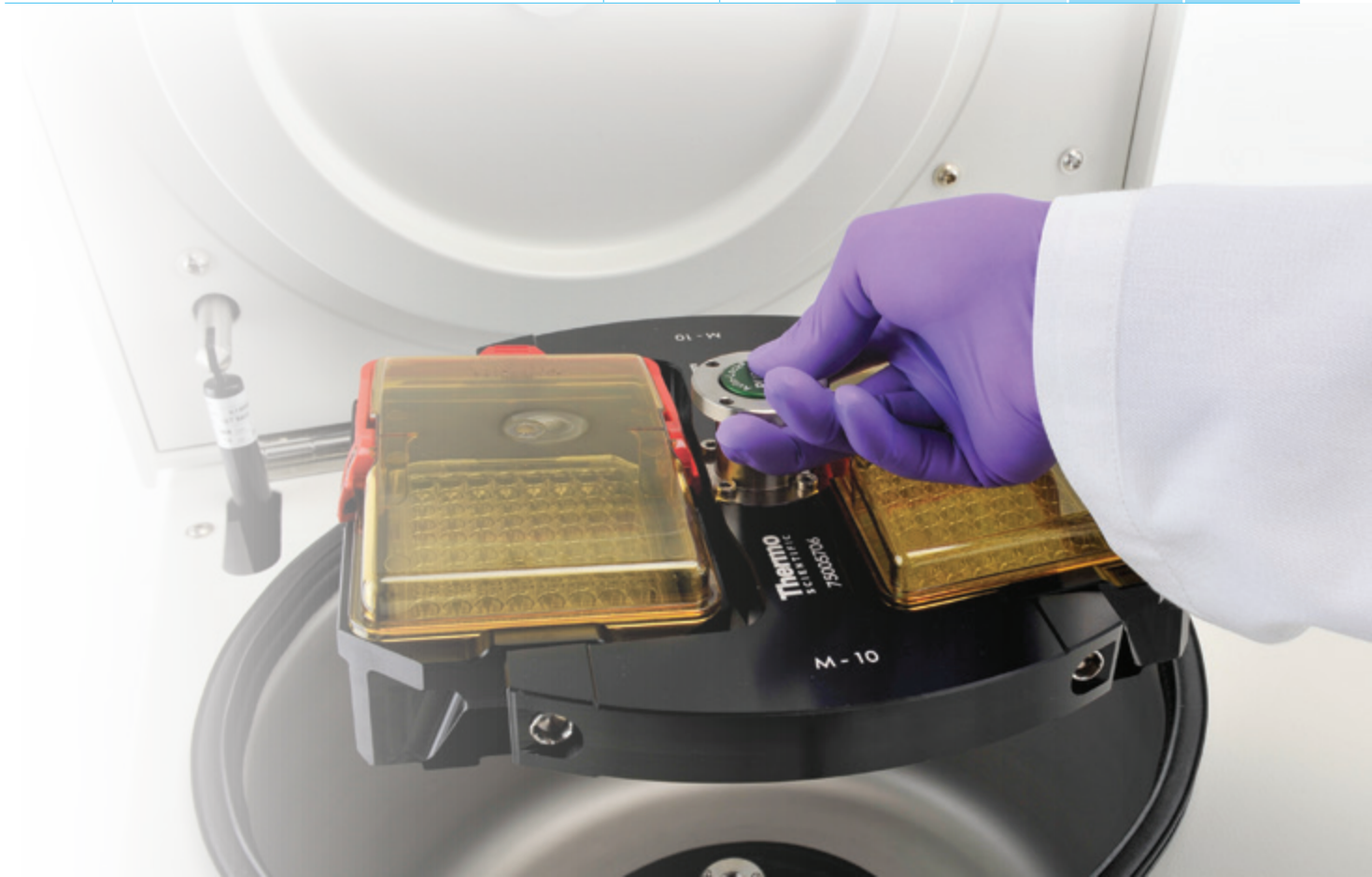

Fixed Angle Rotors



HIGHConic III Rotor



8 x 50 mL Rotor



CLINIConic Rotor



MicroClick 18 x 5 Rotor

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Max Tube Dimensions (O x L, mm)	Max Speed (rpm)	Max RCF (x g)	Max Speed (rpm)	Max RCF (x g)
Fixed Angle Rotors				VENTILATED		REFRIGERATED	
 75005709	HIGHConic III Rotor, 45°, R_{max} 120 mm	6 x 50	30 x 121	8,700	10,155	9,500	12,108
75005731	Replacement Lid (each)						
75003058	Replacement Sealing O-rings (set of 2 with grease)						
Adapters for HIGHConic III Rotor (set of 6)							
75005745	38 mL Round Bottom Tube	6 x 38	25.5 x 110				
75005746	16 mL Round Bottom Tube	6 x 16	18 x 123				
75005755	15 mL Conical Tube	6 x 15	17 x 123				
75005747	12 mL Round Bottom Tube	6 x 12	16 x 95				
75005748	6.5 mL Round Bottom Tube	6 x 6.5	13.5 x 114				
75005770	5 mL Conical Microtube	6 x 5	17 x 100				
75005749	3.5 mL Round Bottom Tube	12 x 3.5	11 x 100				
75005750	1.5/2 mL Microtube	12 x 2	11 x 40				
 75003694	8 x 50 mL Individually Sealed Rotor 25°, R_{max} 143 mm	8 x 50 Sealed	30 x 121	n/a	n/a	5,600	5,014
75003011	Replacement Biocontainment Lids (set of 2)						
75003689	Replacement O-ring Kit						
Adapters for 8 x 50 mL Individually Sealed Rotor (set of 8)							
75005755	15 mL Conical Tube	8 x 50	17 x 121				
75005747	10 mL Blood Collection Tube (16 x 100 mm)	8 x 10	16 x 110				
75005748	7 mL Blood Collection Tube (13 x 100 mm)	8 x 7	13 x 110				
75005749	3.5 mL Blood Collection Tube	16 x 3.5	11 x 75				
75003623	CLINIConic Rotor, 37°, R_{max} 140 mm	30 x 15	16.5 x 131	4,400	3,030	4,400	3,030
Adapters for CLINIConic Rotor (each)							
no adapter needed	15 mL Round/Conical Bottom Tube	30 x 15	16.5 x 131				
75003702	10 mL Round Bottom Tube	30 x 10	16.5 x 95				
11172596	5/7 mL BD Hemogard/ BD Vacutainer Tube	30 x 5/7	13 x 106				
11172595	5 mL BD Hemogard Tube	30 x 5	13 x 75				
 75005765	MicroClick 18 x 5 Rotor, R_{max} 102 mm	18 x 5	17 x 62	n/a	n/a	14,000	22,351
Adapters for MicroClick 18 x 5 (set of 2)							
75005756	1.5/2 mL Microtube	18 x 1.5/2	11 x 45				

Fixed Angle Rotors



MicroClick 30 x 2 Rotor



MicroClick 24 x 2 Rotor



8 x 8 PCR Strip Rotor



Hematocrit Rotor

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Max Tube Dimensions (O x L, mm)	Max Speed (rpm)	Max RCF (x g)	Max Speed (rpm)	Max RCF (x g)
Fixed Angle Rotors				VENTILATED		REFRIGERATED	
 75005719	MicroClick 30 x 2 Rotor, 45°, R_{max} 99 mm	30 x 1.5/2	11 x 45	14,000	21,694	14,000	21,694
75005730	Replacement ClickSeal Biocontainment Lid (each)						
75005726	Replacement O-ring for Lid (each)						
 75005715	MicroClick 24 x 2 Rotor, 45°, R_{max} 85 mm	24 x 1.5/2	11 x 45	16,000	24,328	17,850	30,279
75005725	Replacement ClickSeal Biocontainment Lid (each)						
75003405	Replacement O-ring for Lid (each)						
Adapters for MicroClick 30 x 2 Rotor and MicroClick 24 x 2 Rotor (sets of 30)							
75005752	0.2 mL PCR Tube	30/24 x 0.2	6.5 x 20				
75005753	0.5 mL Microtube	30/24 x 0.5	8 x 44				
75005754	0.25 mL Microtube	30/24 x 0.25	6 x 46				
 75005720	8 x 8 PCR Strip Rotor, 45°, R_{max} 71 mm	8 x 8 PCR strips	8 x 8	15,000	17,860	15,000	17,860
75005730	Replacement ClickSeal Biocontainment Lid (each)						
75005726	Replacement O-ring for Lid (each)						
75005733	Hematocrit Rotor, 90°, R_{max} 85 mm, includes Linear Scale	24 Capillaries	n/a	13,300	16,810	13,300	16,810
76000923	Capillaries (pack of 100)						
75000964	Sealing Kit						
76000938	Replacement Linear Scale Reading Graph						



Thermo Scientific Heraeus Megafuge 8 Centrifuge Series



Specifications	Heraeus Megafuge 8	Heraeus Megafuge 8R
Swing Bucket Rotors		
Max Capacity	4 x 145 mL (with TX-150 rotor)	4 x 145 mL (with TX-150 rotor)
Max Speed	4,500 rpm (with TX-150 rotor)	4,500 rpm (with TX-150 rotor)
Min Speed	300 rpm (with TX-150 rotor)	300 rpm (with TX-150 rotor)
Max RCF	3,260 x g (with TX-150 rotor)	3,260 x g (with TX-150 rotor)
Fixed Angle Rotors		
Max Capacity	6 x 50 mL (with HIGHConic rotor)	6 x 50 mL (with HIGHConic rotor)
Max Speed	16,000 rpm (with MicroClick 24 rotor)	17,850 rpm (with MicroClick 24 rotor)
Min Speed	300 rpm (with MicroClick 24 rotor)	300 rpm (with MicroClick 24 rotor)
Max RCF	24,328 x g (with MicroClick 24 rotor)	30,279 x g (with MicroClick 24 rotor)
Control System		
Microprocessor		
Drive System		
Direct, brushless induction drive		
Rotor Locking System		
Auto-Lock rotor exchange		
Imbalance Detection System		
Continuous vibration measurement, with rotor mass correction		
Programs		
4 direct program buttons, with alphanumeric program naming		3 direct program buttons, plus 96 additional programs accessible via folder, all with alphanumeric program naming
Pulse (Short) Run		
Yes		
Acceleration / Deceleration Rates		
2 (standard and soft)		
Centrifugation Chamber		
Brushed stainless steel		
Max Timer Range		
99h 59min, continuous		
Temperature Range		
n/a		-10°C to 40°C
Pre-Cooling Function		
n/a		Yes, with direct button
Refrigeration System		
n/a		CFC Free
Sound Level		
<58 dBA (with TX-150 rotor); <61 dBA (with MicroClick 24 rotor)		<56 dBA with all rotors; <52 dBA (with MicroClick 24 rotor)
Other Features		
Certified ¹ ClickSeal biocontainment lids Multilingual selection: English, Dutch, French, German, Italian, Russian and Spanish		
Max Power Consumption		
310 W		750 W
Dimensions		
Height (lid open)	67 cm (26.8 in)	70 cm (27.5 in)
Height (lid closed)	31 cm (12.2 in)	32 cm (12.6 in)
Width	37 cm (14.6 in)	46 cm (18.1 in)
Depth	48 cm (18.9 in)	67 cm (26.4 in)
Weight		
35 kg (77 lb)		71 kg (156 lb)
Standards & Certifications		
Standards	IEC 61010-1, IEC 61010-2-020, IEC 61010-2-101, EN 61326-1	
Certifications	UL Listed / CE marked / IVD compliant / Certified biocontainment ¹	

Ordering Information

75007211	Heraeus Megafuge 8 centrifuge, ventilated, 230 V, 50/60 Hz
75007210	Heraeus Megafuge 8 centrifuge, ventilated, 120 V, 60 Hz
75007214	Heraeus Megafuge 8R centrifuge, refrigerated, 230 V, 50/60 Hz
75007213	Heraeus Megafuge 8R centrifuge, refrigerated, 120 V, 60 Hz

¹ Biocontainment certification by the Public Health England, Porton Down, UK



Fits in.
Stands out.

Fits in. More centrifuge solutions for your lab

Perfect Fit. It's simple. From 1 Liter bottles, to 15 and 50 mL conical tubes and microplates, the versatile selection of **Thermo Scientific™ Nalgene™ and Nunc™ centrifugation products** work seamlessly with your complete centrifuge and rotor system, bringing together outstanding quality and performance.



Thermo Scientific Heraeus Benchtop Centrifuges

Accelerate your research with high speeds, rotor versatility and outstanding safety. From intuitive controls, to the certified¹ Thermo Scientific ClickSeal biocontainment rotor lids that “click” in place, Thermo Scientific Heraeus benchtop centrifuges deliver the perfect blend of capabilities to support micro-volume protocols, a wide range of high volume swing-out applications and high speed processing.

Thermo Scientific™ Heraeus™ Pico™ and Fresco™ microcentrifuge series

Select from two available speed options, up to 17,000 x g for routine laboratory work, or 21,000 x g for high-end research applications, and our extensive rotor portfolio, including:

- Standard 24-place rotor supports 1.5 to 2.0 mL tubes and spin columns.
- Unique dual row rotor holds 18 x 1.5/2.0 mL or 18 x 0.5 mL tubes, eliminating the need for adapters.

Thermo Scientific™ Heraeus™ Multifuge™ general purpose centrifuges

Improve performance and endurance with innovative technologies with 1.6 or 4 Liter capacity in a compact footprint.

- Support multiple users with the flexibility to adapt with evolving clinical and research needs, up to 196 x 5/7 mL blood tubes or 40 x 50 mL conical tubes in one spin.
- Lightweight Thermo Scientific™ Fiberlite™ carbon fiber rotors provide outstanding ergonomics and performance.

¹  Biocontainment certification by the Public Health England, Porton Down, UK



thermoscientific.com/smallbenchcentrifuge

© 2014 Thermo Fisher Scientific Inc. All rights reserved.

HERAEUS is a registered trademark of Heraeus Holding GmbH licensed to Thermo Fisher Scientific. Greiner is a registered trademark of Greiner. Corex is a registered trademark of Corning Inc. Kimble is a registered trademark of Kimble Glass. BD Hemogard and BD Vacutainer are registered trademarks of BD Biosciences. All other trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

Australia +61 39757 4300
Austria +43 1 801 40 0
Belgium +32 53 73 42 41
China +800 810 5118 or
+400 650 5118
France +33 2 2803 2180
Germany national toll free 0800 1 536 376
Germany international +49 6184 90 6000

India toll free 1800 22 8374
India +91 22 6716 2200
Italy +39 02 95059 552
Japan +81 3 5826 1616
Netherlands +31 76 579 55 55
New Zealand +64 9 980 6700
Nordic/Baltic/CIS countries
+358 9 3291 0200

Russia +7 812 703 42 15
Spain/Portugal +34 93 223 09 18
Switzerland +41 44 454 12 22
UK/Ireland +44 870 609 9203
USA/Canada +1 866 984 3766

Other Asian countries +852 2885 4613
Countries not listed +49 6184 90 6000

Thermo
SCIENTIFIC
Part of Thermo Fisher Scientific