
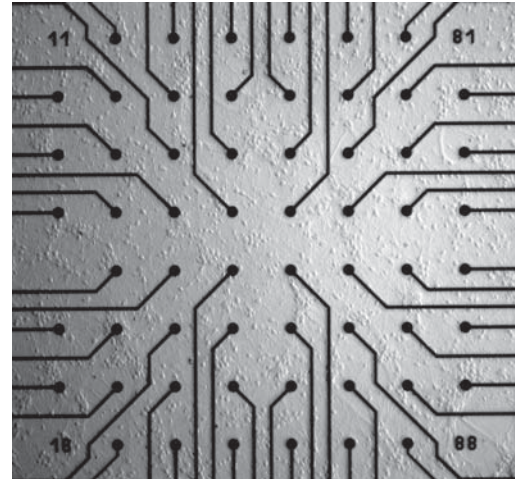


# Standard 60MEA

60MEA200/10iR-ITO, 60MEA200/30iR-ITO,  
60MEA100/10iR-ITO, 60MEA100/10-ITO,  
60MEA100/10iR-Ti, 60MEA200/30-Ti,  
60MEA200/10iR-Ti, 60MEA200/30iR-Ti

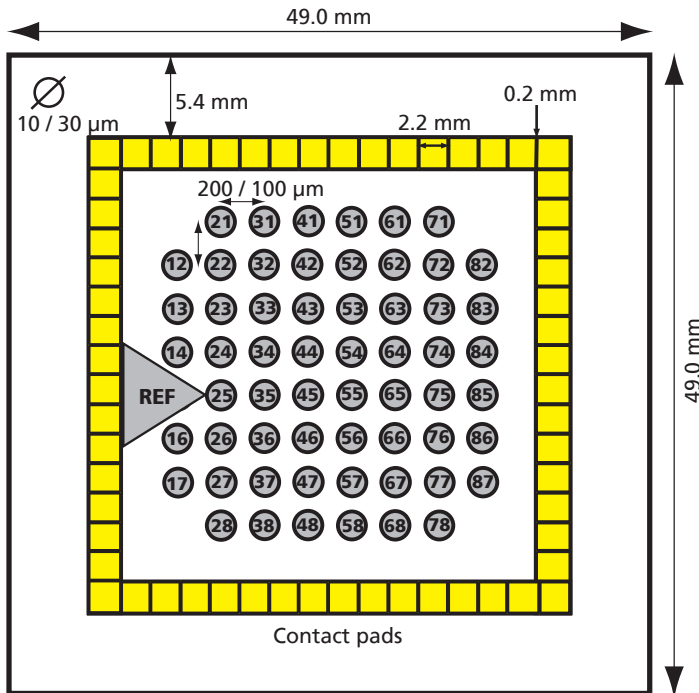

Standard electrode layout grid 8 x 8

## Technical Specifications Standard 60MEA

Temperature compatibility

Dimension (W x D x H)

Base material

Track material

Contact pads

Electrode diameter

Interelectrode distance (centre to centre)

Electrode height

Electrode type

Isolation type

Electrode impedance

Electrode layout grid

Number of recording electrodes

Number of reference electrodes

0 - 125 °C

49 mm x 49 mm x 1 mm

Glass

ITO (Indium tin oxide) or Ti (Titanium)

ITO (Indium tin oxide) or TiN (Titanium nitride)

10 or 30 µm

100 or 200 µm

Planar

TiN (Titanium nitride)

Silicon nitride 500 nm (PEVCD)

30 - 50 kΩ for 30 µm electrodes, 250 - 400 kΩ for 10 µm electrodes

8 x 8

59 (with iR) or 60 (without iR)

1 internal reference electrode (iR) or without internal reference

MC\_Rack

Source layout in "Data Source Setup"

Channel map

2 dim. (MEA) or Configuration

Default

MEA perfusion chamber

(w/o)

(gr)

(pr)

(pr-T)

Without ring


Glass ring: ID +/- 19 mm, OD 24 mm, height 6 / 12 mm

Plastic ring without thread: ID 26.5 mm, OD 30 mm, height 6 / 3 mm

Plastic ring with thread: ID 26 mm, OD 30 mm, height 6 / 15 mm

# Standard 60MEA

## Standard electrode layout grid 8 x 8

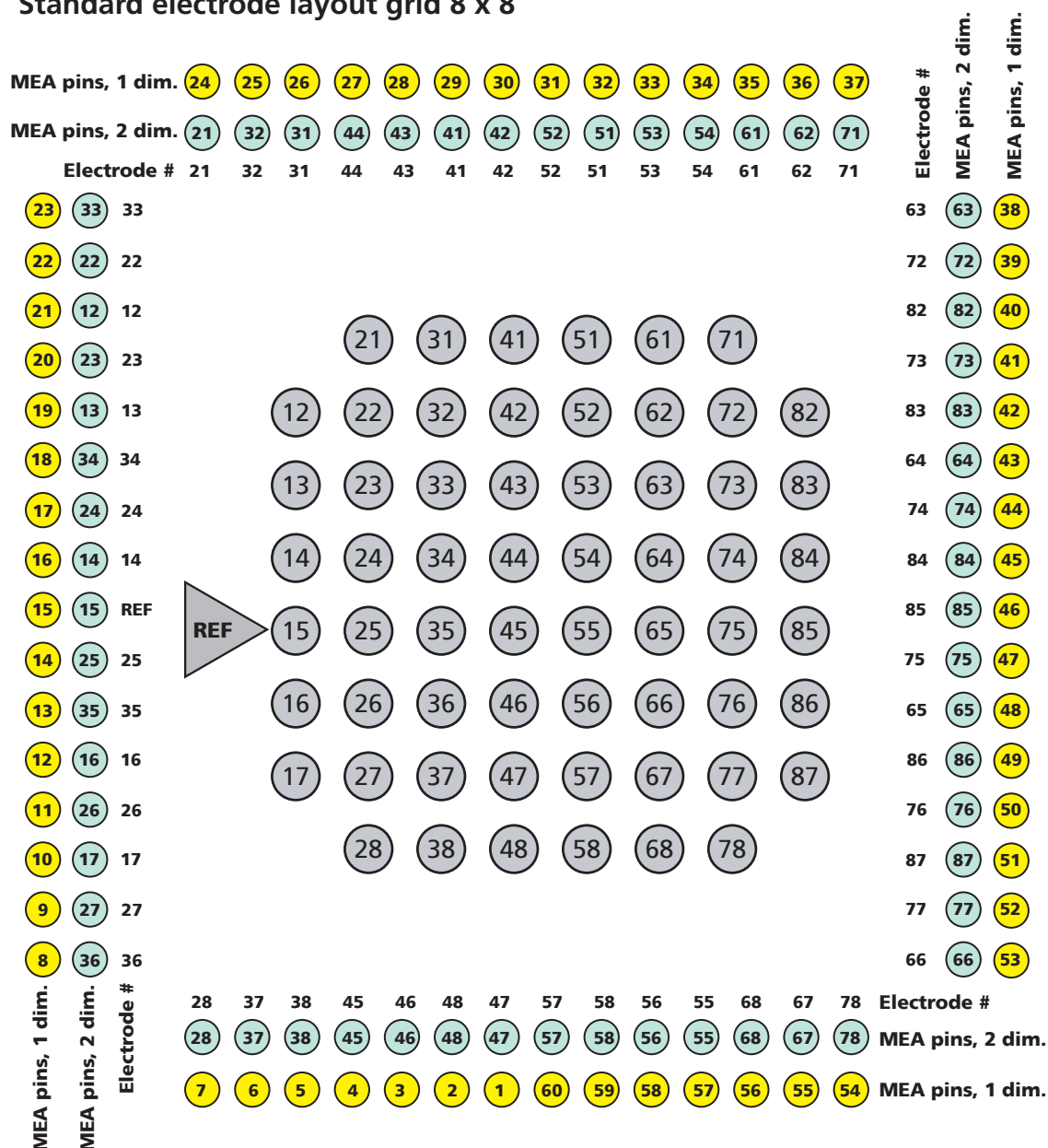

The numbering of MEA electrodes in the 8 x 8 grid follows the standard numbering scheme for square grids: The first digit is the column number, and the second digit is the row number. For example, electrode 23 is positioned in the third row of the second column.

The specified MEA pin numbers (1 dim. or 2 dim.) are the channel numbers that are used in the data acquisition program, when using the 1 dimensional layout or the 2 dimensional layout (or Configuration) in the "Data Source Setup". The electrode 15 is missing in MEAs with internal reference electrode. It is replaced by a big internal reference electrode, connected to pin 15 of the amplifier.