

Bioneer

Bioneer Corporation
is Korea's leading
biotech company.

Exicycler™ 96

Superior 5-color Real-Time Quantitative PCR System

Real-Time Quantitative PCR System
Superlative optics for superior results




CE-IVD

8-11 Munpyeongseo-ro, Daedeok-gu, Daejeon 306-220 Republic of Korea
Tel: (Korea) 1588-9788 (International) +82-42-930-8777 Email: sales@bioneer.com


BIONEER
Innovation • Value • Discovery

Superior Optic Module

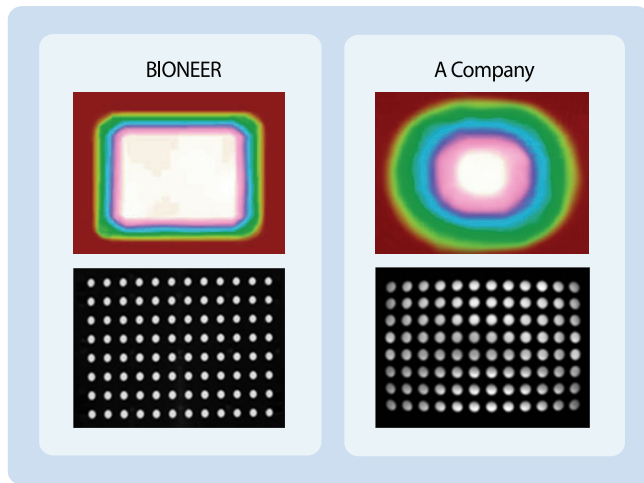




- Homogenous illumination with our exclusive Light Tunnel Technology
- Eliminate intensity differences between wells
- No need to normalize by using a reference dye
- All 5 channels can be used for actual experiments

Advanced Optical Module and Detection System

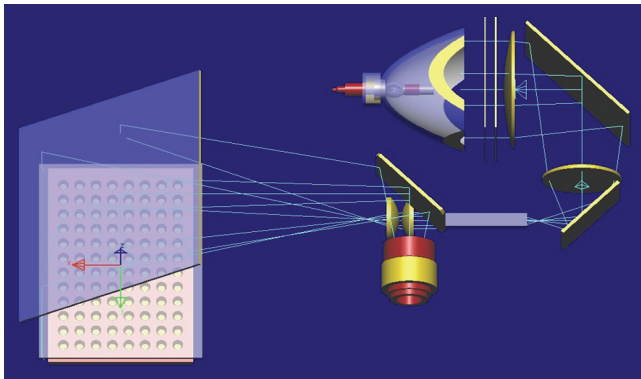


- Obtain fluorescence data from all 96 wells at once with our highly sensitive 2-dimensional CCD
- No time-lag between wells
- Bright white light from our Short Arc lamp provides uniform intensities for all colors
- Light Path Mask eliminates non-well light contribution
- Much longer lamp life relative to Halogen bulbs

Reproducible Cts & Wide Linear Dynamic Range

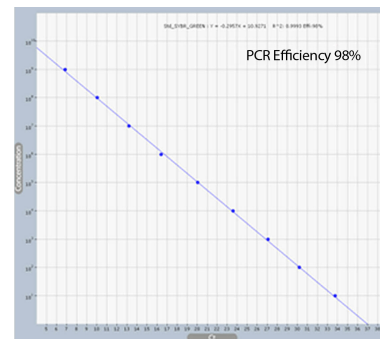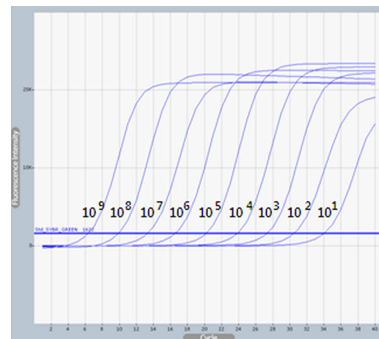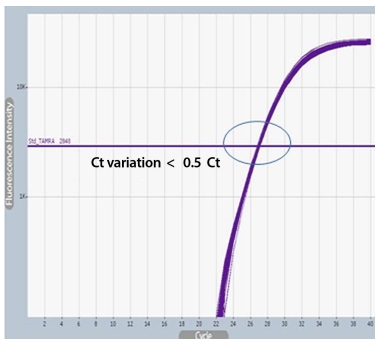

Real 5-color Multiplex PCR

Filter	Excitation	Emission	Fluorescence dye
1	490 nm	520 nm	FAM, SYBR Green I
2	520 nm	550 nm	JOE, TET
3	550 nm	580 nm	TAMRA, Cy3
4	580 nm	610 nm	Texas Red, ROX, Red610
5	630 nm	680 nm	Cy5, Red670

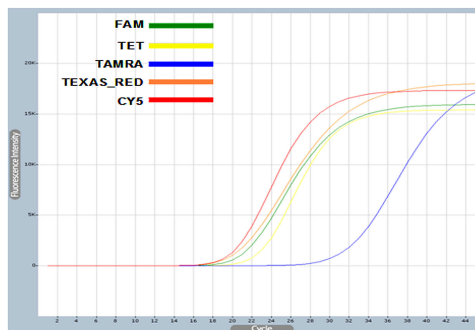

- No need to reserve one channel for reference dye, 5-color multiplexing is available
- Filter 4 is freed for use with Texas Red, Red610 or ROX
- Individual filters for each excitation wavelength results in maximum fluorescence for each dye
- Eliminate the fluorescence overlap between dyes when designing multiplex experiments
- Use fluorescence dyes with long excitation wavelength such as Cy5

Powerful Data Acquisition and Processing Algorithms

Artifacts are removed from the raw fluorescence data by:

- Well Quantitation Algorithm
- Fluorescence Intensity Normalization Algorithm
- Background Subtraction Algorithm
- Cross-Talk Compensation Algorithm

All the basic data analysis steps are automatically done by the following carefully designed algorithms:


- **Amplification Success/Fail Decision Algorithm**
Corrects errors when any well with an amplified sample is mistakenly assigned as empty
- **Baseline Decision Algorithm**
Determines the proper baseline for any possible type of amplification plot
- **Threshold Decision Algorithm**
Determines the proper Ct value using the latest methods

Final results of data analysis are obtained by the software with the following characteristics:

- Individual probe based analysis providing flexibility in experimental design
- Statistically sound and automatic decision for all analysis modules
- Core parameters and options are user adjustable for the fine tuning of analysis results
- A flexible area detection method for Melting Curve Analysis (SYBR Green I)

Full Featured Analysis Modules with User-Friendly GUI

Absolute Quantification

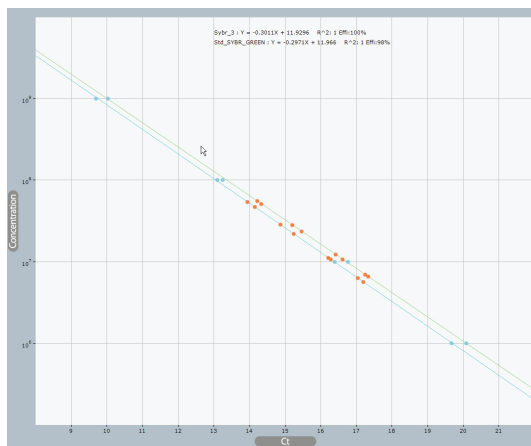


The Ct value of samples are plotted on a standard curve, allowing absolute quantification of unknown samples.

Relative Quantification




The relative expression levels of a target gene can be compared among different samples.

Existence/Nonexistence Assay

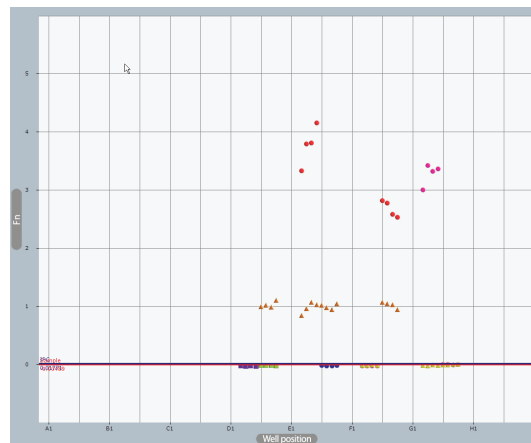

The existence and nonexistence of, for example, pathogenic viruses can be determined.


SNP Genotyping

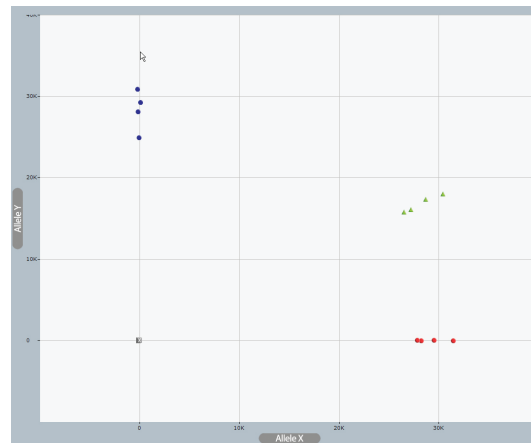

Homozygous and heterozygous alleles can be determined for a SNP site.



Well-to-well Uniformity

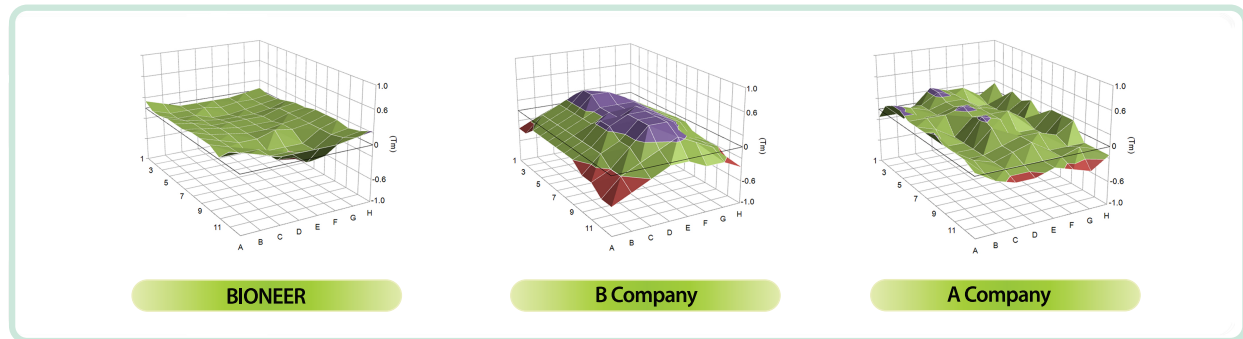

Other Convenient and Powerful Function

- **Self Diagnosis**
Detects hardware and software issues through a self-diagnosis protocol each time the system is turned on
- **Motorized Loading Tray**
Supports automation using robotics for large scale experiments
- **Post-Run Assignment**
The data from all 96 wells are always acquired and kept
- **Thermal Gradient**
For ease of experimental optimization of qPCR experiments
- **Time and Temperature Increment**
- **Uniform Block Heating**
Heat sink with thermal-tunnel guarantees uniform block heating
- **Standard Format Consumables**
- **No moving parts except for the loading tray and filter wheel**
Reliable, quiet operation and low maintenance

Superlative Optics for Superior Results

Related products:

- Reagents
- Plasticware
- Primer and probe

*Bioneer manufactures high quality consumables optimized for Exicycler™ 96.



Specifications

System specifications	
Dimension (mm)	355(W) x 540(D) x 470(H)
Weight (kg)	30 kg
Sample capacity size	96-well plate / 0.2 ml micro tubes
Sample volume	50 μ l recommended
Power consumption	100 ~ 240 VAC, 50 / 60 Hz 850 Watts
Operating temperature	15 ~ 30 °C
Operating humidity	20 ~ 80 %, no condensation
Thermo module specifications	
Method of heating / cooling	Peltier
Temperature range	4.0 °C ~ 99.9 °C
Temperature accuracy	± 0.3 °C
Temperature uniformity	± 0.5 °C
Heating and cooling rate	Max. 2.5 °C/sec
Temperature increment range	0.1 °C ~ 9.9 °C
Time increment range	1 sec ~ 60 sec
Computer specifications	
Operating system	Windows XP & Windows 7 (32 - bit OS only, S/W version 3.54.4 or later)
Processor speed	Intel Dual Core E2160 (1.8 GHz) or higher
Memory	1GB or higher
Communication port	USB 2.0 high speed
Screen resolution	1280 x 1024 or higher
Optics specifications	
Light source	Short arc lamp (120W)
Lamp life time	3,000 hours
Sensor	16 - Bit 2D CCD
Excitation filter / Emission filter	5 Sets

Ordering Information

Cat. No.	Product Description
A-2060	Exicycler™ 96 Real-Time Quantitative PCR System
K-6100	AccuPower® DualStar™ qPCR PreMix, 20 ul reaction, 12 x 8-strip tubes (96 rxns), optical sealing film included
K-6103	AccuPower® DualStar™ qPCR PreMix, 20 ul reaction, 1 x 96-well plate, optical sealing film included
K-6110	AccuPower® DualStar™ qPCR PreMix, 50 ul reaction, 12 x 8-strip tubes (96 rxns), optical sealing film included
K-6113	AccuPower® DualStar™ qPCR PreMix, 50 ul reaction, 1 x 96-well plate, optical sealing film included
3111-52	Opaque white 96-well semi-skirted PCR plate for Real-Time PCR, 25 plates
3111-50	Opaque white 0.2 ml 8-strip PCR tubes for Real-Time PCR, 250 strips
3111-41	Adhesive Optical Sealing Film for Real-Time PCR, 100 sheets

The specifications of this product can be changed without notice.

Contact Us

Bioneer Corporation 8-11, Munpyeongseo-ro Daedeok-gu, Daejeon 306-220 Republic of Korea Tel: (Korea) 1588 - 9788 (International) + 82 - 42 - 930 - 8777 Fax: + 82 - 42 - 930 - 8600 E-mail: sales@bioneer.com	Bioneer, Inc. 1000 Atlantic Avenue Alameda, CA 94501 USA Toll free: + 1 - 877 - 264 - 4300 Fax: + 1 - 510 - 865 - 0350 E-mail: order.usa@bioneer.us.com	Bioneer Trade(Shanghai) Co., Ltd 403 Room, Building 88, no.887 Zuchongzhi Road, Zhangjiang High Tech Park PuDong New District, Shanghai 201203 China Tel: + 86 - 21 - 5080 - 0969 Fax: + 86 - 21 - 5080 - 1620 E-mail: salescn@bioneer.com	Bioneer Pangyo R&D Center B-702, Korea Bio Park Bldg., 700, Daewangpangyo-ro, Bundang-gu, Seongnam-si, Gyeonggi-do, 463-400, Republic of Korea Tel: + 82 - 31 - 628 - 0500 Fax: + 82 - 31 - 628 - 0555
---	---	---	---