

P-2000 Digital Polarimeter

Ultimate flexibility for a wide range of applications

*A new generation of polarimeters
offering simplicity, reliability and
complete validation protocols*



Jasco

P-2000 Digital Polarimeter

Ultimate flexibility for a wide range of applications ranging from quality control to research and technology



Foods

Pharmaceuticals

Sugars and Sweeteners


Essential oils, Flavors and Fragrances

Chemicals

JASCO, a leading manufacturer of polarimeters since 1967, is proud to introduce the new P-2000 multi-option polarimeter. The P-2000 is designed as a customizable polarimeter with various options for a range of applications and budgetary requirements. The instrument system can also be field upgraded as the application requirements change.

Advanced instrument control

Two graphical user interfaces are available including a newly redesigned intelligent remote module (iRM-800) with a color LCD touch screen and Spectra Manager™ II software, the latest version of JASCO's innovative cross-platform spectroscopy software. Both of these control and analysis interfaces allow full-system control and advanced data processing. 21 CFR part 11 compliant versions for both interfaces can be ordered as an option.

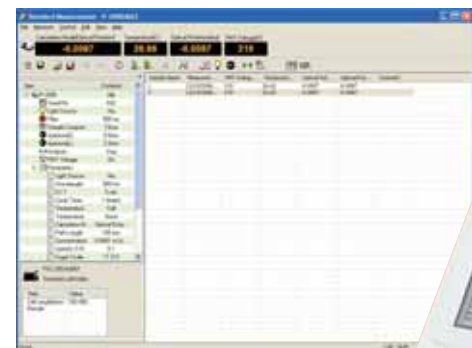

Spectra Manager™ II



iRM-800

The iRM conveniently guides the operator through routines from data acquisition to data processing. The obtained data can be automatically printed to USB printers, or saved to a compact flash memory card for further processing on a PC.

Full line of accessories

In support of diverse applications, a full line of accessories including a variety of cylindrical cells: a demountable stainless steel cell, Peltier thermostatted cell holder, and monochromator options, all offered for functional expansions of the P-2000. A newly developed Peltier cell holder offers accurate temperature control with $\pm 0.1^\circ$ accuracy.

Interference filter



A wide range of interference filters from UV-Vis to NIR

8-position automatic filter changer is optionally available.

Automatic light source switching mirror option

Light source

Up to two light sources can be installed. Available light sources are;

WI (Tungsten-Halogen lamp)
Na (Sodium lamp)
Hg (Mercury lamp)

RSC-200 V-shape cell holder

Temperature control by an optional thermostatted water circulation bath and a water-jacketed cylindrical cell



PTC-203 Peltier thermostatted cell holder

Temperature control by Peltier effect with a SUS demountable cell or a cylindrical cell



SHP-201P Peltier sipper SHP-201 Sipper

For quick measurements of multiple samples



High speed, precision and accuracy

The P-2000 offers a response speed as fast as six degrees per second to provide reproducible data for each measurement with a resolution of 0.0001°. A wide dynamic range of up to $\pm 90^\circ$ enables the system to measure chiral compounds over a broad range of concentrations. The sample chamber is equipped with two integrated temperature probes for cell holder and sample.

Reliability and instrument verification

An automated validation program is available for GLP/cGMP laboratories or those regulated by the FDA. The program can be used for periodic validation of light source energy, zero repeatability, and rotational accuracy/repeatability. Inspected filters, sample cells with certified optical path length, and NIST-traceable rotation plates are also available.

Ultimate flexibility

With its single optical platform, the P-2000 can be customized to provide the optimum performance for every application. By selecting the most suitable combination of optical elements, the instrument can be used in the ultra-violet, visible and/or NIR regions.

Examples of versatile system configurations

A single wavelength system

546 nm with Hg or WI lamp

589 nm with Na or WI lamp

A dual wavelength system

589/546 nm with Na/Hg or WI lamp

5 wavelength system

589/578/546/436/365 nm with Na/Hg lamps

6 wavelength system

633/589/546/436/405/365 nm with WI lamp

7 wavelength system

633/589/578/546/436/405/365 nm with WI lamp

System expansion to UV range

334/313/302/296/280/254 nm with Hg lamp with the optional UV/NIR region wavelength extension unit

334/325 nm with WI lamp with the optional UV/NIR region wavelength extension unit

Multiple wavelength system for UV/Vis range

633/589/578/546/436/405/365/334/325 nm with WI lamp with the optional UV/NIR region wavelength extension unit

System expansion to NIR range

880 nm with WI lamp with the optional UV/NIR region wavelength extension unit

Multiple wavelength system for NIR and Visible range

880/589 nm with WI lamp with the optional UV/NIR region wavelength extension unit

880/589/546 nm with WI lamp with the optional UV/NIR region wavelength extension unit

Hardware Specifications

Principle	Automatic digital polarimeter with symmetric angular oscillation using the optical-null balance method
Light source	Tungsten-Halogen lamp (WI), Sodium lamp (Na), Mercury lamp (Hg) (Up to two light sources can be installed.)
Wavelength	880, 633, 589, 578, 546, 436, 405, 365, 334, 325, 313, 302, 296, 280, 254 nm
Aperture	1.8, 3 and 8 mm diameters
Angular range	±90°
Response speed	6"/sec
Measurement accuracy	±0.002° (up to 1°), ±0.2% (larger than 1°)
Repeatability	0.002°
Resolution	0.0001°
Integration time	1 - 100 sec
Detector	Photomultiplier tube
Automatic recognition	Light source, filter, cell holder
Temperature measurement	Cell holder, cell (sample)
Temperature measurement range	0~40°C (Temperature accuracy: ±0.1°C)
External output	Analog output
Dimensions	653 (W) x 249 (D) x 364 (H) mm
Weight	Approx. 30 kg
Power requirement	AC100~240V ±10%, 50 or 60 Hz, 80 - 300 W

iRM-800 Intelligent Remote Module

Display	320 x 240 pixel color LCD Touch sensitive screen
Measurement modes	Polarimeter measurement Readout modes: Optical rotation, optical specific rotation, concentration, sugar scale Z, Brix purity, optical purity Statistical calculation: Average value, standard deviation, coefficient of variation Time course measurement Readout modes: Optical rotation, optical specific rotation Data processing: Reaction rate calculation
Other standard function	Instrument validation
Data handling	Data format: JASCO format, text format
Data save	Compact flash memory card
Optional	CFR version Printer: USB compatible printer

Spectra Manager II

Measurement modes	Polarimeter measurement Readout modes: Optical rotation, optical specific rotation, concentration, sugar scale Z, Brix purity, optical purity Statistical calculation: Average value, standard deviation, coefficient of variation Time course measurement Readout modes: Optical rotation, optical specific rotation Data processing: Reaction rate calculation
Other standard function	Instrument validation
Optional	CFR version

Lamps and filters

Filter	Light Source		
	WI (Tungsten-Halogen lamp)	Na (Sodium lamp)	Hg (Mercury lamp)
880 mm	UV/NIR region wavelength extension unit		
633 mm			
589 mm			
578 mm			
546 mm			
436 mm			
405 mm			
365 mm			
334 mm	UV/NIR region wavelength extension unit		UV/NIR region wavelength extension unit
325 mm	UV/NIR region wavelength extension unit		
313 mm			UV/NIR region wavelength extension unit
302 mm			UV/NIR region wavelength extension unit
296 mm			UV/NIR region wavelength extension unit
280 mm			UV/NIR region wavelength extension unit
254 mm			UV/NIR region wavelength extension unit



● Specifications are subject to change without notice.

JASCO INTERNATIONAL CO., LTD.

4-21, Sennin-cho 2-chome, Hachioji, Tokyo 193-0835, Japan
Tel: +81-426-66-1322 Fax: +81-426-65-6512 Internet: <http://www.jascoint.co.jp/english/>
Australia, China, Hong Kong, India, Indonesia, Iran, Korea, Malaysia, New Zealand,
Pakistan, Philippines, Russia, Singapore, South Africa, Taiwan, Thailand

JASCO INCORPORATED

8649 Commerce Drive, Easton, Maryland 21601-9903, U.S.A
Tel: +1-800-333-5272 Tel: +1-410-822-1220 Fax: +1-410-822-7526 Internet: <http://www.jascoinc.com>
Canada, Costa Rica, Mexico, Puerto Rico, Argentina, Brazil, Chile, Colombia, Paraguay, Peru, Uruguay

JASCO EUROPE s.r.l.

Via Confalonieri 25, 23894 Cremella (Lc), Italy
Tel: +39-039-956439 Fax: +39-039-958642 www.jasco-europe.com
JASCO Deutschland www.jasco.de, JASCO UK www.jasco.co.uk, JASCO France www.jascofrance.fr,
JASCO Benelux www.jasco.nl, JASCO Spain www.jasco-spain.com, JASCO Scandinavia www.jascoscandinavia.se
Austria, Finland, Greece, Hungary, Poland, Portugal, Romania, Switzerland, Algeria, Cyprus,
Egypt, Israel, Jordan, Kuwait, Lebanon, Morocco, Saudi Arabia, Syria, Tunisia, Turkey, U.A.E.



Serving the global marketplace
with Analytical Instrumentation

JASCO Corporation

2967-5, Ishikawa-cho, Hachioji, Tokyo 192-8537, Japan
<http://www.jasco.co.jp>

BP01-0702 Printed in Japan